



US006364880B1

(12) **United States Patent**  
**Michelson**

(10) **Patent No.:** **US 6,364,880 B1**  
(45) **Date of Patent:** **\*Apr. 2, 2002**

(54) **SPINAL IMPLANT WITH BONE SCREWS**

4,059,115 A 11/1977 Jumashv et al.

(76) **Inventor:** **Gary Karlin Michelson**, 438 Sherman Canal, Venice, CA (US) 90291

(List continued on next page.)

**FOREIGN PATENT DOCUMENTS**

(\*) **Notice:** This patent issued on a continued prosecution application filed under 37 CFR 1.53(d), and is subject to the twenty year patent term provisions of 35 U.S.C. 154(a)(2).

DE 106 101 10/1942  
DE 3132520 A1 6/1982

(List continued on next page.)

**OTHER PUBLICATIONS**

Subject to any disclaimer, the term of this patent is extended or adjusted under 35 U.S.C. 154(b) by 0 days.

Rathke, F.W., et al.; Surgery of the Spine; Atlas of Orthopaedic Operations, vol. 1, p. 137, W.B. Saunders Co., Philadelphia (1979).

(List continued on next page.)

(21) **Appl. No.:** **09/563,705**

*Primary Examiner*—Henry J. Recla

(22) **Filed:** **May 2, 2000**

*Assistant Examiner*—(Jackie)Tan-Uyen T. Ho

(74) *Attorney, Agent, or Firm*—Martin & Ferraro, LLP

**Related U.S. Application Data**

(57) **ABSTRACT**

(63) Continuation of application No. 09/126,585, filed on Jul. 31, 1998, now Pat. No. 6,136,001, which is a continuation of application No. 08/926,334, filed on Sep. 5, 1997, now Pat. No. 6,120,503, which is a continuation of application No. 08/589,787, filed on Jan. 22, 1996, now abandoned, which is a continuation of application No. 08/219,626, filed on Mar. 28, 1994, now abandoned.

A spinal fixation device for stabilizing one or more segments of the human spine and for preventing the dislodgement of intervertebral spinal fusion implants, which remains permanently fixated once applied. The spinal fixation device of the present invention comprises of a staple member made of material appropriate for human surgical implantation which is of sufficient length to span the disc space between two adjacent vertebrae and to engage, via essentially perpendicular extending projections, the vertebrae adjacent to that disc space. A portion of the staple of the spinal fixation device interdigitates with an already implanted intervertebral spinal fusion implant which itself spans the disc space to engage the adjacent vertebrae, and the spinal fixation device is bound to the spinal fusion implant by a locking means. The spinal fixation device of the present invention is of great utility in restraining the vertebrae adjacent to the spinal fusion implant from moving apart as the spine is extended and also serves as an anchor for a multi-segmental spinal alignment means for aligning more than one segment of the spine.


(51) **Int. Cl.**<sup>7</sup> ..... **A61B 17/56**  
(52) **U.S. Cl.** ..... **606/61**  
(58) **Field of Search** ..... 606/61, 60, 72;  
623/17, 16

(56) **References Cited**

**U.S. PATENT DOCUMENTS**

350,420 A 10/1886 Dillon  
2,065,659 A 12/1936 Cullen  
2,181,746 A 11/1939 Siebrandt  
2,774,350 A 12/1956 Cleveland, Jr.  
3,298,372 A 1/1967 Feinberg  
3,426,364 A 2/1969 Lumb  
3,875,595 A 4/1975 Froning  
3,916,907 A 11/1975 Peterson  
4,047,524 A 9/1977 Hall

**203 Claims, 10 Drawing Sheets**

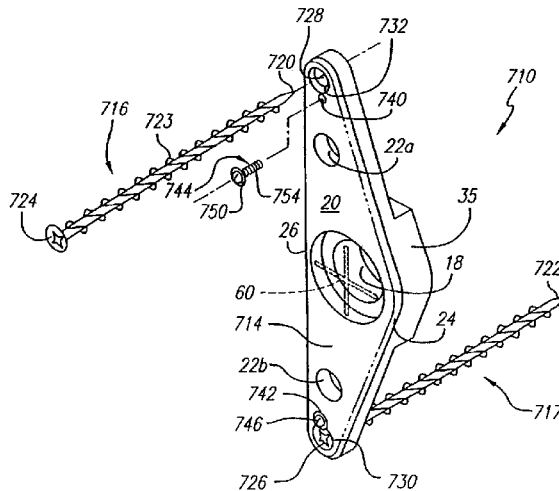

U.S. PATENT DOCUMENTS

4,271,832 A 6/1981 Evans et al.  
 4,289,123 A 9/1981 Dunn  
 4,294,753 A 10/1981 Urist  
 4,401,112 A 8/1983 Rezaina  
 4,439,152 A 3/1984 Small  
 4,455,256 A 6/1984 Urist  
 4,526,909 A 7/1985 Urist  
 4,554,914 A 11/1985 Kapp et al.  
 D281,814 S 12/1985 Pratt et al.  
 4,563,489 A 1/1986 Urist  
 4,570,623 A 2/1986 Ellison et al.  
 4,570,624 A 2/1986 Wu  
 4,592,346 A 6/1986 Jurgutis  
 4,596,574 A 6/1986 Urist  
 4,599,086 A 7/1986 Doty  
 4,604,995 A 8/1986 Stephens et al.  
 4,611,581 A 9/1986 Steffee  
 4,619,989 A 10/1986 Urist  
 4,696,290 A 9/1987 Steffee  
 4,743,256 A 5/1988 Brantigan  
 4,743,260 A 5/1988 Burton  
 4,759,769 A 7/1988 Hedman et al.  
 4,761,471 A 8/1988 Urist  
 4,789,732 A 12/1988 Urist  
 4,795,804 A 1/1989 Urist  
 4,805,602 A 2/1989 Puno et al.  
 4,857,456 A 8/1989 Urist  
 4,877,020 A 10/1989 Vich  
 4,878,915 A 11/1989 Brantigan  
 4,903,882 A 2/1990 Long  
 4,913,144 A 4/1990 Del Medico  
 4,960,420 A 10/1990 Goble et al.  
 5,015,247 A 5/1991 Michelson  
 5,026,373 A \* 6/1991 Ray et al. .... 606/61  
 5,084,048 A 1/1992 Jacob et al.  
 5,092,893 A 3/1992 Smith

5,108,422 A 4/1992 Green et al.  
 5,258,031 A 11/1993 Salib et al.  
 5,314,427 A 5/1994 Goble et al.  
 5,352,229 A 10/1994 Goble et al.  
 5,364,399 A \* 11/1994 Lowery et al. .... 606/60  
 5,395,372 A 3/1995 Holt et al.  
 5,397,364 A 3/1995 Kozak et al.  
 5,534,031 A \* 7/1996 Matsuzaki et al. .... 606/61  
 5,556,431 A 9/1996 Buttner-Janz  
 5,562,663 A 10/1996 Wisniewski et al.


FOREIGN PATENT DOCUMENTS

EP 0 077 159 4/1983  
 EP 0 179 695 4/1986  
 EP 0 307 241 3/1989  
 EP 0 551 187 A1 7/1993  
 EP 0 552 109 A1 7/1993  
 EP 0 552 109 7/1993  
 ES 2 048 671 3/1994  
 FR 2 709 410 A1 3/1995  
 FR 2 709 410 3/1995  
 GB 2082754 A 3/1982  
 JP 60-31706 2/1985  
 SU 1107854 8/1984

OTHER PUBLICATIONS



Morscher, E., et al.; Die vordere Verplattung der Halswirbelsäule mit dem Hohlschrauben-Plattensystem aus Titanium, *Der Chirurg*, vol. 57, pp. 702-707 (1986) with English Translation.  
 European Search Report dated Jul. 27, 1998 from corresponding European Application No. EP 95 91 4197.  
 PCT International Search Report dated Jul. 3, 1995 from corresponding International Application No. PCT/US95/03770.

\* cited by examiner

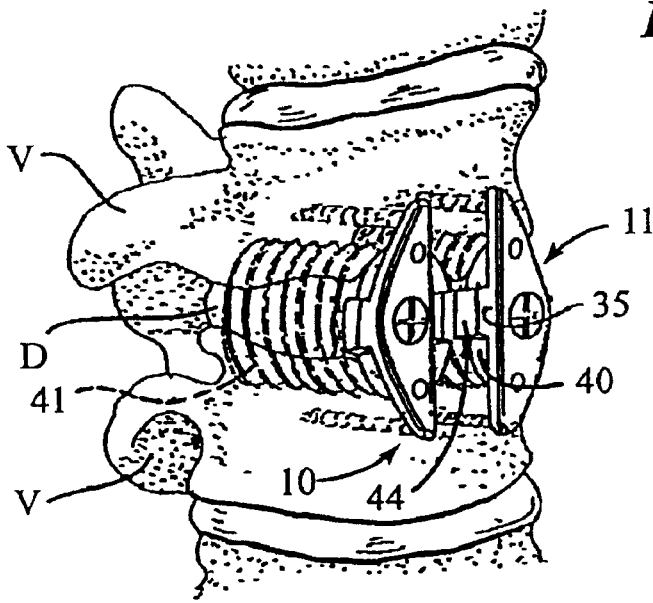


*Fig. 1*

*Fig. 2*  
PRIOR ART

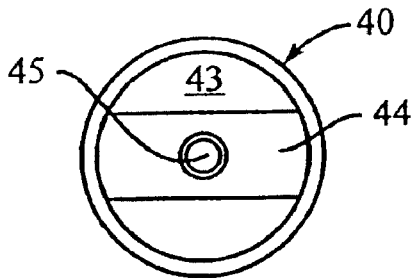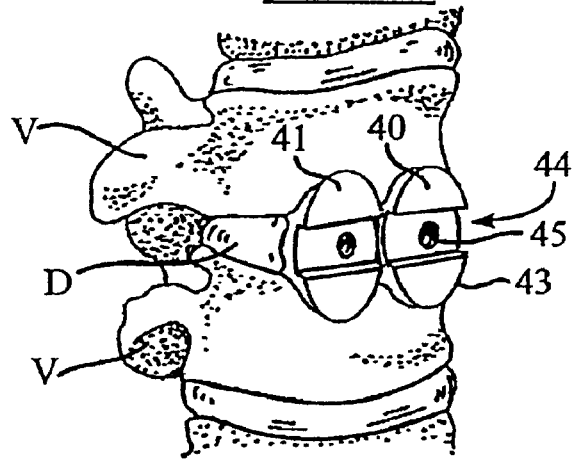

*Fig. 4*  
PRIOR ART

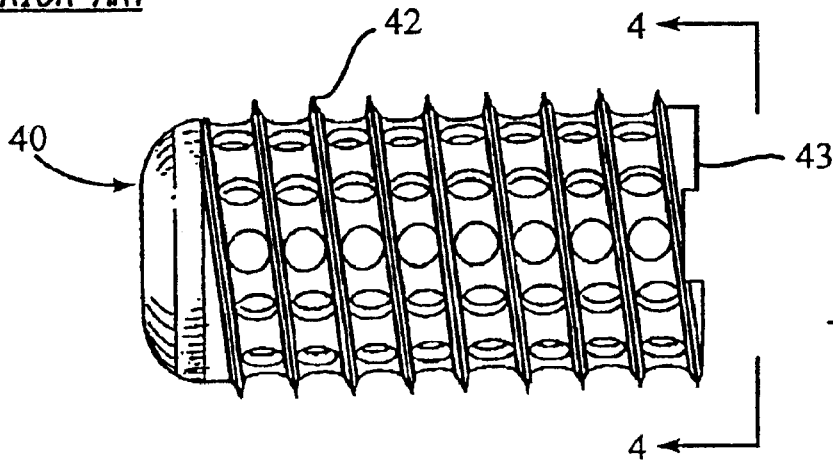

*Fig. 3*  
PRIOR ART

Fig. 5

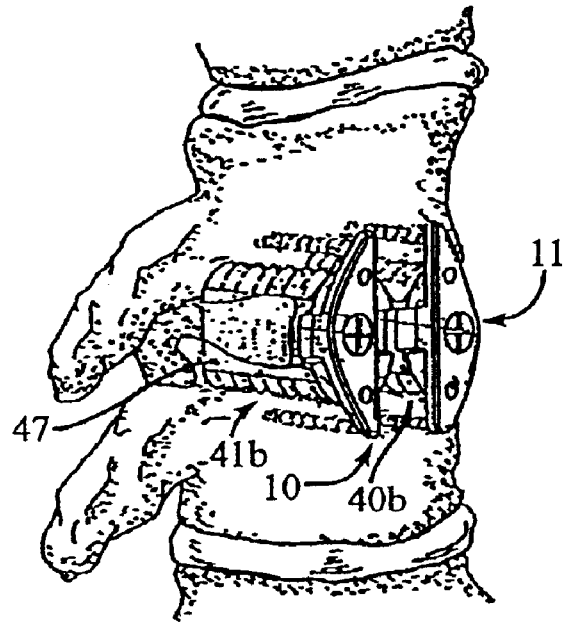
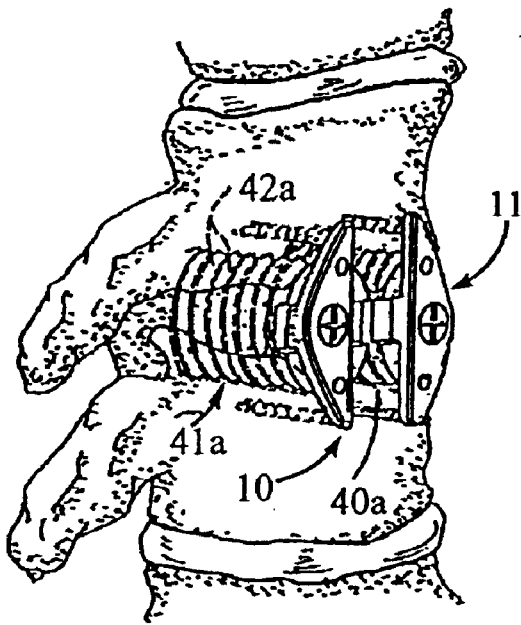


Fig. 6



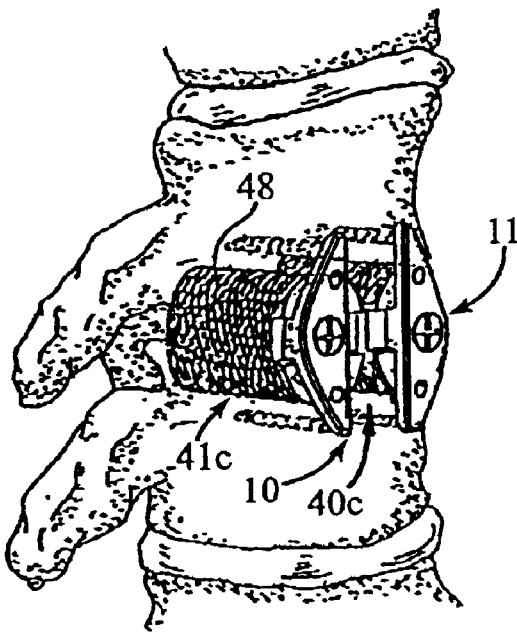



Fig. 7

Fig. 12

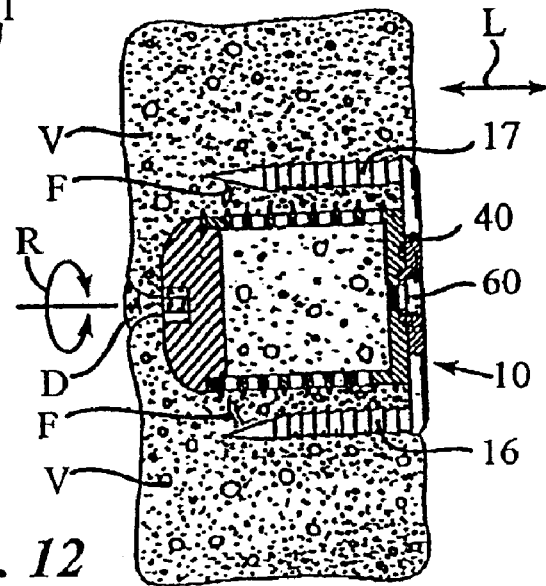


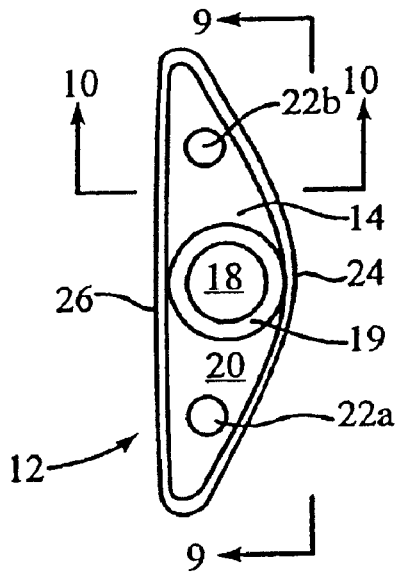


**Fig. 8**

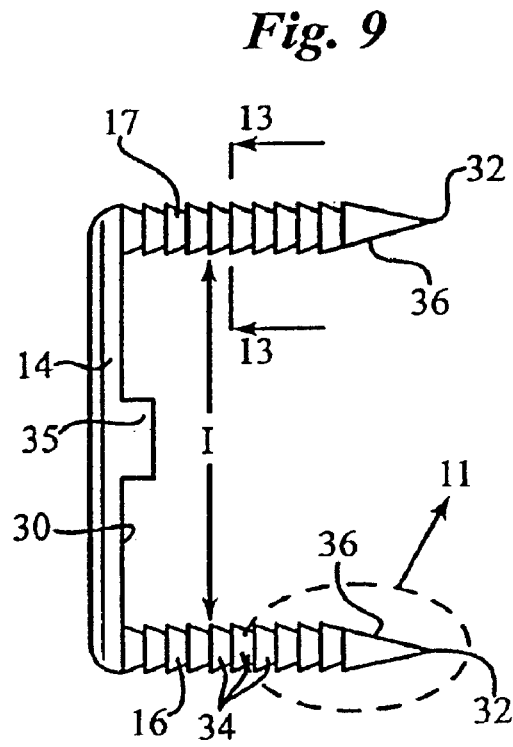


**Fig. 9**

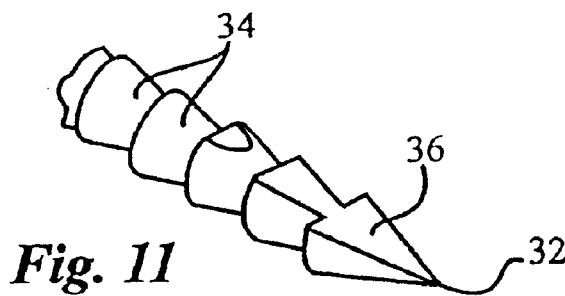


**Fig. 11**

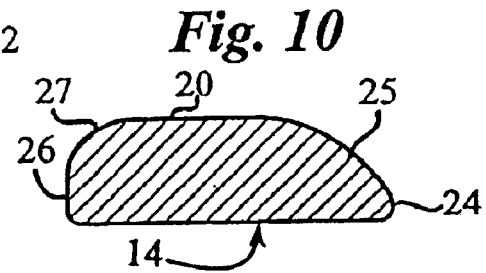


**Fig. 10**

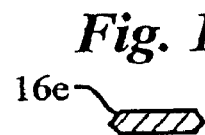

**Fig. 13F**

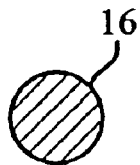

**Fig. 13A**

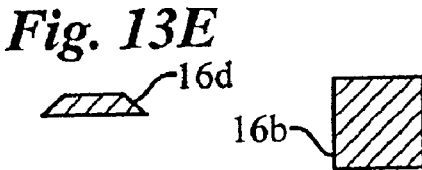

**Fig. 13E**




**Fig. 13B**


**Fig. 13C**




**Fig. 13D**

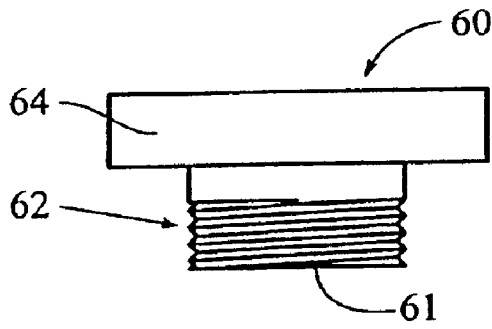


*Fig. 14*

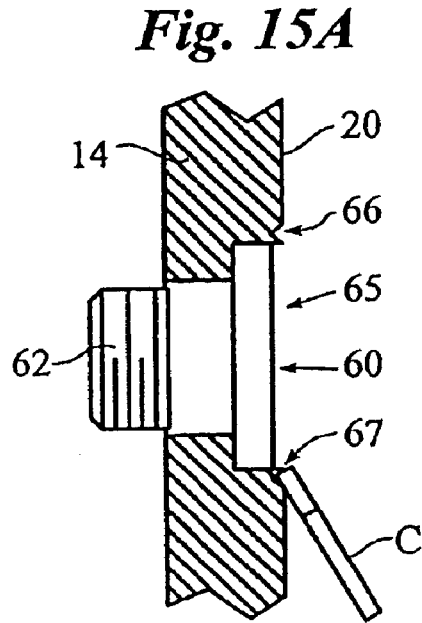


*Fig. 15A*

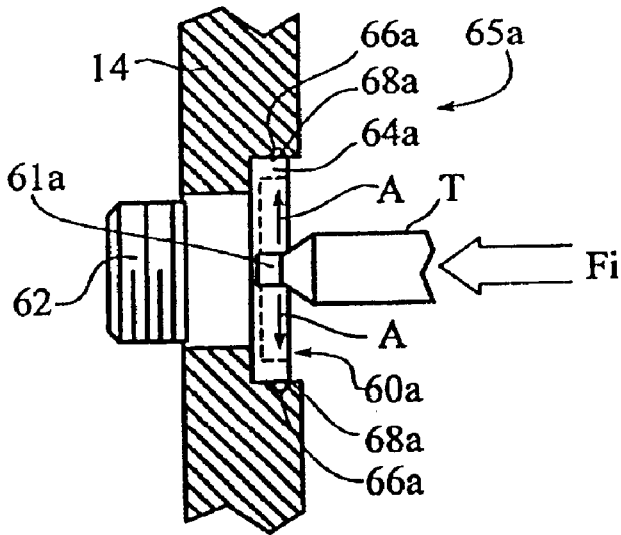


*Fig. 15B*

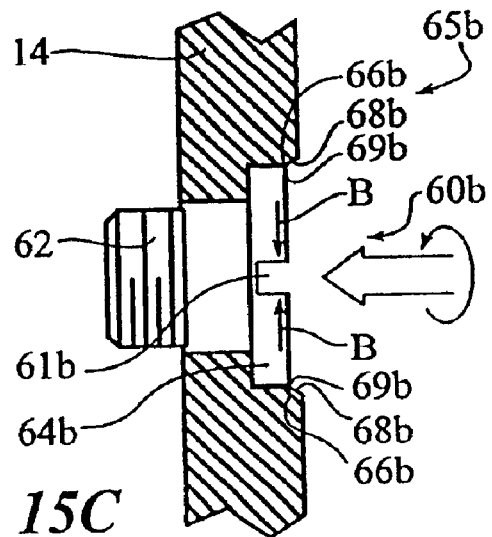


*Fig. 15C*

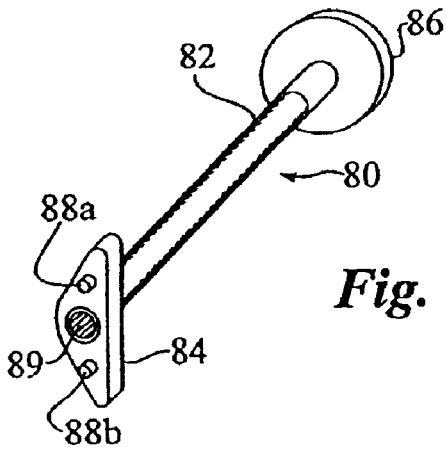


**Fig. 16A**

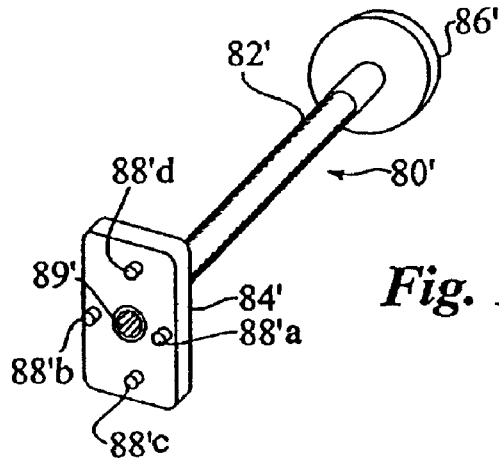


**Fig. 16B**

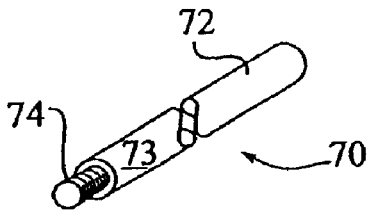

**Fig. 17A**

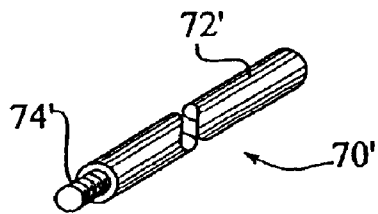

**Fig. 17B**

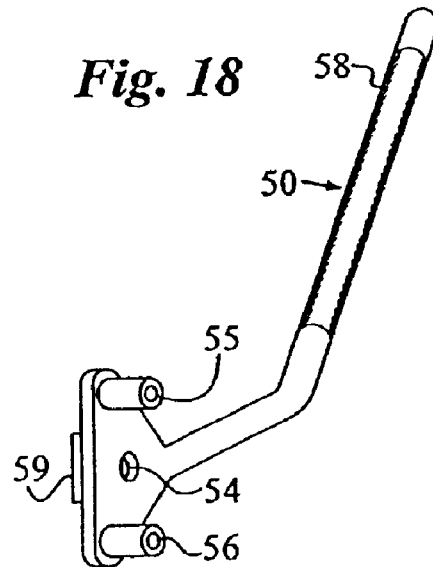

**Fig. 18**

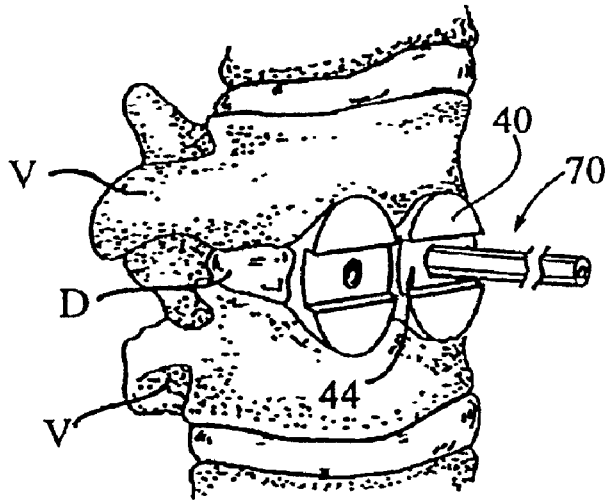


Fig. 19

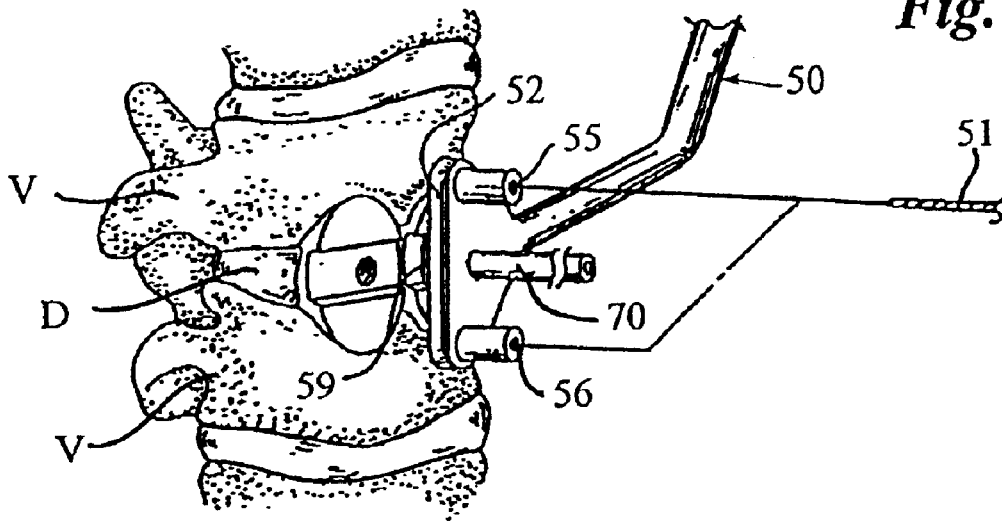


Fig. 20

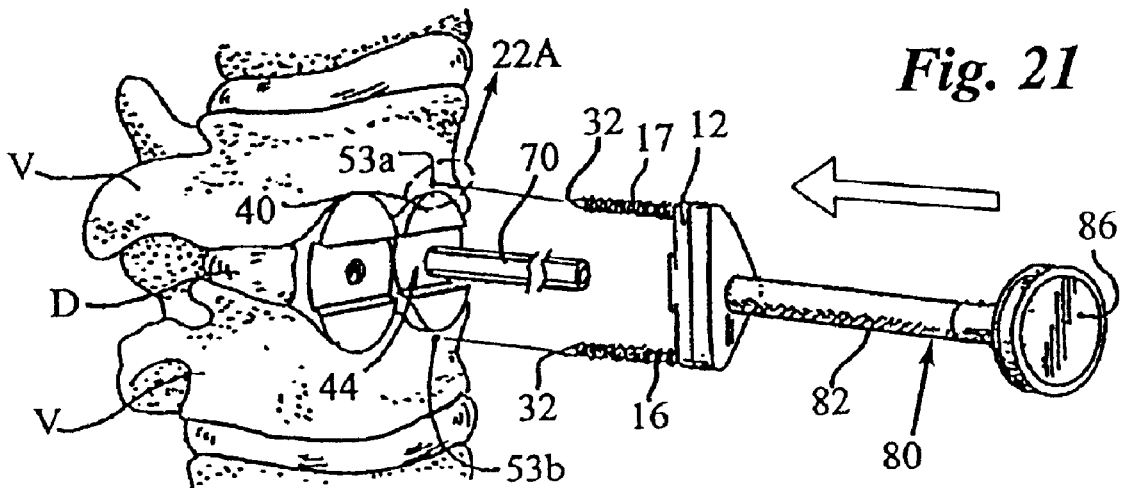


Fig. 21



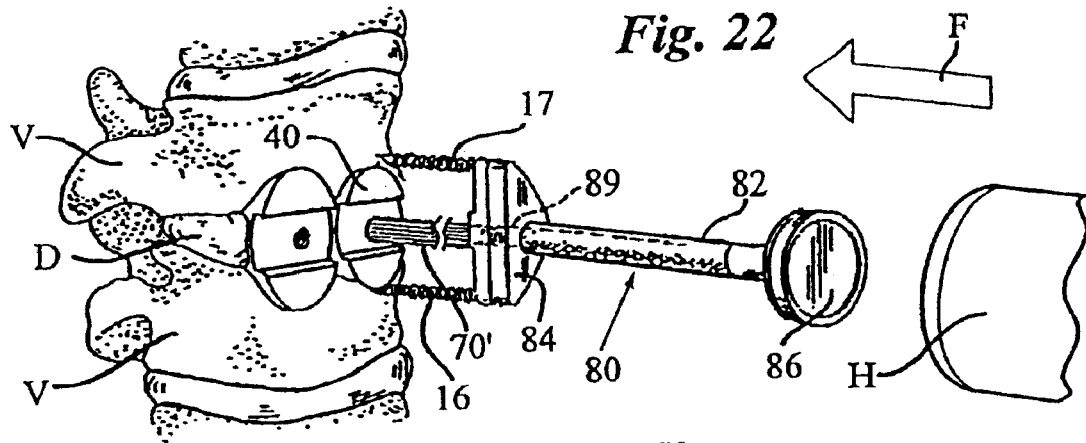


Fig. 22

Fig. 22A

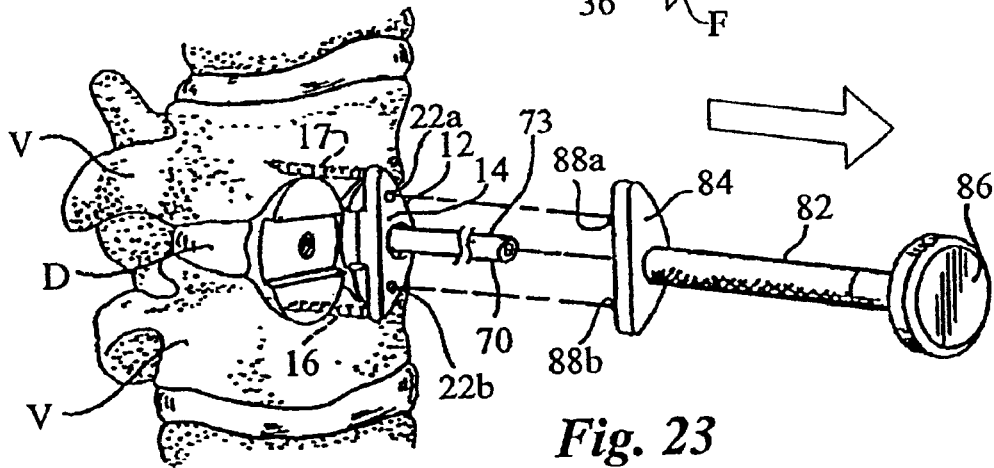
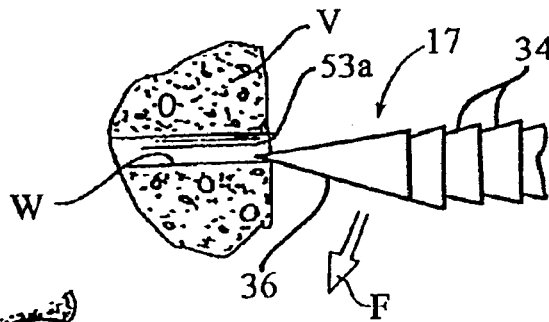


Fig. 23








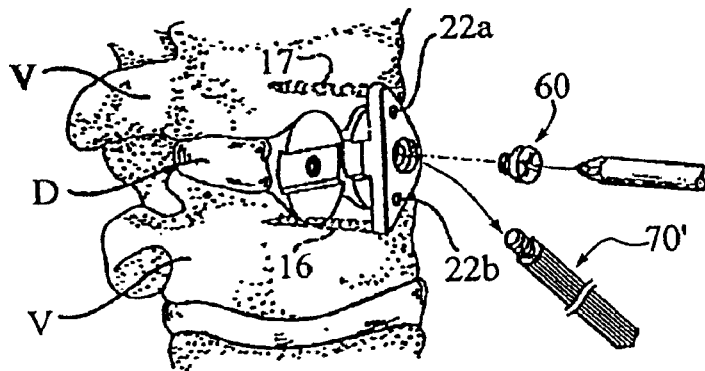


Fig. 24

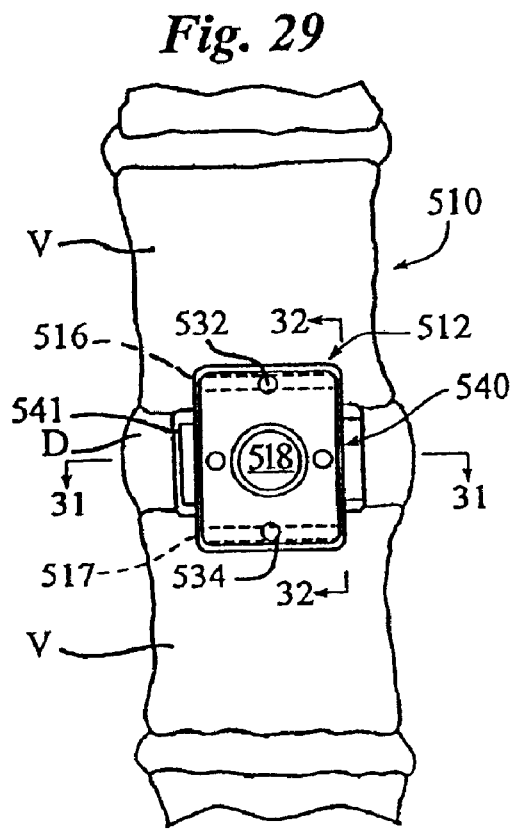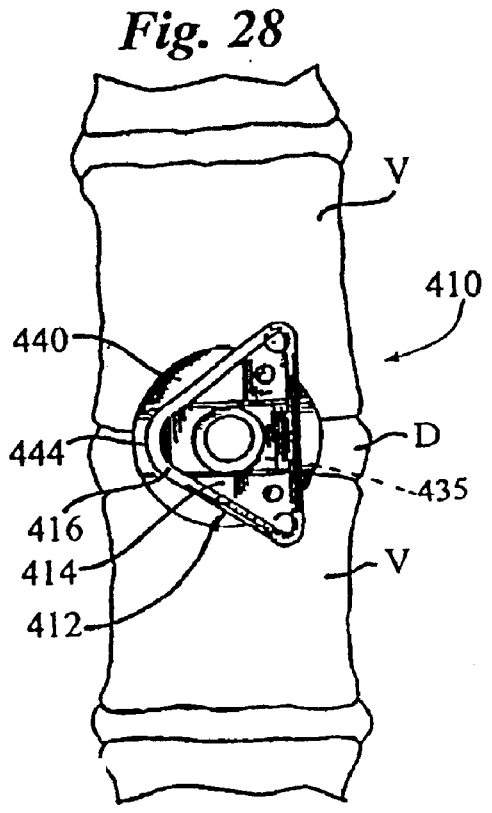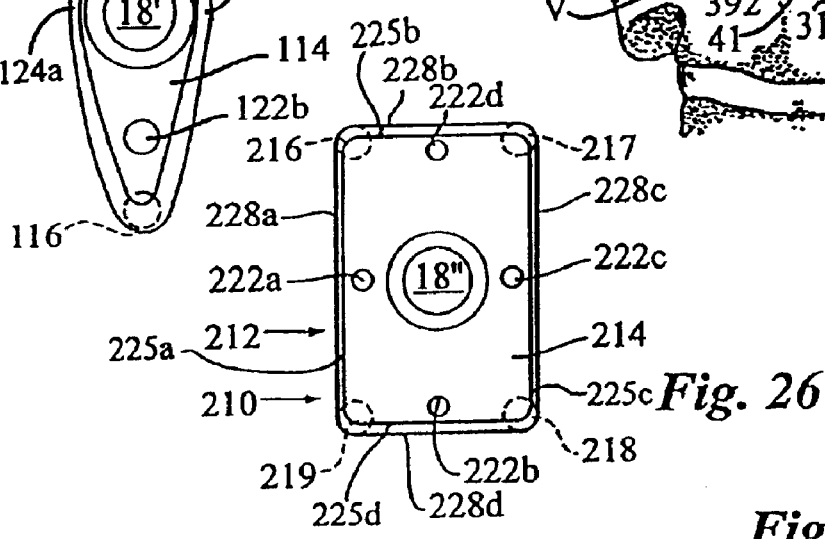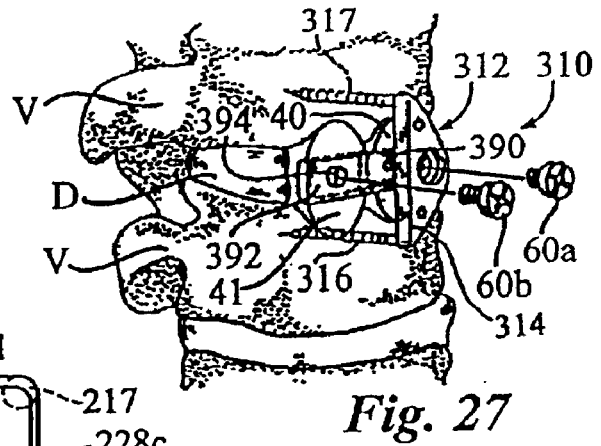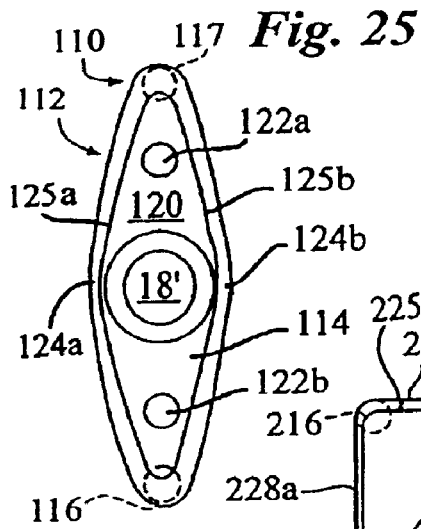

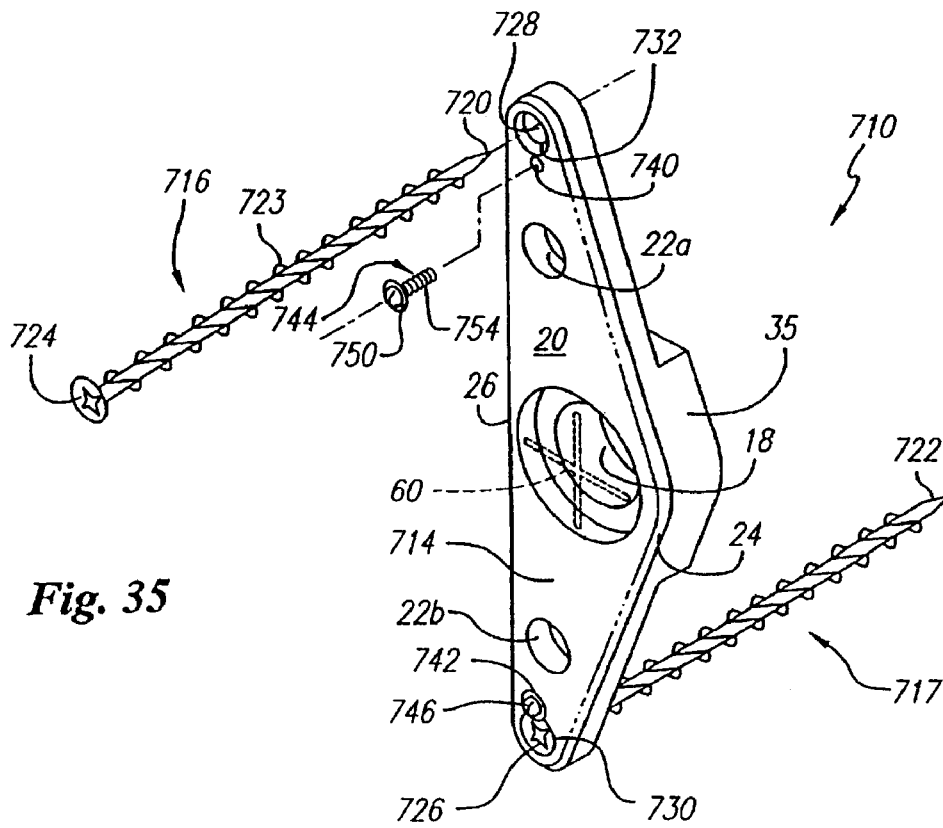

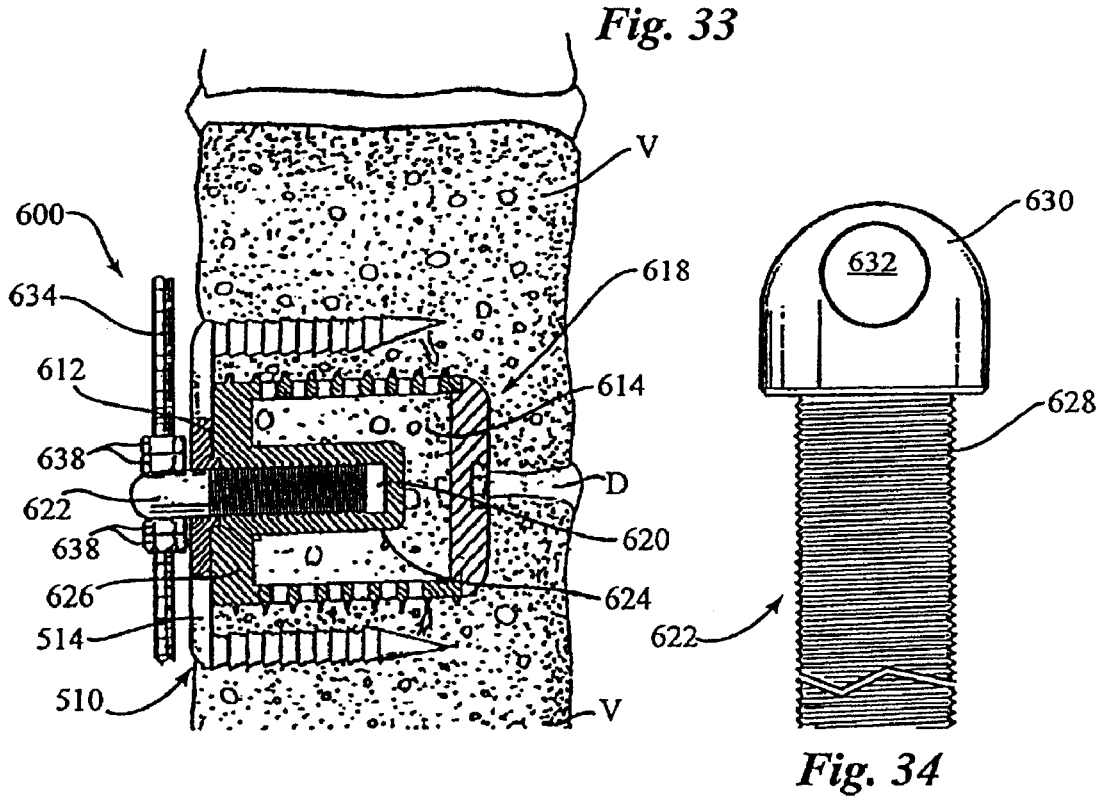

**SPINAL IMPLANT WITH BONE SCREWS**

This application is a continuation of prior application Ser. No. 09/126,585, filed Jul. 31, 1998, now U.S. Pat. No. 6,136,001, which is a continuation of application Ser. No. 08/926,334, filed Sep. 5, 1997, now U.S. Pat. No. 6,120,503 which is a continuation of application Ser. No. 08/589,787, filed Jan. 22, 1996, now abandoned which is a continuation of application Ser. No. 08/219,626, filed Mar. 28, 1994, now abandoned.

**BACKGROUND OF THE INVENTION****1. Field of the Invention**

This invention relates to surgical interbody fixation devices and in particular to a surgically implantable device for the stabilization of adjacent vertebrae of the human spine undergoing spinal arthrodesis and for the prevention of the dislodgement of spinal fusion implants used in the fusion process.

**2. Description of the Related Art**

When a segment of the human spine degenerates, or otherwise becomes diseased, it may become necessary to surgically remove the affected disc of that segment, and to replace it with bone for the purpose of obtaining a spinal fusion by which to restore more normal, pre-morbid, spatial relations, and to provide for enhanced stability across that segment. Performing such surgery of the spine from an anterior (front) approach offers the great advantage of avoiding the spinal cord, dural sac, and nerve roots. Unfortunately, in entering the disc space anteriorly a very important band-like structure called the anterior longitudinal ligament, is violated. This structure physiologically acts as a significant restraint resisting the anterior displacement of the disc itself and acting as a tension band binding the front portions of the vertebrae so as to limit spinal hyperextension.

Historically, various devices have been utilized in an attempt to compensate for the loss of this important stabilizing structure. These devices have assumed the form of blocks, bars, cables, or some combination thereof, and are bound to the vertebrae by screws, staples, boats, or some combination thereof. The earliest teachings are of a metal plate attached to adjacent vertebrae with wood-type screws. Dwyer teaches the use of a staple-screw combination. Brantigan U.S. Pat. No. 4,743,256 issued on May 10, 1988, teaches the use of a block inserted to replace the disc, affixed to a plate then screwed to the vertebrae above and below. Raezian U.S. Pat. No. 4,401,112 issued on Aug. 30, 1993, the use of a turnbuckle affixed to an elongated staple such that at least one entire vertebral body is removed, the turnbuckle portion is placed within the spine, and the staple extends both above and below the turnbuckle and engages the adjacent vertebrae to the one removed.

Unfortunately, both staples and screws have quite predictably demonstrated the propensity to back out from the vertebrae. This is quite understandable as any motion, either micro or macro, tends to stress the interface of the metallic implant to the bone, and in doing so causes the bone to relieve the high stress upon it by resorbing and moving away from the metal. This entropic change is universally from the more tightened and thus well-fixated state, to the less tightened and less fixated state. For a staple, this is specifically from the more compressed and engaged state, to the less compressed and disengaged state. Similarly, screws in such a dynamic system loosen and back out.

The potential consequences of such loosening and consequent backing out of the hardware from the anterior aspect

of the vertebral column may easily be catastrophic. Because of the proximity of the great vessels, aortic erosions and perforations of the vena cava and iliac vessels have usually occurred with unfortunate regularity and have usually resulted in death.

Therefore, the need exists for a device which is effective in restoring stability to a segment of the spine such as, but not limited to, the anterior aspect of the human spine and which will without danger remain permanently fixated once applied.

**SUMMARY OF THE INVENTION**

The present invention is directed to a spinal fixation device for stabilizing a segment of the human spine and for preventing the dislodgement of intervertebral spinal fusion implants, which remains permanently fixated to the spine once applied. The spinal fixation device of the present invention comprises a staple member made of a material appropriate for human surgical implantation and which is of sufficient length to span the disc space between two adjacent vertebrae. The staple member engages, via essentially perpendicular extending projections, the vertebrae adjacent to that disc space. The projections are sharpened and pointed so as to facilitate their insertion into the vertebrae and are segmented or ratcheted to prevent the staple member from disengaging and backing out once inserted.

In the preferred embodiment of the spinal fixation device of the present invention, a portion of the staple member interdigitates with an already implanted intervertebral spinal fusion implant and the staple member is bound to the spinal fusion implant by a locking mechanism such as a screw with a locking thread pattern. The anchoring of the staple member via a locking mechanism to a spinal fusion implant protects the patient from the danger of the staple member itself disengaging and backing out. Further, if the spinal fusion implant is externally threaded, such as the spinal fusion implant taught by Michelson, U.S. Pat. No. 5,015,247 issued on May 14, 1991, then the staple member could only back out if the spinal fusion implant were free to rotate. However, the rotation of the spinal fusion implant in this instance is blocked by its connection to the staple member which is fixated across the disc space in such a way as to be incapable of rotation. Thus, the staple member is made safe against dislodgement by attachment to the spinal fusion implant and the stability of the spinal fusion implant is assured as it is also stabilized by the staple member and each works in connection with the other to remove the only remaining degree of freedom that would allow for the disengagement of either.

The spinal fixation device of the present invention is broadly applicable to the anterior, posterior and lateral aspects of the spinal column, be it the cervical, thoracic or lumbar area. In particular, the use of a staple member spanning the disc space and engaging the adjacent vertebrae which is applied to the anterior aspect of the spine is of great utility in restraining those vertebral bodies from moving apart as the spine is extended and thus is effective in replacing the anterior longitudinal ligament of the patient.

The spinal fixation device of the present invention provides the advantage of facilitating cross vertebral bony bridging (fusion via immobilization) which when achieved relieves all of the forces on the inserted spinal fusion implants. The spinal fixation device of the present invention may be coated with materials, to promote bone fusion and thus promote the incorporation and ultimate entombment of the spinal fixation device into the bone fusion mass. The use

of a bone fusion promoting material results in a speedier vertebra to vertebra fusion as bone may grow along the coated spinal fixation device bridging the two vertebrae so that the spinal fixation device acts as a trellis and supplies essential chemical elements to facilitate the bone fusion process.

Another advantage provided by the spinal fixation device of the present invention is that as it is inserted it compresses the adjacent vertebrae together, thus increasing the compressive load on the spinal fusion implants or implants within the disc space, such compression being beneficial to fusion and further stabilizing the spinal fusion implants.

A further advantage of the spinal fixation device of the present invention is that it may be used as an anchor such that a multiplicity of spinal fixation devices may then be interconnected via a cable, rod, bar, or plate, so as to achieve or maintain a multi-segmental spinal alignment.

Alternatively, the spinal fixation device of the present invention could be made of resorbable materials, such as biocompatible resorbable plastics, that resorb at an appropriate rate such that once the spinal fixation device is no longer needed (i.e. when spinal fusion is complete) the body would resorb the spinal fixation device. The spinal fixation device could be only in part resorbable such that the projections of the staple member would be non-resorbable and would remain incarcerated in the vertebrae and sealed off once the resorbable portion of the staple is resorbed by the body.

As a further alternative, the spinal fixation device of the present invention could be made wholly of in part of ceramic and more particularly made of or coated with a ceramic such as hydroxyapatite that would actively participate in the fusion process.

#### OBJECTS OF THE PRESENT INVENTION

It is an object of the present invention to provide a spinal fixation device having a staple member spanning the disc space and engaging two adjacent vertebrae of the spine to restrain the vertebrae from moving apart as the spine is extended;

It is another object of the present invention to provide a spinal fixation device that is effective in replacing the function of the anterior longitudinal ligament of a patient;

It is a further object of the present invention to provide a means for protecting the patient from the danger of the spinal fixation device itself disengaging and backing out by its being anchored to an intervertebral spinal fusion implant;

It is still another object of the present invention to provide a spinal fixation device that blocks the rotation of an intervertebral spinal fusion implant by its connection to the staple member which is fixated across the disc space in such a way as to be incapable of rotation thereby preventing the spinal fusion implant from backing out;

It is yet another object of the present invention to provide a spinal fixation device that is broadly applicable to the anterior aspect of the spinal column, be it the cervical, thoracic or lumbar area;

It is another object of the present invention to provide a spinal fixation device which may be applied longitudinally at any point about the circumference of the anterior aspect of the spine;

It is also another object of the present invention to provide a spinal fixation device that stabilizes a surgically implanted spinal fusion implant and works in connection with the spinal fusion implant to prevent disengagement of either;

It is another object of the present invention to provide a spinal fixation device that achieves cross vertebral bony bridging (fusion) which ultimately relieves all of the forces on intervertebral spinal fusion implants inserted within the disc space between two adjacent vertebrae, and provides for a permanently good result;

It is another object of the present invention to provide a spinal fixation device that serves as an anchor, such that a multiplicity of these anchors may then be interconnected via a cable, rod, bar, or plate, so as to achieve or maintain a multi-segmental spinal alignment; and

It is a further object of the present invention to provide a spinal fixation device that directly participates in the bony bridging of two adjacent vertebrae and participates in the spinal fusion process across those vertebrae.

These and other objects of the present invention will become apparent from a review of the accompanying drawings and the detailed description of the drawings.

#### BRIEF DESCRIPTION OF THE DRAWINGS

FIG. 1 is a perspective side view of a segment of the spinal column having two spinal fusion implants shown partially in hidden line inserted across the disc space between two adjacent vertebrae with each spinal fusion implant having a spinal fixation device of the present invention shown partially in hidden line secured thereto, spanning across the disc space and inserted into the vertebrae.

FIG. 2 is a perspective side view of a segment of the spinal column having two spinal fusion implants inserted across the disc space between two adjacent vertebrae.

FIG. 3 is an elevational side view of a cylindrical threaded spinal fusion implant.

FIG. 4 is an end view of the cylindrical threaded spinal fusion implant along lines 4—4 of FIG. 3.

FIG. 5 is a perspective side view of a segment of the spinal column having two non-threaded spinal fusion implants with external ratchetings, shown in hidden line, inserted across the disc space between two adjacent vertebrae with each spinal fusion implant having a spinal fixation device of the present invention, shown partially in hidden line, coupled thereto, spanning across the disc space and inserted into the vertebrae.

FIG. 6 is a perspective side view of a segment of the spinal column having two spinal fusion implants having truncated sides with external ratchetings shown in hidden line inserted across the disc space between two adjacent vertebrae with each spinal fusion implant having a spinal fixation device of the present invention shown partially in hidden line coupled thereto, spanning across the disc space and inserted into the vertebrae.

FIG. 7 is a perspective side view of a segment of the spinal column having two spinal fusion implants having a knurled external surface shown in hidden line inserted across the disc space between two adjacent vertebrae with each spinal fusion implant having a spinal fixation device of the present invention shown partially in hidden line coupled thereto, spanning across the disc space and inserted into the vertebrae.

FIG. 8 is a top plan view of the spinal fixation device of the present invention.

FIG. 9 is a side view of the spinal fixation device of the present invention along lines 9—9 of FIG. 8.

FIG. 10 is a cross sectional view taken along lines 10—10 of FIG. 8 showing the top member of the spinal fixation device of the present invention.

FIG. 11 is an enlarged fragmentary perspective side view of a projection of the spinal fixation device of the present invention taken along line 11 of FIG. 9.

FIG. 12 is a cross sectional view of the spinal fixation device of the present invention inserted into the vertebrae and secured to the spinal fusion implant with the arrows showing the forces exerted, the rotational axis and the longitudinal axis of the spinal fusion implant.

FIG. 13A is a cross sectional view along line 13—13 of FIG. 9 of the preferred embodiment of the projections of the present invention.

FIGS. 13B, 13C, 13D, 13E, and 13F are cross sectional views taken along line 13—13 of FIG. 9 showing alternative embodiments of the projections of the spinal fixation device of the present invention.

FIG. 14 is an enlarged elevational side view of the locking screw used to secure the spinal fixation device of the present invention to a spinal fusion implant.

FIG. 15A is a cross sectional view of a securing means for locking the locking screw of the present invention.

FIG. 15B is a cross sectional view of a first alternative embodiment of the securing means for locking the locking screw of the present invention.

FIG. 15C is a cross sectional view of a second alternative embodiment of the securing means for locking the locking screw of the present invention.

FIG. 16A is a perspective side view of the instrumentation used for driving the spinal fixation device of the present invention into the vertebrae.

FIG. 16B is a perspective side view of a first alternative embodiment of the instrumentation used for driving the spinal fixation device of the present invention into the vertebrae.

FIG. 17A is a perspective side view of an alignment rod used to align the spinal fixation device of the present invention.

FIG. 17B is a perspective side view of an alternative embodiment of the alignment rod having splines used to align the spinal fixation device of the present invention.

FIG. 18 is a front perspective view of the drill template instrument.

FIG. 19 is a perspective side view of the alignment rod attached to a spinal fusion implant inserted in the disc space between two adjacent vertebrae.

FIG. 20 illustrates the step of drilling guide holes in the vertebrae adjacent to the spinal fusion implant with the drill template instrument of FIG. 18.

FIG. 21 illustrates a step of the method of inserting the spinal fixation device of the present invention with the alignment rod attached to the spinal fusion implant and the spinal fixation device placed on the driver instrumentation.

FIG. 22 illustrates a step of the short method of inserting the spinal fixation device of the present invention with the driver instrument engaging the splined alignment rod and a hammer for applying an impaction force and driving the driver instrument.

FIG. 22A is an enlarged fragmentary view of a projection being inserted into an insertion hole drilled within a vertebra shown in cross section taken along line 22A of FIG. 21.

FIG. 23 illustrates another step of the method of inserting the spinal fixation device of the present invention in which the spinal fixation device has been driven into the vertebrae and the driver instrumentation has been removed.

FIG. 24 illustrates another step of the method of inserting the spinal fixation device of the present invention with the

splined alignment rod being removed from the spinal fusion implant and the locking screw being inserted and secured the spinal fixation device to the spinal fusion implant.

FIG. 25 is a top plan view of a first alternative embodiment of the spinal fixation device of the present invention.

FIG. 26 is a top plan view of a second alternative embodiment of the spinal fixation device of the present invention.

FIG. 27 is a perspective side view of a third alternative embodiment of the spinal fixation device of the present invention coupled to two spinal fusion implants and inserted in adjacent vertebrae of the spinal column.

FIG. 28 is a top plan view of a fourth alternative embodiment of the spinal fixation device of the present invention inserted into the vertebrae of the spinal column having a spinal fusion implant inserted in the disc space.

FIG. 29 is a top plan view of a fifth alternative embodiment of the spinal fixation device of the present invention inserted into the vertebrae of the spinal column having a spinal fusion implant inserted in the disc space.

FIG. 30 is a perspective bottom view of the fourth alternative embodiment of the spinal fixation device of the present invention.

FIG. 31 is a cross sectional view along lines 31—31 of FIG. 29 showing the fifth alternative embodiment of the spinal fixation device of the present invention inserted into the adjacent vertebrae and coupled to a spinal fusion implant.

FIG. 32 is a cross sectional view along lines 32—32 of FIG. 29 showing the projections of the fifth alternative embodiment of the present invention with respect to a spinal fusion implant inserted within the disc space.

FIG. 33 is a cross sectional view of a spinal fixation device of the present invention engaging two adjacent vertebrae and being attached to a spinal fusion implant, shown being used as an anchor for a multi-segmental spinal alignment means.

FIG. 34 is an enlarged elevational side view of a threaded post used to connect the spinal fixation device of the present invention to a multi-segmental spinal alignment means.

FIG. 35 is an exploded perspective view of a sixth alternative embodiment of the spinal fixation device of the present invention having independent projection members that are screws.

#### DETAILED DESCRIPTION OF THE DRAWINGS

Referring to FIGS. 1 and 2, two identical spinal fixation devices of the present invention, each being generally referred to by the numerals 10 and 11, respectively, are shown inserted into two vertebrae V adjacent to a disc D of a segment of the human spine. Each spinal fixation device 10 and 11 is shown coupled to identical spinal fusion implants 40 and 41 that have been surgically implanted in the disc space between adjacent vertebrae V. In this manner, the spinal fixation devices 10 and 11 stabilize a segment of the spine, prevent the dislodgement of the spinal fusion implant 40, and remain permanently fixated to the spine once applied. The spinal fixation devices 10 and 11 are identical such that the description of one is equally applicable to the other. Thus, the description that follows will be directed to spinal fixation device 10.

Referring to FIGS. 3—4, the spinal fusion implant 40 such as, but not limited to, the spinal fusion implant described by Michelson, U.S. Pat. No. 5,015,247 issued on May 14, 1991, is shown. The spinal fusion implant 40 is cylindrical in shape

and has external threads **42** at its outer perimeter for engaging the bone of the vertebrae V adjacent to the disc D. The spinal fusion implant **40** has an insertion end **43** having a depression **44** and a threaded aperture **45** for engaging a portion of the spinal fixation device **10** and also for engaging a portion of an instrument used to insert the spinal fixation device **10** into the vertebrae V.

Referring to FIGS. 5–7, it is appreciated that the spinal fixation devices **10** and **11** of the present invention are not limited in use with a threaded spinal fusion implant **40** and **41**, but may be used with different types of spinal fusion implants. For example, the spinal fixation devices **10** and **11** may be coupled to spinal fusion implants **40a** and **41a**, respectively, each having external ratchetings **42a** instead of external threads **42** as shown in FIG. 5. Alternatively, the spinal fixation devices **10** and **11** may be coupled to spinal fusion implants **40b** and **41b**, respectively, each having a partially cylindrical shape with at least one truncated side **47** as shown in FIG. 6. As a further alternative, the spinal fixation devices **10** and **11** may be coupled to spinal fusion implants **40c** and **41c**, respectively, each having a knurled external surface **48** as shown in FIG. 7. It is also appreciated that the spinal fixation devices may be used with a variety of other bone fusion implants without departing from the scope of the present invention.

Referring to FIGS. 8–9, in the preferred embodiment, the spinal fixation device **10** of the present invention comprises a staple member **12** having a substantially planar top member **14** which is of sufficient length to span one intervertebral disc D and to engage, via a plurality of essentially perpendicular extending projections **16** and **17**, the vertebrae V adjacent to that disc D. The top member **14** has a central opening **18** within a concentric, countersunk recess **19** for receiving therethrough a screw or similar coupling means for coupling the spinal fixation device **10** to the spinal fusion implant **40**. The top member **14** has an upper surface **20** having a pair of openings **22a** and **22b** for receiving the posts **88a** and **88b** of a driving instrument **80** which is described in greater detail below in reference to FIGS. 16A and 16B.

Referring to FIG. 10, a cross sectional view of the top member **14** is shown. In the preferred embodiment, the top member **14** is generally triangularly shaped and is radiused along curved side **24** and straight side **26**. The curved side **24** of the top member **14** is radiused at its upper edge **25** and at the upper edge **27** of straight side **26** to conform to the external curvature of the vertebrae V. In this manner, smooth surfaces are created at the upper edges **25** and **27** of the top member **14** that are contoured to the shape of the external curvature of the vertebrae V when the staple member **12** is in place. The smooth contoured surface of the upper edges **25** and **27** of the top member **14** prevent aortic erosions and perforations of the vessels proximate the vertebral column such as the vena cava and the iliac vessels which might otherwise result from friction.

In the preferred embodiment of the spinal fixation device **10**, the top member **14** has a width ranging from 6.0 mm to 28.0 mm, with 10.0 mm being the preferred width, and having a thickness in the range of 2.0 mm to 4.0 mm, with 3.0 mm being the preferred thickness. The staple member **12** is made of material appropriate for human surgical implantation including all surgically appropriate metals such as but not limited to, titanium, titanium alloy, chrome molybdenum alloys, stainless steel; or non-metallic materials including permanent or resorbable substances or composites, carbon fiber materials, resins, plastics, ceramics or others.

Further, the staple member **12** of the present invention may be treated with, or even composed of, materials known

to participate in or promote in the fusion process or bone growth. The spinal fixation device **10** may be coated with materials to promote bone fusion and thus promote the incorporation and ultimate entombment of the spinal fixation device **10** into the bone fusion mass. The use of a bone fusion promoting material such as, but not limited to hydroxyapatite, hydroxyapatite tricalcium phosphate or bone morphogenic protein, results in a speedier vertebra V to vertebra V fusion as bone may grow along the coated spinal fixation device **10** bridging the two vertebrae V so that the spinal fixation device **10** acts as a trellis and supplies essential chemical elements to facilitate the bone fusion process.

Referring again to FIG. 9, the projections **16** and **17** are positioned at opposite ends of the top member **14** and depend downwardly and extend perpendicularly from the bottom surface **30** of the top member **14**. The projections **16** and **17** each terminate in a distal end **32** that is pointed and sharpened to facilitate the insertion of the projections **16** and **17** into the vertebrae V.

The staple member **12** is most effective when the inter-projection distance I between projections **16** and **17** is at least 4.0 mm and preferably 6.0 mm greater than the diameter of the particular spinal fusion implant **40** for which the spinal fixation device **10** is being used so that at least 2.0 mm and preferably 3.0 mm of bone from the vertebrae V will be present between the spinal fusion implant **40** and each of the projections **16** and **17**. Typically, intervertebral spinal fusion implants have a diameter that ranges from 12.0 mm to 28.0 mm, therefore, the interprojection distance I typically will range from 18.0 mm to 34.0 mm for most applications.

In the preferred embodiment, the projections **16** and **17** comprise a series of segmented and ratcheted portions **34**. The segmented and ratcheted portions **34** provide for a “one way” insertion of the staple member **12** to prevent the backing-out of the projections **16** and **17** once they are inserted into the bone of the vertebrae V. In the preferred embodiment, each segmented and ratcheted portion **34** of the projections **16** and **17** is conical in shape and the diameter of each segmented and ratcheted portion **34** increases in the direction from the distal end **32** toward the top member **14** so that the projections **16** and **17** resemble a stack of cones. The segmented and ratcheted portions **34** are spaced approximately 2.0 mm to 4.0 mm apart, with 3.0 mm being the preferred distance between each segmented and ratcheted portion **34**.

Referring to FIG. 11–12, in the preferred embodiment of the spinal fixation device **10**, in order to further facilitate the insertion of the projections **16** and **17** into the vertebrae V, the distal end **32** of each projection **16** has an eccentric, incline-planed inner surface **36** as shown in FIG. 11. The eccentric, incline-planed inner surface **36** of each of the projections **16** and **17** create a force F which pushes the bone of the vertebrae V toward the spinal fusion implant **40** as the staple member **12** is inserted into each of the vertebrae V as shown in FIG. 12.

Referring to FIGS. 13A–13F, in the preferred embodiment of the spinal fixation device **10**, the projections **16** and **17** are cylindrical in shape having a circular cross section as shown for projection **16** in FIG. 13A. Alternatively, the projection **16a** may have a triangular cross section as shown in FIG. 13B; the projection **16b** may have a square cross section as shown in FIG. 13C; the projection **16c** may have a rectangular cross section as shown in FIG. 13D; the projection **16d** may have a trapezoidal cross section as



shown in FIG. 13E; or the projection 16e may have a cross section with a configuration as shown in FIG. 13F.

In the preferred embodiment, the projections 16 and 17 each have a diameter of approximately 2.0 mm to 4.0 mm, with 3.0 mm being the preferred diameter at the widest point. The projection 16 and 17 each have a length ranging from 16.0 mm to 28.0 mm, with 22.0 mm being the preferred length when the spinal fixation device 10 is implanted in the direction of the anterior aspect of the vertebra V to the posterior aspect of the vertebrae V. Alternatively, it is appreciated that the projections 16 and 17 each could have a longer length depending on the diameter of the vertebrae V in which the projections 16 and 17 are implanted.

Referring again to FIG. 9, the top member 14 of the staple member 12 has a central bar 35 extending from the center of its bottom surface 30, for interdigitating and mating to an already implanted intervertebral spinal fusion implant 40. In the preferred embodiment, the central bar 35 has a thickness in the range of 0.5 mm to 1.5 mm, with 0.5 mm being the preferred thickness.

Referring to FIG. 1, the central bar 35 is configured so that it complements and engages the depression 44 at the insertion end 43 of the spinal fusion implant 40. Once engaged to the depression 44, the bar 35 interdigitates with the depression 44 of the spinal fusion implant 40 to lock and prevent the rotation of the spinal fusion implant 40.

Referring to FIG. 14, in the preferred embodiment, the staple member 12 is secured to the spinal fusion implant 40 by a screw 60 having threaded end 61 with a locking thread pattern 62 and screw head 64. The locking thread pattern 62 has a reduced pitch at the bottom of the threaded end 61 such that the screw 60 is self-locking. However, it is appreciated that the threaded pattern 62 may be any of the means for locking a screw well known by those skilled in the art.

Referring to FIGS. 2 and 8, the threaded end 61 of the screw 60 passes through the central opening 18 of the top member 14 and the threaded pattern 62 threads into the threaded aperture 45 of the spinal fusion implant 40. The screw head 64 fits within the countersunk recess 19 of the top member 14 such that the screw head 64 is at or below the plane of the upper surface 20 of the top member 14. In the preferred embodiment, the central opening 18 has a diameter ranging from 4.5 mm to 5.5 mm, with 5.0 mm being the preferred diameter. The countersunk recess 19 has a diameter in the range of 6.0 mm to 8.0 mm with 7.0 mm being the preferred diameter.

Referring to FIGS. 15A, 15B, and 15C, an enlarged cross sectional view of three different embodiments of a securing means 65 for locking the screw 60 once it is threaded to the spinal fusion implant 40 are shown. In FIG. 15A, the securing means 65 comprises a notch 66 in the surface 20 of the top member 14 which is preferably made of metal. Once the screw 60 is threaded and securely tightened to the spinal fusion implant 40, a chisel C is used to bend a portion 67 of the top member 14 into the central opening 18 and against the screw head 64 so as to prevent the outward excursion and any unwanted loosening of the screw 60.

In FIG. 15B, a second embodiment of the securing means 65a is shown comprising a central score 66a concentric with the central opening 18. A screw 60a having a slot 61a in the screw head 64a is threaded and securely tightened to the spinal fusion implant 40. An instrument T is partially inserted into slot 61a after which an impaction force  $F_i$  is applied to the instrument T to spread apart the screw head 64a in the direction of the arrows A so that the screw head 64a becomes deformed from the impaction force  $F_i$  and fits

within the central score 66a. Once the screw head 64a is in the central score 66a, the outward excursion of the screw 60a is prevented by the top lip 68 of the central score 66a.

In FIG. 15C, a third embodiment of the securing means 65b is shown comprising a screw 60b having a screw head 64b with a slightly flanged portion 69b near the top and a slot 61b. The central opening 18 has along its circumference a recess 66b for receiving the flanged portion 69b of the screw head 64b. The securing means 65b relies on the natural resiliency of the metal screw head 64b such that when the screw 60b is being driven by a screw driver, the screw head 64b flexes in the direction of the arrows B. In this manner, the flanged portion 69b of the screw head 64b slides along the interior of the central opening 18 so that the screw head 64b is below the top lip 68b of the recess 66b. Once the screw driver is removed from the screw 60b, the screw head 64b returns to its natural state in the direction opposite to the arrows B so that the flanged portion 69b is within the recess 66b. The outward excursion of the screw 60 is thus prevented by the top lip 68b which blocks the screw head 64b by catching the flanged portion 69b.

FIGS. 16A–18 show the instrumentation used for installing the spinal fixation device 10. Referring to FIG. 16A, a driving instrument 80 used for inserting the spinal fixation device 10 into the vertebrae V is shown having a hollow tubular shaft 82 which terminates at one end to a bottom flat member 84 and terminates to a top flat member 86 at the other end. The bottom flat member 84 is preferably configured so that it conforms to the shape of the top member 14 of the staple member 12.

The driving instrument 80 has a pair of short posts 88a and 88b extending from the bottom flat member 84. The posts 88a and 88b are oriented on the bottom flat member 84 so as to correspond to the position of the openings 22a and 22b in the upper surface 20 of the top member 14 of the staple member 12. Each of the posts 88a and 88b fit into each of the openings 22a and 22b and keep the staple member 12 aligned on the bottom flat member 84 of the driving instrument 80. It is appreciated that the openings 22a and 22b in the top member 14 may be depressions within the surface 20 of the top member 14 or may be holes that pass through the top member 14. In the preferred embodiment, the openings 22a and 22b gave a diameter ranging from 1.5 mm to 3.5 mm, with 2.5 mm being the preferred diameter.

Referring to FIG. 16B, an alternative embodiment of the driving instrument 80' which is, used for inserting into the vertebrae V the spinal fixation device 210, described in detail below in reference to FIG. 26, is shown having a hollow tubular shaft 82' which terminates at one end to a bottom flat member 84' and terminates to a top flat member 86' at the other end. The bottom flat member 84' is rectangular in shape so that it conforms to the shape of the top member 214 of the spinal fixation device 210.

The driving instrument 80' has a pair of short posts 88'a, 88'b, 88'c and 88'd extending from the bottom flat member 84'. The posts 88'a–88'd are oriented on the bottom flat member 84' so as to correspond to the position of the openings 222a–222d of the spinal fixation device 210. Each of the and keep the spinal fixation device 210 aligned on the bottom flat member 84' of the driving instrument 80'.

Referring to FIG. 17A, an alignment rod 70 comprising a cylindrical shaft 72 having a smooth exterior surface 73 and a threaded end 74 may, be threadably attached to the threaded aperture 45 of the spinal fusion implant 40 is shown. The alignment rod 70 fits through the central opening 18 of the spinal fixation device 10 and is used to properly

align the projections 16 and 17 on each side of the spinal fusion implant 40 prior to engaging the vertebrae V. Further, the alignment rod 70 also serves as a guide post for the drilling template instrument 50 described in greater detail below.

Referring to FIG. 17B, as an alternative embodiment of the alignment rod 70, a splined alignment rod 70' that has a finely splined surface 72' along its longitudinal axis and a threaded end 74' that may be attached to the threaded aperture 45 of the spinal fusion implant is shown.

Referring to FIG. 18, a drilling template instrument 50 for creating a pair of insertion holes 53a and 53b in each of the vertebrae V for receiving each of the projection 16 and 17 respectively is shown. The drilling template instrument 50 has a template 52 with a central aperture 54 therethrough and guide passages 55 and 56 for guiding a drill bit 51 of a drilling tool. Attached to the template 52 is a handle 58 which angles away from the template 52 so as not to obstruct the line of sight of the surgeon and to allow easy access to the template 52 and easy access to the guide holes 55 and 56 for the drill bit 51. Extending from the center of the bottom surface of the template 52 is a central member 59 (similar in structure and function to the central bar 35) for mating to an already implanted intervertebral spinal fusion implant 40. The central member 59 interdigitates with the depression 42 of the spinal fusion implant 40 so that the template 52 is properly oriented about the spinal fusion implant 40 and the guide holes 55 and 56 are properly oriented with respect to the vertebrae V adjacent to the spinal fusion implant 40. The alignment rod 70 serves as a guide post for the drill template instrument 50 as it fits through the central aperture 54 of the template 52 and aligns the template 52 with respect to the spinal fusion implant 40 and insures that it is coaxial. The central aperture 54 of the drilling template instrument 50 is smooth so that if it is placed over a splined alignment rod 70' the drilling template instrument 50 may be easily rotated about the splined alignment rod 70' into position such that the central member 59 is able to mate and interdigitate with the depression 44 of the spinal fusion implant 40.

Referring to FIGS. 19–24, the spinal fixation device 10 of the present invention is inserted in the following manner: At least one spinal fusion implant 40 is surgically implanted so that it is substantially within the disc space between two adjacent vertebrae V and engages at least a portion, of each of the two adjacent vertebrae V. Once the spinal fusion implant 40 is in place, the alignment rod 70 is attached to the threaded aperture 45 of the spinal fusion implant 40. The alignment rod 70 serves as a guide post for the drilling template instrument 50 as it fits through the central aperture 54 of the template 52 and aligns the template 52 coaxially with respect to the spinal fusion implant 40.

Referring to FIG. 20, once the template 52 is properly aligned and the drilling template instrument 50 is seated so that the central member 59 interdigitates with the spinal fusion implant 40, the insertion holes 53a and 53b are drilled in each of the adjacent vertebrae V with a drilling instrument having a drill bit, 51 with a diameter that is substantially smaller than the diameter of each the projections 16 and 17 of the staple member 12.

Once the drilling of the insertion holes 53a and 53b is completed, the drill template instrument 50 is removed from the spinal fusion implant 40 and from the alignment rod 70. The alignment rod 70 is left in place attached to the threaded aperture 45 of the spinal fusion implant 40.

Referring to FIG. 21, the staple member 12 is placed onto the driving instrument 80 used for driving and fixing the

staple member 12 into the vertebrae V so that the bottom flat member 84 and the posts 88a and 88b are aligned with the top member 14 and the depressions 22a and 22b of the top member 14. The alignment rod 70 serves as a guide post for the staple member 12 as it fits through the central opening 18 of the staple member 12 and aligns the staple member 12 coaxially with respect to the spinal fusion implant 40.

Referring to FIG. 22, once the staple member 12 is properly placed onto the bottom flat member 84 of the driving instrument 80, the staple member 12 and the driving instrument 80 are aligned with respect to the alignment rod 70 so that the alignment rod 70 passes through the central opening 18 of the staple member 12 and is inserted into the central hollow portion 89 of the driving instrument 80. The staple member 12 and the driving instrument 80 are then lowered along the alignment rod 70 so that the sharp distal end 32 of each of the projections 16 and 17 comes into contact with the external surface of the vertebrae V and is aligned with the previously drilled insertion holes 53a and 53b.

As shown in FIG. 22A, it is preferred that the insertion holes 53a and 53b be drilled so that when the projections 16 and 17 are inserted into the holes 53a and 53b, the incline planed inner surface 36 of each of the projections 16 and 17 contacts the inner wall W of the insertion holes 53a and 53b that is closest to the spinal fusion implant 40. In this manner a compression force F is created as each of the projections 16 and 17 of the staple member 12 is inserted into insertion holes 53a and 53b, respectively, compressing the bone of the vertebrae V toward the spinal fusion implant 40.

Referring to FIG. 23, the staple member 12 is then driven into the vertebrae V by applying a high impaction force to the driving instrument 80 with a hammer H or other impacting means against the top flat member 86 of the driving instrument 80. The staple member 12 is driven into the vertebrae V such that the projections 16 and 17 are moved forward into the insertion holes 53a and 53b, respectively, until the bottom surface 30 of the top member 14 of the staple member 12 comes to rest against the surface of the vertebrae V.

Referring to FIGS. 23–24, the driving instrument 80 is lifted away from the alignment rod 70 so that the alignment rod 70 is no longer within the central hollow portion 89 of the driving instrument 80. The alignment rod 70 is unthreaded from the threaded aperture 45 and is removed from the spinal fusion implant 40. The staple member 12 is secured to the spinal fusion implant 40 with the locking screw 60 which has a threaded pattern 62 with a reduced pitch. The reduced pitch of the locking screw 60 locks the locking screw 60 to the spinal fusion implant 40 with minimal turning of the locking screw 60 and prevents any unwanted loosening. Further, any of the three embodiments of the securing means 65, 65a or 65b described above in reference to FIGS. 15A–15C may be used to further prevent any unwanted loosening and outward excursion of the screw 60.

Referring, back to FIG. 12, once the staple member 12 is driven into the vertebrae V and is secured to the spinal fusion implant 40, the spinal fusion implant 40 is prevented from rotating along its rotational axis R by its connection to the staple member 12 which is fixated across the disc space between the vertebrae V. The staple member 12 is prevented from backing out from the vertebrae V along the longitudinal axis L by its connection to the spinal fusion implant 40 and by the segmented and ratcheted portions 34 of the projections 16 and 17. In this manner, the staple member 12

and the spinal fusion implant **40** interact to prevent the dislodgement of each other from the vertebrae **V** in which they are implanted. Thus, the staple member **12** is made safe against dislodgement by attachment to the spinal fusion implant **40** and the stability of the spinal fusion implant **40** is assured as it is also stabilized by the staple member **12** and each works in connection with the other to remove the only remaining degree of freedom that would allow for the disengagement of either. In addition, the inclined inner surface **36** at the distal end **32** of the projections **16** and **17** forces bone toward the spinal fusion implant **40** along force lines **F** to further secure the spinal fusion implant **40** and further prevent the dislodgement of the spinal fusion implant **40**.

It is appreciated by those skilled in the art that when the bone of the vertebrae **V** is sufficiently soft, a shorter method (hereinafter referred to as the "Short Method") of inserting the spinal fixation device **10** is possible by omitting the steps of drilling the insertion holes **53a** and **53b** prior to inserting the staple member **12** into the vertebrae **V**.

Referring to FIG. **22**, in the Short Method, the splined alignment rod **70'** that is finely splined along its longitudinal axis is used instead of the alignment rod **70**. Once the splined alignment rod **70'** has been attached to the spinal fusion implant **40**, the staple member **12** may be placed over the splined alignment rod **70'** so that the splined alignment rod **70'** passes through the aperture **18** and into the central aperture **89** of the driving instrument **80**. The central aperture **89** of the driving instrument **80** is correspondingly splined to the splines of the splined alignment rod **70'** so that the staple member **12** can be aligned with respect to the spinal implant **40**. The alignment of the staple member **12** and the driving instrument **80** is maintained as the corresponding splines of the central aperture **89** interdigitate with the splines of the splined alignment rod **70'** and prevent the rotation of the staple member **12** about the splined alignment rod **70'**. The prevention of rotation about the splined alignment rod **70'** is especially important when the Short Method is used to insert the spinal fixation device **10**, as no insertion holes **53a** and **53b** have been drilled in the vertebrae **V**. The staple **12** can be driven directly into the vertebrae **V** by the application of a high impact force to the driving instrument **80** as described above and shown in FIG. **22**.

Once the staple member **12** is driven into the vertebrae **V**, the steps of the longer method described above are used to secure the spinal fixation device to the spinal fusion implant **40** are the same. The Short Method of inserting the staple member **12** reduces the amount of time, required to insert and secure the spinal fixation device **10** of the present invention and thus reduces the overall duration of the spinal fixation surgical procedure.

While the present invention has been described with respect to its preferred embodiment, it is recognized that alternative embodiments of the present invention may be devised without departing from the inventive concept.

For example, referring to FIG. **25**, a first alternative embodiment of a spinal fixation device **110** having a staple member **112** with a top member **114** generally in the shape of an elongated oval having two curved sides **124a** and **124b** is shown. In this alternative embodiment, the curved sides **124a** and **124b** have upper edges **125a** and **125b**, respectively, that are radiused to conform to the external curvature of the vertebrae **V** thereby creating smooth contoured surfaces as described above for the spinal fixation device **10**, the preferred embodiment of the present invention. The top member **114** has openings **122a** and **122b** in the

upper surface **120** of the top member **114** and has two projections **116** and **117** depending downwardly from the bottom surface **130** of the top member **114** at opposite ends of the staple member **112**. The projections **116** and **117** are the same as the projections **16** described above for the preferred embodiment.

Referring to FIG. **26**, a second alternative embodiment of the spinal fixation device **210** having a staple member **212** is shown with a top member **214** that is generally rectangular in shape and has an upper surface **220** with openings **222a**, **222b**, **222c**, and **222d**. The top member **214** has four projections **216**, **217**, **218**, and **219** depending from its bottom surface at each of its corners. The projections **216–217** are the same as the projections **16** and **17** described above in the preferred embodiment. The top member **214** has four straight sides **228a**, **228b**, **228c**, and **228d** having upper edges **225a**, **225b**, **225c**, and **225d**, respectively, that are radiused to conform to the external curvature of the vertebrae **V** create a smooth surface as described above for the preferred embodiment. The driving instrument **80'** shown in FIG. **16B** is used to insert the spinal fixation device **210**.

Referring to FIG. **27**, a third alternative embodiment of the spinal fixation device **310** having a staple **312** with a top member **314** that is generally triangular is shown. The top member **314** has two projections **316** and **317** depending from the bottom surface of the top member **314** that engage the vertebrae **V**. Extending from the center of the bottom surface of the top member **314** is a central member **390** which is similar to the central bar **35** of the preferred embodiment of the spinal fixation device **10** in that the central member **390** interdigitates with the depression **44a** of the spinal fusion implant **40**. However, the central bar **390** also has an extension arm **392** that extends laterally from the top member **314** to span the diameter of an adjacent spinal fusion implant **41**. The extension arm **392** interdigitates with the depression **44** of the spinal implant **41**. The extension arm **392** has a central aperture **374** for receiving a screw **60b** used to couple the extension arm **392** to the spinal fusion implant **41**. In this manner, a single spinal fixation device **310** is capable of interdigitate with two adjacent spinal fusion implants **40** and **41** to lock and prevent the rotation and any excursion of the spinal fusion implants **40** and **41**. The fixation of two spinal fusion implants **40** and **41** is possible while leaving no protruding metal, such as the top member **314**, on the side of the spine where the vessels are located in close approximation to the vertebrae as is the case with the  $L_4$  and  $L_5$  vertebrae where the vessels are located over the left side of those vertebrae. It is appreciated that any of the securing means **65–65b**, described above may be used to lock the screw **60b** to the extension arm **392**.

Referring to FIG. **28**, a fourth alternative embodiment of the spinal fixation device **410** having a staple member **412** with a top member **414** that is generally triangular in shape is shown in the installed position. The top member **414** is wider and larger than top member **14** as it is used with an implant **440** having a large diameter in the range of 22.0 mm to 28.0 mm. The top member **414** needs to be wider when used with implant **440** in order to provide a central bar **435** of sufficient length to interdigitate and mate with the depression **444** of the implant **440** in order to prevent its rotation. Further, the top member **414** is tapered at portion **416** so as not to cause erosion or pressure against the vessels that may be present in the area of the spine adjacent to the portion **416** of the top member **414**.

Referring to FIGS. **29–32**, a fifth alternative embodiment of the spinal fixation device **510** with a staple member **512**

having a generally rectangular top member **514** is shown. The staple member **512** is similar in structure to the staple **212** described above except that the top member **514** has multipronged projection blades **516** and **517** depending from its lower surface **530** as shown in FIG. **30**. The multipronged projection blades **516** and **517** have the same function and similar structure as the projections **16** and **17** described above and include segmented and ratcheted portions **534** which are similar in design are function to segmented and ratcheted portions **34**. The multipronged blade projections **516** and **517** offer the added advantage of increasing the strength and stability of the staple member **514** once it is inserted into the bone of the vertebrae V providing a greater area of engagement of the staple member **512** to the vertebrae V.

The lower surface **530** has knobs **532** and **534** extending therefrom for engaging and interdigitating with a spinal implant **540** having an insertion end **541** with openings **542** and **544** for receiving knobs **532** and **534** respectively.

Referring to FIGS. **31** and **32**, the spinal fusion implant **540** is shown inserted within the disc space between two adjacent vertebrae V. The spinal implant **540** is generally rectangular in shape. The multiprong blade projections **516** and **517** have a width that is approximately equal or slightly less than the width of the spinal fusion implant **540**. Once inserted, the spinal fixation device, **510** compresses the bone of the vertebrae V towards the spinal fusion implant **540** as discussed above in reference to FIG. **12**. The spinal fixation device **510** may be secured to the spinal fusion implant **540** with a screw **60** as discussed above.

The spinal fixation device **510** having a staple member **512** is the preferred embodiment of the present invention for use with a multi-segmental spinal alignment means **600** described in greater detail below in that the staple **512** provides a more solid anchoring means that can resist greater torsion forces resulting from the application of the multi-segmental spinal alignment means **600** to align the spine.

Alternatively, for all of the embodiments described above, the spinal fixation device **10** of the present invention could be made of resorbable materials, such as bio-compatible resorbable plastics, that resorb at an appropriate rate such that once the spinal fixation device **10** is no longer needed (i.e. when spinal fusion is complete) the body would resorb the spinal fixation device **10**. One such resorbable material is polygalactone, however any other resorbable plastic or other material safely usable within the human body are also within the scope of the present invention.

Further, the spinal fixation device could be only in part resorbable such that the projections **16** and **17** of the staple member **12** would be non-resorbable and would remain incarcerated in the vertebrae V and sealed off once the resorbable portion of the staple is resorbed by the body.

Referring to FIGS. **33** and **34**, as a further application, the spinal fixation device **510** of the present invention may be used as an anchor for a multi-segmental spinal alignment means **600**, such that a multiplicity of spinal fixation devices may then be interconnected via a cable, rod, bar, or plate, so as to achieve or maintain any desired multi-segment spinal alignment. In the preferred embodiment, the multi-segmental spinal alignment means **600** comprises more than one spinal fixation device **510** of the present invention placed in series along the spine such that each spinal fixation device **510** spans one disc D and engages two adjacent vertebrae V. The spinal fixation device **510** is preferred over the other embodiments of the present invention in that it has

a greater area of engagement with the vertebrae V so as to provide a solid anchoring means for the multi-segmental spinal alignment means **600**. However, it is appreciated that other embodiments including but not limited to those described herein may be utilized as anchoring means for the multi-segmental spinal alignment means **600**.

When used as an anchor, each spinal fixation device **510** interdigitates with and is connected to a spinal fusion implant **610** having an insertion end **612**, an interior chamber **614** and is inserted in the disc space between the two adjacent vertebrae. The spinal fusion implant **610** has a threaded blind hole **620** for receiving a threaded post **622** therein. The blind hole **620** has a casing that is made of strong surgically, implantable material such as, but not limited to titanium. The casing **624** extends from the insertion end **612** of the spinal fusion implant **610** into the interior central chamber **614**. The insertion end **612** has a rigid construction that is capable of withstanding high torsion forces resulting from the tensioning of the multi-segmental spinal alignment means to align segments of the spine. In the preferred embodiment, the insertion end **612** of the spinal fusion implant has an end portion **626** that closes the insertion end **612**. The end portion is substantially thicker than the rest of the spinal fusion implant **610** and in the preferred embodiment, the end portion **626** has thickness ranging from 1.5 mm to 4.0 mm, with 2.5 mm being the preferred thickness.

Referring to FIG. **34**, the threaded post **622** has a threaded end **628** with a locking thread pattern that is substantially longer than the locking thread pattern **62** of the screw **60** described above and a head portion **630** having a hole **632** for receiving a rod **634** or a cable therethrough. The head portion **630** has a rounded exterior surface to prevent any damage such as aortic erosion to the vessels in the area adjacent to the spine. In the preferred embodiment the threaded post has a diameter ranging from 3.0 mm to 6.0 mm, with 4.5 mm being the preferred diameter and has a length ranging from 15.0 mm to 25.0 mm, with 20.0 mm being the preferred length. The head portion **630** extends at a height above the top member **514** of the spinal fixation device **510** of approximately 8.0 mm to 16.0 mm, with 12.0 mm being the height preferred once it is threadably attached the spinal fusion implant **610** such that it does not significantly protrude from the spinal column into the tissue and vessels adjacent thereto.

Once the threaded post **622** is attached to the spinal fusion implant **610**, the head portion **630** of each threaded post **622** are connected to one another by the rod **634** having a sufficient diameter to fit through the hole **632** of each head portion **630**. The rod **634** has at least a portion thereof that is threaded so that a plurality of lock nuts **638** may be used to secure the rod **634** to the head portions **630**. The lock nuts **638** may also be used as length adjusting means to adjust the length of the rod **634** between head portions **630** so that segmental portions of the spine may be held closer together or held further apart for the purposes of aligning the spine. It is appreciated that a plurality of multi-segmental spinal alignment means **600** may be placed in series either on one side or on opposite sides of the spine, such that one side of the spine may be extended while the other side may be held stationary or may be compressed in order to achieve proper spinal alignment. The multi-segment spinal alignment may be maintained by keeping the rod tensioned with the lock nuts **638** or by any other means well known by those skilled in the art. It is also appreciated that in place of a rod **634** a cable, a plate or any other means well known by those skilled in the art may be used to interconnect the multi-segmental spinal alignment means.

Referring to FIG. 35, a sixth alternative embodiment of the spinal fixation device of the present invention is shown and generally referred to by the numeral 710. The spinal fixation device 710 comprises a top member 714 that is similar to the top member 14 described above, except that it does not have projections 16 and 17 extending from the bottom surface. Like numbers are being used to designate identical features of the top members 14 and 714.

In the top member 714, instead of having projections 16 and 17, independent projection members 716 and 717 in the form of screws are used to secure the top member 714 of the spinal fixation device 710 to the vertebrae V of the spine. The projection screw members 716 and 717 each terminate in a sharp distal end 720 and 722 respectively, have a threaded portion 723, and have screw heads 724 and 726 for engaging a screw driver or similar driving instrument.

The top member 714 has a hole 728 on one end and a hole 730 at its other end through which each of the projection screw members 716 and 717 respectively, may pass. The projections screw members 716 and 717 pass through the holes 728 and 730 to engage the vertebrae V. Each of the holes 728 and 730 has a concentric counter sunk recess 732 for receiving and seating the screw heads 724 and 726 of the projection screw members 716 and 717 so that the screw heads 724 and 726 are flush or below the top surface 20 of the top member 714 once inserted into the vertebrae V.

As the projection screw members 716 and 717 are threaded, they can be rotationally advanced into the vertebrae instead of by way of an impaction force such that the potential for damage to the vertebrae V is reduced. The threads of the threaded portion 723 follow one another as the projection screw members 716 and 717 are being screwed into the bone such that the integrity of the vertebrae V is preserved. Also, as the projection screw members 716 and 717 are independent from the top member 714, the penetration depth of the spinal fixation device 710 into the bone of the vertebrae V may be easily altered by selecting different sized projection screw members 716 and 717 appropriate for the particular vertebrae being fused. Further, it is possible to configure the holes 728 and 730 in the top member 714 such that the projection screw members 716 and 717 may be inserted into the vertebrae V from a number of different angles relative to the top member 714.

Adjacent and proximate to each of the holes 728 and 730 are threaded openings 740 and 742, respectively, for receiving locking screws 744 and 746 respectively. Each of the locking screws 744 and 746 have a head portion 750 and a locking thread portion 754 for threadably and lockably engaging the threaded openings 740 and 742. The locking screws 744 and 746 are attached to the top member 714 after the projection screw members 716 and 717 have been inserted into the vertebrae V. At least a part of the head portion 750 and 752 blocks and preferably makes contact with the screw projections 716 and 717 to prevent any unwanted loosening and outward excursion of the screw projections 716 and 717.

It is appreciated that the projection members 716 and 717, instead of being threaded screws, may have a number of other configurations such as, but not limited to, the configurations of the projections described above for the various embodiments of the present invention. If the projection members 716 and 717 are ratcheted instead of being threaded, they can be driven into the vertebrae V with a driving instrument and impaction force as described above for the method of the present invention.

While the present invention has been described with respect to its preferred embodiment and a number of alter-

native embodiments, it is recognized that additional variations of the present invention may be devised without departing from the inventive concept and scope of the present invention.

What is claimed is:

1. An apparatus comprising:

an interbody spinal fusion implant for surgical implantation within a disc space between two adjacent vertebral bodies in a segment of a human spine having a longitudinal axis, said implant comprising upper and lower portions for contacting each of the adjacent vertebral bodies when positioned therein, each of said upper and lower portions having at least one opening adapted to communicate with one of the adjacent vertebral bodies, said openings of said upper and lower portions being in communication with one another and adapted for permitting the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant, a hollow interior for holding bone growth promoting material, said hollow interior being in communication with at least one opening in each of said upper and lower portions, said implant having an insertion end for entry into the spine, a trailing end opposite said insertion end, and a mid-longitudinal axis passing through said implant from said insertion end to said trailing end; and

two opposed bone screws adapted for placement one each into each of the adjacent vertebral bodies adjacent the disc space to be fused and into which said implant is adapted to be positioned, each of said two opposed bone screws having a proximal end, a distal end, and a threaded shaft, said two opposed bone screws being connected to said implant proximate said trailing end so that a substantial and continuous length of said threaded shaft of each of said two opposed bone screws in a direction from said distal end toward said proximal end of said two opposed bone screws are spaced from said implant, respectively, said mid-longitudinal axis and said proximal ends of said two opposed bone screws being in a plane and when in use the plane being aligned with the longitudinal axis of the spine.

2. The apparatus of claim 1, wherein at least a portion of said upper and lower portions are arcuate along at least a portion of their length.

3. The apparatus of claim 1, wherein said upper and lower portions further comprise a protrusion for engaging the adjacent vertebral bodies.

4. The apparatus of claim 3, wherein said protrusion is a thread.

5. The apparatus of claim 1, wherein at least one of said insertion and trailing ends is open for loading bone growth promoting material into said hollow interior.

6. The apparatus of claim 5, further comprising an end cap for closing said open end.

7. The apparatus of claim 1, wherein said hollow interior is a chamber and the bone growth promoting material includes a bone graft.

8. The apparatus of claim 1, wherein said implant is configured for implantation across the disc space in the thoracolumbar region of the human spine.

9. The apparatus of claim 1, wherein said spinal implant includes an artificial material other than bone.

10. The apparatus of claim 1, further comprising a connector having openings for receiving said bone screws, said trailing end having an aperture for receiving said connector.

11. The apparatus of claim 10, wherein said aperture in said trailing end is generally aligned with the longitudinal axis of the spine.

19

12. The apparatus of claim 10, wherein said connector has a shank adapted to cooperatively engage said aperture.
13. The apparatus of claim 12, wherein said shank has threads adapted to cooperatively engage said aperture.
14. The apparatus of claim 1, wherein said bone screws are appropriately sized and configured to function for their intended purpose.
15. The apparatus of claim 1, wherein said bone screws are sized and configured for use in the cervical spine.
16. The apparatus of claim 1, wherein said bone screws are sized and configured for use in the lumbar spine.
17. The apparatus of claim 1, wherein said bone screws are sized and configured for use in the anterior aspect of the spine.
18. The apparatus of claim 1, wherein said bone screws have a sharp distal end.
19. The apparatus of claim 1, wherein said bone screws have a head opposite said distal end for engaging a driving instrument.
20. The apparatus of claim 1, wherein said bone screws having different lengths.
21. The apparatus of claim 1, further comprising a locking device for locking each of said bone screws into connected position with said implant.
22. The apparatus of claim 21, wherein said locking device is a locking screw.
23. The apparatus of claim 1, wherein said implant is made of an artificial material that is stronger than bone.
24. The apparatus of claim 1, wherein said implant is made of an artificial material that is harder than bone received.
25. The apparatus of claim 1, wherein said implant comprises harvested bone.
26. The apparatus of claim 1, wherein said implant is in combination with bone growth promoting material.
27. The apparatus of claim 26, wherein said bone growth promoting material is selected from one of bone morphogenetic protein, hydroxyapatite, and genes coding for the production of bone.
28. The apparatus of claim 1, wherein said implant is treated with a bone growth promoting substance.
29. The apparatus of claim 1, wherein said implant is a source of osteogenesis.
30. The apparatus of claim 1, wherein said implant is at least in part bioabsorbable.
31. The apparatus of claim 1, wherein said implant comprises metal.
32. The apparatus of claim 31, wherein said metal includes titanium.
33. The apparatus of claim 1, wherein said implant comprises a plastic material.
34. The apparatus of claim 1, wherein said implant comprises a ceramic material.
35. The apparatus of claim 1, wherein said implant is formed of a porous material.
36. The apparatus of claim 1, wherein said implant is formed of a material that intrinsically participates in the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant.
37. The apparatus of claim 1, wherein said implant has an interior surface and a hollow defined therein, said hollow being capable of containing bone growth promoting material.
38. The apparatus of claim 31, wherein said bone growth promoting material is selected from one of bone morphogenetic protein, hydroxyapatite, and genes coding for the production of bone.

20

39. The apparatus of claim 1, wherein said at least one opening is adapted to retain fusion-promoting materials.
40. The apparatus of claim 1, wherein at least a portion of said implant is treated to promote bone ingrowth between said implant and said adjacent vertebral bodies.
41. The apparatus of claim 1, wherein said implant is in combination with harvested bone.
42. An apparatus comprising:  
an interbody spinal fusion implant for surgical implantation within a disc space between two adjacent vertebral bodies in a segment of the spine having a longitudinal axis, said implant comprising upper and lower portions for contacting each of the adjacent vertebral bodies when positioned therein, each of said upper and lower portions having at least one opening adapted to communicate with one of the adjacent vertebral bodies, said openings of said upper and lower portions being in communication with one another and adapted for permitting for the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant, a hollow interior for holding bone growth promoting material, said hollow interior being in communication with at least one opening in each of said upper and lower portions, said implant having an insertion end for entry into the spine, a trailing end opposite said insertion end, and a mid-longitudinal axis passing through said implant from said insertion end to said trailing end; and  
two opposed bone screws adapted for placement one each into each of the adjacent vertebral bodies adjacent the disc space to be fused and into which said implant is adapted to be positioned, each of said two opposed bone screws having a proximal end, a distal end, and a threaded shaft, said two opposed bone screws being attached to said implant proximate said trailing end so that a substantial and continuous length of said threaded shaft of each of said two opposed bone screws in a direction from said distal end toward said proximal end of said two opposed bone screws are spaced from said implant, respectively, said mid-longitudinal axis and said proximal ends of said two opposed bone screws being in a plane and when in use the plane being aligned with the longitudinal axis of the spine.
43. The apparatus of claim 42, wherein at least a portion of said upper and lower portions are arcuate along at least a portion of their length.
44. The apparatus of claim 42, wherein said upper and lower portions further comprise a protrusion for engaging the adjacent vertebral bodies.
45. The apparatus of claim 44, wherein said protrusion is a thread.
46. The apparatus of claim 44, wherein at least one of said insertion and trailing ends is open for loading bone growth promoting material into said hollow interior.
47. The apparatus of claim 46, further comprising an end cap for closing said open end.
48. The apparatus of claim 42, wherein said hollow interior is a chamber and the bone growth promoting material includes a bone graft.
49. The apparatus of claim 42, wherein said implant is configured for implantation across the disc space in the thoracolumbar region of the human spine.
50. The apparatus of claim 42, wherein said spinal implant includes an artificial material other than bone.
51. The apparatus of claim 42, further comprising a connector having openings for receiving said bone screws, said trailing end having an aperture for receiving said connector.

52. The apparatus of claim 51, wherein said aperture in said trailing end is generally aligned with the longitudinal axis of the spine.

53. The apparatus of claim 51, wherein said connector has a shank adapted to cooperatively engage said aperture.

54. The apparatus of claim 53, wherein said shank has threads adapted to cooperatively engage said aperture.

55. The apparatus of claim 42, wherein said bone screws are appropriately sized and configured to function for their intended purpose.

56. The apparatus of claim 42, wherein said bone screws are sized and configured for use in the cervical spine.

57. The apparatus of claim 42, wherein said bone screws are sized and configured for use in the lumbar spine.

58. The apparatus of claim 42, wherein said bone screws are sized and configured for use in the anterior aspect of the spine.

59. The apparatus of claim 42, wherein said bone screws have a sharp distal end.

60. The apparatus of claim 42, wherein said bone screws have a head opposite said distal end for engaging a driving instrument.

61. The apparatus of claim 42, wherein said bone screws having different lengths.

62. The apparatus of claim 42, further comprising a locking device for locking each of said bone screws into connected position with said implant.

63. The apparatus of claim 62, wherein said locking device is a locking screw.

64. The apparatus of claim 31, wherein said implant is made of an artificial material that is stronger than bone.

65. The apparatus of claim 31, wherein said implant is made of an artificial material that is harder than bone.

66. The apparatus of claim 31, wherein said implant comprises harvested bone.

67. The apparatus of claim 31, wherein said implant is in combination with bone growth promoting material.

68. The apparatus of claim 67, wherein said bone growth promoting material is selected from one of bone morphogenetic protein, hydroxyapatite, and genes coding for the production of bone.

69. The apparatus of claim 42, wherein said implant is treated with a bone growth promoting substance.

70. The apparatus of claim 42, wherein said implant is a source of osteogenesis.

71. The apparatus of claim 42, wherein said implant is at least in part bioabsorbable.

72. The apparatus of claim 42, wherein said implant comprises metal.

73. The apparatus of claim 72, wherein said metal includes titanium.

74. The apparatus of claim 42, wherein said implant comprises a plastic material.

75. The apparatus of claim 42, wherein said implant comprises a ceramic material.

76. The apparatus of claim 42, wherein said implant is formed of a porous material.

77. The apparatus of claim 42, wherein said implant is formed of a material that intrinsically participates in the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant.

78. The apparatus of claim 42, wherein said implant has an interior surface and a hollow defined therein, said hollow being capable of containing bone growth promoting material.

79. The apparatus of claim 78, wherein said bone growth promoting material is selected from one of bone morpho-

genetic protein, hydroxyapatite, and genes coding for the production of bone.

80. The apparatus of claim 42, wherein said at least one opening is adapted to retain fusion-promoting materials.

81. The apparatus of claim 42, wherein at least a portion of said implant is treated to promote bone ingrowth between said implant and said adjacent vertebral bodies.

82. The apparatus of claim 42, wherein said implant is in combination with harvested bone.

83. An apparatus comprising:  
an interbody spinal fusion implant for surgical implantation within a disc space between two adjacent vertebral bodies in a segment of the spine having a longitudinal axis, said implant comprising upper and lower portions for contacting each of the adjacent vertebral bodies when positioned therein, each of said upper and lower portions having at least one opening adapted to communicate with one of the adjacent vertebral bodies, said openings of said upper and lower portions being in communication with one another and adapted for permitting for the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant, a hollow interior for holding bone growth promoting material, said hollow interior being in communication with at least one opening in each of said upper and lower portions, said implant having an insertion end for entry into the spine, a trailing end opposite said insertion end, and a mid-longitudinal axis passing through said implant from said insertion end to said trailing end; and

two opposed bone screws adapted for placement one each into each of the adjacent vertebral bodies adjacent the disc space to be fused and into which said implant is adapted to be positioned, each of said two opposed bone screws having a proximal end, a distal end, and a threaded shaft, said two opposed bone screws being coupled to said trailing end of said implant so that a substantial and continuous length of said threaded shaft of each of said two opposed bone screws in a direction from said distal end toward said proximal end of said two opposed bone screws are spaced from said implant, respectively, said mid-longitudinal axis and said proximal ends of said two opposed bone screws being in a plane and when in use the plane being aligned with the longitudinal axis of the spine.

84. The apparatus of claim 83, wherein at least a portion of said upper and lower portions are arcuate along at least a portion of their length.

85. The apparatus of claim 83, wherein said upper and lower portions further comprise a protrusion for engaging the adjacent vertebral bodies.

86. The apparatus of claim 85, wherein said protrusion is a thread.

87. The apparatus of claim 83, wherein at least one of said insertion and trailing ends is open for loading bone growth promoting material into said hollow interior.

88. The apparatus of claim 87, further comprising an end cap for closing said open end.

89. The apparatus of claim 83, wherein said hollow interior is a chamber and the bone growth promoting material includes a bone graft.

90. The apparatus of claim 83, wherein said implant is configured for implantation across the disc space in the thoracolumbar region of the human spine.

91. The apparatus of claim 83, wherein said spinal implant includes an artificial material other than bone.

92. The apparatus of claim 83, further comprising a connector having openings for receiving said bone screws, said trailing end having an aperture for receiving said connector.

93. The apparatus of claim 92, wherein said aperture in said trailing end is generally aligned with the longitudinal axis of the spine.

94. The apparatus of claim 92, wherein said connector has a shank adapted to cooperatively engage said aperture.

95. The apparatus of claim 94, wherein said shank has threads adapted to cooperatively engage said aperture.

96. The apparatus of claim 83, wherein said bone screws are appropriately sized and configured to function for their intended purpose.

97. The apparatus of claim 83, wherein said bone screws are sized and configured for use in the cervical spine.

98. The apparatus of claim 83, wherein said bone screws are sized and configured for use in the lumbar spine.

99. The apparatus of claim 83, wherein said bone screws are sized and configured for use in the anterior aspect of the spine.

100. The apparatus of claim 83, wherein said bone screws have a sharp distal end.

101. The apparatus of claim 83, wherein said bone screws have a head opposite said distal end for engaging a driving instrument.

102. The apparatus of claim 83, wherein said bone screws having different lengths.

103. The apparatus of claim 83, further comprising a locking device for locking each of said bone screws into connected position with said implant.

104. The apparatus of claim 103, wherein said locking device is a locking screw.

105. The apparatus of claim 83, wherein said implant is made of an artificial material that is stronger than bone.

106. The apparatus of claim 83, wherein said implant is made of an artificial material that is harder than bone.

107. The apparatus of claim 83, wherein said implant comprises harvested bone.

108. The apparatus of claim 83, wherein said implant is in combination with bone growth promoting material.

109. The apparatus of claim 108, wherein said bone growth promoting material is selected from one of bone morphogenetic protein, hydroxyapatite, and genes coding for the production of bone.

110. The apparatus of claim 83, wherein said implant is treated with a bone growth promoting substance.

111. The apparatus of claim 83, wherein said implant is a source of osteogenesis.

112. The apparatus of claim 83, wherein said implant is at least in part bioabsorbable.

113. The apparatus of claim 83, wherein said implant comprises metal.

114. The apparatus of claim 113, wherein said metal includes titanium.

115. The apparatus of claim 83, wherein said implant comprises a plastic material.

116. The apparatus of claim 83, wherein said implant comprises a ceramic material.

117. The apparatus of claim 83, wherein said implant is formed of a porous material.

118. The apparatus of claim 83, wherein said implant is formed of a material that intrinsically participates in the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant.

119. The apparatus of claim 83, wherein said implant has an interior surface and a hollow defined therein, said hollow being capable of containing bone growth promoting material.

120. The apparatus of claim 119, wherein said bone growth promoting material is selected from one of bone

morphogenetic protein, hydroxyapatite, and genes coding for the production of bone.

121. The apparatus of claim 83, wherein said at least one opening is adapted to retain fusion-promoting materials.

122. The apparatus of claim 83, wherein at least a portion of said implant is treated to promote bone ingrowth between said implant and said adjacent vertebral bodies.

123. The apparatus of claim 83, wherein said implant is in combination with harvested bone.

124. An apparatus comprising:  
an interbody spinal fusion implant for surgical implantation within a disc space between two adjacent vertebral bodies in a segment of the spine having a longitudinal axis, said implant comprising upper and lower portions for contacting each of the adjacent vertebral bodies when positioned therein, each of said upper and lower portions having at least one opening adapted to communicate with one of the adjacent vertebral bodies, said openings of said upper and lower portions being in communication with one another and adapted for permitting for the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant, a hollow interior for holding bone growth promoting material, said hollow interior being in communication with at least one opening in each of said upper and lower portions, said implant having an insertion end for entry into the spine, a trailing end opposite said insertion end, and a mid-longitudinal axis passing through said implant from said insertion end to said trailing end; and

two opposed bone screws adapted for placement one each into each of the adjacent vertebral bodies adjacent the disc space to be fused and into which said implant is adapted to be positioned, each of said two opposed bone screws having a proximal end, a distal end, and a threaded shaft, said two opposed bone screws being linked to said trailing end of said implant so that a substantial and continuous length of said threaded shaft of each of said two opposed bone screws in a direction from said distal end toward said proximal end of said two opposed bone screws are spaced from said implant, respectively, said mid-longitudinal axis and said proximal ends of said two opposed bone screws being in a plane and when in use the plane being aligned with the longitudinal axis of the spine.

125. The apparatus of claim 124, wherein at least a portion of said upper and lower portions are arcuate along at least a portion of their length.

126. The apparatus of claim 124, wherein said upper and lower portions further comprise a protrusion for engaging the adjacent vertebral bodies.

127. The apparatus of claim 126, wherein said protrusion is a thread.

128. The apparatus of claim 124, wherein at least one of said insertion and trailing ends is open for loading bone growth promoting material into said hollow interior.

129. The apparatus of claim 128, further comprising an end cap for closing said open end.

130. The apparatus of claim 124, wherein said hollow interior is a chamber and the bone growth promoting material includes a bone graft.

131. The apparatus of claim 124, wherein said implant is configured for implantation across the disc space in the thoracolumbar region of the human spine.

132. The apparatus of claim 124, wherein said spinal implant includes an artificial material other than bone.

133. The apparatus of claim 124, further comprising a connector having openings for receiving said bone screws, said trailing end having an aperture for receiving said connector.



134. The apparatus of claim 133, wherein said aperture in said trailing end is generally aligned with the longitudinal axis of the spine.

135. The apparatus of claim 133, wherein said connector has a shank adapted to cooperatively engage said aperture. 5

136. The apparatus of claim 135, wherein said shank has threads adapted to cooperatively engage said aperture.

137. The apparatus of claim 124, wherein said bone screws are appropriately sized and configured to function for their intended purpose. 10

138. The apparatus of claim 124, wherein said bone screws are sized and configured for use in the cervical spine.

139. The apparatus of claim 124, wherein said bone screws are sized and configured for use in the lumbar spine.

140. The apparatus of claim 124, wherein said bone screws are sized and configured for use in the anterior aspect of the spine. 15

141. The apparatus of claim 124, wherein said bone screws have a sharp distal end.

142. The apparatus of claim 124, wherein said bone screws have a head opposite said distal end for engaging a driving instrument. 20

143. The apparatus of claim 124, wherein said bone screws having different lengths.

144. The apparatus of claim 124, further comprising a locking device for locking each of said bone screws into connected position with said implant. 25

145. The apparatus of claim 144, wherein said locking device is a locking screw.

146. The apparatus of claim 124, wherein said implant is made of an artificial material that is stronger than bone. 30

147. The apparatus of claim 124, wherein said implant is made of an artificial material that is harder than bone.

148. The apparatus of claim 124, wherein said implant comprises harvested bone. 35

149. The apparatus of claim 124, wherein said implant is in combination with bone growth promoting material.

150. The apparatus of claim 149, wherein said bone growth promoting material is selected from one of bone morphogenetic protein, hydroxyapatite, and genes coding for the production of bone. 40

151. The apparatus of claim 124, wherein said implant is treated with a bone growth promoting substance.

152. The apparatus of claim 124, wherein said implant is a source of osteogenesis. 45

153. The apparatus of claim 124, wherein said implant is at least in part bioabsorbable.

154. The apparatus of claim 124, wherein said implant comprises metal.

155. The apparatus of claim 154, wherein said metal includes titanium. 50

156. The apparatus of claim 124, wherein said implant comprises a plastic material.

157. The apparatus of claim 124, wherein said implant comprises a ceramic material. 55

158. The apparatus of claim 124, wherein said implant is formed of a porous material.

159. The apparatus of claim 124, wherein said implant is formed of a material that intrinsically participates in the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant. 60

160. The apparatus of claim 124, wherein said implant has an interior surface and a hollow defined therein, said hollow being capable of containing bone growth promoting material. 65

161. The apparatus of claim 160, wherein said bone growth promoting material is selected from one of bone

morphogenetic protein, hydroxyapatite, and genes coding for the production of bone.

162. The apparatus of claim 124, wherein said at least one opening is adapted to retain fusion-promoting materials.

163. The apparatus of claim 124, wherein at least a portion of said implant is treated to promote bone ingrowth between said implant and said adjacent vertebral bodies.

164. The apparatus of claim 124, wherein said implant is in combination with harvested bone.

165. An apparatus comprising  
 an interbody spinal fusion implant for surgical implantation within a disc space between two adjacent vertebral bodies in a segment of a human spine, said implant comprising upper and lower portions for contacting each of the adjacent vertebral bodies when positioned therein, each of said upper and lower portions having at least one opening adapted to communicate with one of the adjacent vertebral bodies, said openings of said upper and lower portions being in communication with one another and adapted for permitting for the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant, a hollow interior for holding bone growth promoting material, said hollow interior being in communication with at least one opening in each of said upper and lower portions, said implant having an insertion end for entry into the spine and a trailing end opposite said insertion end, said hollow interior containing bone morphogenetic protein for promoting bone growth; and

bone screws adapted for placement into each of the adjacent vertebral bodies adjacent the disc space to be fused and into which said implant is adapted to be positioned, said bone screws being connected to said implant proximate said trailing end.

166. The apparatus of claim 165, wherein at least a portion of said upper and lower portions are arcuate along at least a portion of their length. 35

167. The apparatus of claim 165, wherein said upper and lower portions further comprise a protrusion for engaging the adjacent vertebral bodies.

168. The apparatus of claim 167, wherein said protrusion is a thread. 40

169. The apparatus of claim 165, wherein at least one of said insertion and trailing ends is open for loading bone growth promoting material into said hollow interior.

170. The apparatus of claim 169, further comprising an end cap for closing said open end. 45

171. The apparatus of claim 165, wherein said hollow interior is a chamber and the bone growth promoting material includes a bone graft.

172. The apparatus of claim 165, wherein said implant is configured for implantation across the disc space in the thoracolumbar region of the human spine.

173. The apparatus of claim 165, wherein said spinal implant includes an artificial material other than bone.

174. The apparatus of claim 165, further comprising a connector having openings for receiving said bone screws, said trailing end having an aperture for receiving said connector. 55

175. The apparatus of claim 174, wherein said aperture in said trailing end is generally aligned with the longitudinal axis of the spine.

176. The apparatus of claim 174, wherein said connector has a shank adapted to cooperatively engage said aperture.

177. The apparatus of claim 176, wherein said shank has threads adapted to cooperatively engage said aperture.

178. The apparatus of claim 165, wherein said bone screws are appropriately sized and configured to function for their intended purpose. 65

27

179. The apparatus of claim 165, wherein said bone screws are sized and configured for use in the cervical spine.

180. The apparatus of claim 165, wherein said bone screws are sized and configured for use in the lumbar spine.

181. The apparatus of claim 165, wherein said bone screws are sized and configured for use in the anterior aspect of the spine.

182. The apparatus of claim 165, wherein said bone screws have a sharp distal end.

183. The apparatus of claim 165, wherein said bone screws have a head opposite said distal end for engaging a driving instrument.

184. The apparatus of claim 165, wherein said bone screws having different lengths.

185. The apparatus of claim 165, further comprising a locking device for locking each of said bone screws into connected position with said implant.

186. The apparatus of claim 165, wherein said locking device is a locking screw.

187. The apparatus of claim 165, wherein said implant is made of an artificial material that is stronger than bone.

188. The apparatus of claim 165, wherein said implant is made of an artificial material that is harder than bone.

189. The apparatus of claim 165, wherein said implant comprises harvested bone.

190. The apparatus of claim 165, wherein said implant is in combination with bone growth promoting material.

191. The apparatus of claim 190, wherein said bone growth promoting material is selected from one of hydroxyapatite and genes coding for the production of bone.

28

192. The apparatus of claim 165, wherein said implant is treated with a bone growth promoting substance.

193. The apparatus of claim 165, wherein said implant is a source of osteogenesis.

194. The apparatus of claim 165, wherein said implant is at least in part bioabsorbable.

195. The apparatus of claim 165, wherein said implant comprises metal.

196. The apparatus of claim 195, wherein said metal includes titanium.

197. The apparatus of claim 165, wherein said implant comprises a plastic material.

198. The apparatus of claim 165, wherein said implant comprises a ceramic material.

199. The apparatus of claim 165, wherein said implant is formed of a porous material.

200. The apparatus of claim 165, wherein said implant is formed of a material that intrinsically participates in the growth of bone from adjacent vertebral body to adjacent vertebral body through said implant.

201. The apparatus of claim 165, wherein said at least one opening is adapted to retain fusion-promoting materials.

202. The apparatus of claim 165, wherein at least a portion of said implant is treated to promote bone ingrowth between said implant and said adjacent vertebral bodies.

203. The apparatus of claim 165, wherein said implant is in combination with harvested bone.

\* \* \* \* \*

UNITED STATES PATENT AND TRADEMARK OFFICE  
**CERTIFICATE OF CORRECTION**

PATENT NO. : 6,364,880 B1  
DATED : April 2, 2002  
INVENTOR(S) : Gary Karlin Michelson

Page 1 of 1

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

Title page.

Item [56], **References Cited**, U.S. PATENT DOCUMENTS, please insert the following:

-- 5,306,307 \* 4/1994 Senter et al.  
5,702,449 \* 12/1997 McKay --.

Item [57], **ABSTRACT**,

Line 17, change "as,the" to -- as the --.

Column 19,

Line 31, delete "received";  
Line 37, change "morphage-" to -- morphoge- --; and  
Line 64, change "31" to -- 37 --;

Column 20,

Line 51, change "44" to -- 42 --;

Column 21,

Line 30, change "31" to -- 42 --;  
Line 32, change "31" to -- 42 --;  
Line 34, change "31" to -- 42 --;  
Line 36, change "31" to -- 42 --;

Column 26,

Line 4, delete "a";  
Line 9, after "comprising" insert -- : --; and  
Line 18, change "165" to -- 185 --.

Signed and Sealed this

Twenty-ninth Day of October, 2002

Attest:



Attesting Officer

JAMES E. ROGAN  
Director of the United States Patent and Trademark Office