

US007993347B1

# (12) United States Patent

#### Michelson

# (10) Patent No.: US 7,993,347 B1 (45) Date of Patent: Aug. 9, 2011

# (54) GUARD FOR USE IN PERFORMING HUMAN INTERBODY SPINAL SURGERY

(75) Inventor: Gary Karlin Michelson, Venice, CA

(US)

(73) Assignee: Warsaw Orthopedic, Inc., Warsaw, IN

(US)

(\*) Notice: Subject to any disclaimer, the term of this

patent is extended or adjusted under 35

U.S.C. 154(b) by 880 days.

- (21) Appl. No.: 09/626,636
- (22) Filed: Jul. 27, 2000

#### Related U.S. Application Data

- (60) Continuation of application No. 08/484,927, filed on Jun. 7, 1995, now Pat. No. 6,096,038, which is a division of application No. 08/074,781, filed on Jun. 10, 1993, now Pat. No. 5,484,437.
- (51) **Int. Cl. A61F 5/04** (2006.01)
- (52) U.S. Cl. ...... 606/90

See application file for complete search history.

#### (56) References Cited

#### U.S. PATENT DOCUMENTS

62,123	Α		2/1867	Elgin	
350,420	Α		10/1886	Dillon	
493,730	Α	*	3/1893	MacKenzie	606/170
640,108	Α	*	12/1899	Dalzell	606/170
1,137,585	Α		4/1915	Craig	
2,065,659	Α		12/1936	Cullen	

2,181,746 A	11/1939	Siebrandt			
2,243,718 A	5/1941	De G. Moreira			
2,372,622 A	3/1945	Fassio			
2,514,665 A	7/1950	Myller			
2,537,070 A	1/1951	Longfellow			
2,543,780 A	3/1951	Hipps et al.			
2,677,369 A	5/1954	Knowles			
2,774,350 A	12/1956	Cleveland			
2,789,558 A	4/1957	Rush			
2,832,343 A	4/1958	Mose			
2,842,131 A	7/1958	Smith			
2,878,809 A	3/1959	Treace			
3,128,768 A	4/1964	Geistauts			
3,298,372 A	1/1967	Feinberg			
3,426,364 A	2/1969	Lumb			
3,486,505 A *	12/1969	Morrison 606/90			
3,604,487 A	9/1971	Gilbert			
(Continued)					

#### FOREIGN PATENT DOCUMENTS

DE 1961531 7/1970 (Continued)

#### OTHER PUBLICATIONS

Otero-Vich, Jose M.; Anterior Cervical Interbody Fusion with Threaded Cylindrical Bone; J. Neurosurg 63:750-753 (Nov. 1985).


(Continued)

Primary Examiner — Ryan Severson (74) Attorney, Agent, or Firm — Martin & Ferraro, LLP

#### (57) ABSTRACT

Apparatus and a method of inserting spinal implants is disclosed in which an intervertebral space is first distracted, a hollow sleeve having teeth at one end is then driven into the vertebrae adjacent that disc space. A drill is then passed through the hollow sleeve removing disc and bone in preparation for receiving the spinal implant which is then inserted through the sleeve.

#### 28 Claims, 24 Drawing Sheets

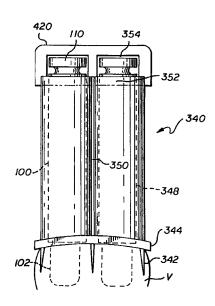

## US 7,993,347 B1

Page 2

LIC DATENT	DOCLIMENTS	4,552,200 A	11/1095	Sinha et al.
U.S. PATENT	DOCUMENTS		11/1985	
3,605,123 A 9/1971				Kapp et al.
3,618,611 A 11/1971				Pratt et al.
	Ames et al.	4,566,466 A	1/1986	Ripple et al.
	Halloran	4,570,623 A	2/1986	Ellison et al.
	Merig, Jr.	4,570,624 A	2/1986	Wu
	Sherwin	4,592,346 A		Jurgutis
	Ma et al.	4,599,086 A	7/1986	
3,855,638 A 12/1974		4,600,000 A		Edwards
	Stubstad et al.	4,602,638 A	7/1986	
3,867,932 A 2/1975		4,604,995 A		Stephens
-,,	Fischell	4,608,052 A		Van Kampen et al.
3,875,595 A * 4/1975	Froning 623/17.12	4,611,581 A	9/1986	
	Fischell		10/1986	
	Neufeld	, ,	12/1986	
3,905,047 A 9/1975	Long	4,634,720 A 4,636,217 A		Dorman et al. Ogilvie et al.
3,915,151 A 10/1975	Kraus	4,636,526 A		Dorman et al.
3,916,907 A 11/1975	Peterson	4,645,503 A		Lin et al.
3,918,440 A 11/1975	Kraus	4,653,486 A	3/1987	
3,942,535 A 3/1976	Schulman	4,653,509 A		Oloff et al.
	Zaffaroni	4,655,777 A	4/1987	
	Bokros et al.	4,661,536 A		Dorman et al.
	Scharbach et al.	4,664,567 A		Edwards
	Kawahara et al.	4,665,920 A		Campbell
	Sawyer et al.	4,677,883 A	7/1987	
D245,259 S 8/1977		4,677,972 A	7/1987	
4,051,905 A 10/1977		4.693.721 A		Ducheyne
	Jumashev et al.	4,696,290 A	9/1987	
	Entherly et al.		10/1987	Dorman et al.
	Mann et al.	4,710,075 A	12/1987	Davison
	Kawahara et al.	4,711,609 A	12/1987	Seefluth
	Berner et al. Stravropoulos et al.	4,713,004 A	12/1987	Linkow et al.
	Broemer et al.		12/1987	
4,168,326 A 9/1979 4,175,555 A 11/1979		4,721,103 A		Freedland
	Grell et al.	4,733,995 A	3/1988	
	Holmes	4,736,738 A		Lipovsek et al.
	Schulman et al.	4,232,679 A		Schulman
4,206,516 A 6/1980		4,743,256 A		Brantigan
	Tomonaga et al.	4,743,260 A	5/1988	
D257,511 S 11/1980		4,759,766 A		Buettner-Janz et al.
	Schulman	4,759,769 A		Hedman et al.
	Jones et al.	4,769,881 A		Pedigo et al.
	Sutherland			Kees, Jr. et al.
	Hirabayashi et al.		11/1988	
4,262,369 A 4/1981			12/1988	
	Evans et al.	4,805,602 A 4,820,305 A		Puno et al.
D260,525 S 9/1981	Lassiter	4,830,000 A	5/1989	Harms et al.
4,289,123 A 9/1981	Dunn	4,834,757 A		Brantigan
4,293,962 A 10/1981	Fuson	4,848,327 A		Perdue
4,309,777 A 1/1982	Patil	4,851,008 A		Johnson
4,328,593 A 5/1982	Sutter et al.	4,863,476 A		Shepperd
4,333,469 A 6/1982	Jeffcoat et al.	4,863,477 A	9/1989	Monson
	Perrett	4,865,603 A	9/1989	
4,349,921 A 9/1982			10/1989	
	Guillemin et al.		10/1989	
	Reimels et al.			Brantigan 623/17.11
	Rezaian	4,903,882 A	2/1990	Long
	Cosentino	4,904,260 A	2/1990	Ray et al.
	Carson 600/138	4,904,261 A	2/1990	Dove et al.
	Hirshorn et al.	4,911,718 A	3/1990	Lee et al.
	Otte et al.	4,913,143 A *	4/1990	Oloff et al 606/170
	Pratt et al.	4,913,144 A	4/1990	Del Medico
4,439,152 A 3/1984 4,450,834 A 5/1984	Fischer	4,936,848 A	6/1990	
	Sutter et al.	4,943,291 A	7/1990	
	Belykh et al.	4,955,885 A	9/1990	
	Nicholson et al.	4,955,908 A		Frey et al.
4,501,269 A 2/1985		4,957,495 A	9/1990	2
	Kambara et al.			Goble et al.
RE31,865 E 4/1985				Ray et al.
	Duarte			Hergenroeder
	Anderson et al.	4,969,888 A	11/1990	Scholten et al.
	Ashman et al.	4,987,904 A	1/1991	
	Rowe, Jr. et al.	4,997,434 A	3/1991	Seedhom et al.
	Jacobson	5,015,247 A	5/1991	Michelson
	Ashman et al.	5,015,255 A	5/1991	Kuslich
	Brighton et al.	5,026,373 A		Ray et al.
•	~			•

D318.629 S						
Seq. 22.573   A 7,1991   Chow   FP   0,162,095   11/1985	D318,629 S	7/1991	Michelson	EP	0 077 159	4/1983
South   Sout				EP		
South   Sout						
Sol.						
Soin						
Sol.						
Soil   Association   Company   Com						
Soft 436 A   12 991 Jang						
Solid						
5084,050   A   11992   Draenert   FP   0421485   A1   101998   S1,010340   A   41992   Wollschildunger et al.   FR   238778   5,1983   S1,010340   A   41992   Wollschildunger et al.   FR   2,205   279   7,1976   S1,010340   A   41992   Wollschildunger et al.   FR   2,205   279   7,1976   S1,010340   A   41992   Great at al.   FR   2,205   S1,010340   A   5,1992   Carboral et al.   FR   2,703   S80   101,1994   A   5,1922   Carboral et al.   FR   2,703   S80   101,1994   A   5,1922   Carboral et al.   FR   2,703   S80   101,1994   A   5,1922   Carboral et al.   FR   2,703   S80   101,1994   A   5,1922   Carboral et al.   FR   2,703   S80   101,1994   A   5,192,192   A   4,1925   A   5,192,192   A   4,1925   A   5,192,192   A   4,1925	5,071,406 A				0577179 A1	1/1994
5,102,414   A   41992   Kirsch   FR   239378   51,1985   5,108,1995   A   41992   Lauran   606.86 B   FR   0,179 0.95   41,1986   5,108,429   A   41992   Center et al.   FR   2,195 729   71,1976   5,108,429   A   41992   Center et al.   FR   2,195 729   71,1976   5,116,304   A   5,1992   Carbwell   GB   1,049,090   11,1977   5,123,026   A   61,1992   Fisharodi   GB   1,492,990   11,1977   5,116,304   A   9,1993   Boyd   GB   1,513,187   11,1978   5,116,304   A   9,1993   Sind principle at al.   606.172   GB   2,006,675   A   2,1981   5,103,237   A   2,1993   Wilk   Sind principle at al.   606.172   GB   2,006,675   A   2,1982   5,193,577   A   3,1993   Bantigm   GB   2,104,277   A   3,1986   5,197,571   A   3,1993   Bantigm   GB   2,104,277   A   3,1986   5,197,571   A   3,1993   Bantigm   GB   2,104,277   A   3,1986   5,197,571   A   3,1993   Santigm   FR   2,1982   5,246,478   A   9,1993   Grabam   FR   60-3,1766   2,1985   5,246,478   A   9,1993   Grabam   FR   60-4,3984   3,1985   5,246,478   A   1,1995   Wilk   Sind principle at al.   Sind	5,071,437 A	12/1991	Steffee		0627204 A2	12/1994
5,103,414   A 41992   Kirsch   ES 283078   51985   5,103,419 A 41992   Laurain   606.86 B   FR 0 179.05   41986   5,103,329 A 0 41992   Carlovell   GB 179.05   41986   5,103,420 A 14992   Carlovell   GB 199.05   41986   5,122,130 A 61992   Keller   GB 1291470   101972   5,123,026 A 61992   Pisharodi   GB 1492990   11/1977   5,116,304 A 21993   Boyd   GB 1492990   11/1977   5,171,278 A 121992   Pisharodi   GB 1492990   11/1977   5,171,278 A 121993   Boyd   GB 1531487   11/1978   5,195,771 A 31993   Bantigan   GB 2076657 A 121/981   5,195,771 A 31993   Bantigan   GB 2164277 A 31/986   5,195,771 A 31993   Bantigan   GB 2164277 A 31/986   5,195,771 A 31993   Bantigan   GB 2164277 A 31/986   5,195,771 A 31993   Salie et al.   P 60-31706   21/985   5,242,443 A 91993   Kambin   P 60-31706   21/985   5,246,458 A 91993   Graham   P 60-13784   21/982   5,246,458 A 11/1993   Wilk   SI 106397   12/1983   5,246,577 A 12/1981   Banmann et al.   SI 1107854   41/986   5,229,225 A 11/994   Bammann et al.   SI 1107854   41/994   5,229,225 A 11/994   Bammann et al.   SI 1107854   41/994   5,243,257 A 1 61985   Shapiro   WO 93001771   21/1993   5,330,310 A 1 11/995   Kanish et al.   Si 1107854   41/996   5343,258 A 61/994   Gobie et al.   WO 93001771   21/993   5,330,317 A 1 11/995   Kanish et al.   Si 1106397   11/996   5,343,538 A 11/995   Kanish et al.   Si 1106397   11/996   5,343,538 A 11/995   Kanish et al.   Si 1106397   11/996   5,343,538 A 11/995   Kanish et al.   Si 1106397   11/996   5,343,548 A 11/995   Kanish et al.   Si 1107854   41/996   Shapiro   WO 930007377   11/993   Shapiro	5,084,050 A	1/1992	Draenert	EP	0421485 A1	10/1998
5,108,319 A * 441992   Green et al.   FR   295 729   71,1976     5,108,422 A   441992   Green et al.   FR   281 336   11,1986     5,112,130 A   5,1092   Kerolin et al.   FR   281 336   11,1986     5,112,130 A   5,1092   Kerolin et al.   FR   281 336   11,1986     5,112,130 A   6,1992   Kerolin et al.   FR   281 336   11,1986     5,132,232 A   6,1992   Kerlen et al.   GB   12,944 0   10,1972     5,132,232 A   6,1992   Boyd   GB   1531487   11,1978     5,131,278 A   12,1992   Boyd   GB   1531487   11,1978     5,132,237 A   12,1992   Boyd   GB   2082754   A   31,0982     5,122,337 A   1,3193   Brantigan   GB   2082754   A   31,0982     5,122,335 A   6,1993   Brantigan   GB   2164277   A   31,986     5,246,458 A   9,1993   Kambin   PR   60,43084   31,985     5,246,458 A   9,1993   Kambin   PR   60,43084   31,985     5,236,351 A   11,1993   Salib et al.   PR   60,43084   31,985     5,230,351 A   11,1993   Bagby   SE   106 101   7,1939     5,263,354 A   11,1994   Bagby   SE   106 101   7,1939     5,263,354 A   11,1994   Bagby   SE   106 101   7,1939     5,263,354 A   11,1994   Bagby   SE   106 101   1,1948     5,304,131 A   12,1981   Bammigan   SI   12,12734   31,1986     5,304,131 A   12,1981   Bammigan   SI   12,12734   31,1986     5,304,336 A   1,1995   Kagna et al.   WO   91,06,566   51,1991     5,433,541 A   8,1995   Seffee et al.   WO   91,06,566   51,1991     5,433,541 A   8,1995   Seffee et al.   WO   91,06,566   51,1991     5,433,541 A   1,1996   Kagna et al.   WO   91,06,566   51,1991     5,433,541 A   1,1996   Kagna et al.   WO   91,06,566   51,1991     5,443,541 A   8,1995   Seffee et al.   WO   91,06,566   51,1991     5,443,541 A   8,1995   Seffee et al.   WO   91,06,566   51,1991     5,443,541 A   1,1996   Kagna et al.   WO   91,06,566   51,1991     5,443,541 A   8,1995   Seffee et al.   WO   91,06,566   51,1991     5,443,541 A   1,1996   Kagna et al.   WO   91,06,566   51,1991     5,443,541 A   1,1996   Kagna et al.   WO   91,06,566   51,1991     5,443,541 A   1,1996   Kagna et al.   WO   91,06,566			Kirsch	ES	283078	5/1985
5,108,395 A * 4/1992   Laurain						
5,118,429						
S.112,336 A						
5,116,304   A 5/1992   Cadwell   GB						
5,122,139						
5,123,926 A						
5,146,933						
5.171,278 A   21/1992   Pishanodi   GB   2076657 A   21/1981     5.182,471 A   21/1993   Wilk   GB   2082754 A   31/1982     5.192,233 A   31/1993   Bantigan   GB   2164077 A   31/1986     5.192,337 A   31/1993   Bantigan   GB   2164077 A   31/1986     5.292,335 A   61/1993   Bantigan   PP   57-29348   21/1982     5.292,424 A   64.093   Anthon   PP   50-31766   21/1985     5.282,835 A   11/1993   Bantigan   PP   60-31766   21/1985     5.282,835 A   21/1994   Bantannan et al.   PP   61-12859   71/1987     5.262,972 A   21/1993   Bantigan   PP   61-12859   71/1987     5.282,835 A   21/1994   Ledergerber   SU   10/1985   81/1984     5.292,232 A   31/1994   Ledergerber   SU   11/19784   81/1984     5.292,232 A   31/1994   Vickerson et al.   SU   12/12774   31/1986     5.304,332 A   51/1994   Vickerson et al.   SU   12/12774   31/1986     5.304,332 A   51/1994   Vickerson et al.   SU   12/12774   31/1986     5.304,332 A   51/1994   Vickerson et al.   SU   12/12774   31/1986     5.305,307 A   21/1994   Vickerson et al.   SU   12/12774   31/1986     5.305,307 A   21/1994   Vickerson et al.   SU   12/12774   31/1986     5.305,307 A   21/1994   Vickerson et al.   SU   12/12774   31/1986     5.305,307 A   21/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986     5.306,309 A   41/1994   Vickerson et al.   SU   12/12774   31/1986   SU   12/1278   31/1986   SU   12/1278   31/1986   SU   12/1278   31/1986   SU   12/1278   31/1986   SU   12/1986						
S.183,471   A   2/1993   Wilk	5,146,933 A	9/1992	Boyd	GB	1531487	11/1978
5.192,237 A 3/1993   Brantigan   GB 2126094 A 3/1984     5.197,971 A 3/1993   Brantigan   GB 216477 A 3/1986     5.197,971 A 3/1993   Brantigan   GB 216477 A 3/1986     5.219,355 A 61993   Brantigan   P 60-43794   2/1982     5.242,443 A 9/1993   Sambin   P 60-43794   3/1985     5.264,684 A 9/1993   Graham   P 60-43984   3/1985     5.286,884 A 9/1993   Graham   P 61-122889   6/1986     5.286,884 A 9/1993   Graham   P 61-122889   6/1986     5.286,884 A 9/1993   Wife   SE 100-100     5.279,292 A 1/1994   Brantigan   SU 1107854   8/1984     5.282,885 A 2/1994   Ledergerber   SU 112460   1/1984     5.282,885 A 2/1994   Ledergerber   SU 1124760   1/1984     5.392,225 A 3/1994   Vickerson et al.   SU 1217374   3/1986     5.304,312 A 4/1994   Wagner et al.   WO 8401298   4/1984     5.313,962 A 5/1994   Obenchain   WO W0 9000007   1/1990     5.314,427 A 5/1995   Gobbe et al.   WO 9106266   5/1991     5.352,229 A 10/1994   Shapino   Mahews   128/898     5.370,697 A 12/1994   Bumgarther   Su 12/1944   Shapino   Sheridan   Su 12/1944   Shapino   Sheridan   Su 12/1944   Shapino   Sheridan   Su 12/1945   Sheridan   Su 12/1945   Sheridan   Su 12/1945   Sheridan   Su 12/1946   Shapino   Sheridan   Su 12/1946   Shapino   Sheridan   Su 12/1947   Sheridan   Su 12/1948   Sheridan   Su 12/1948   Sheridan   Su 12/1948   Sheridan   Su 12/1949   Sheridan   Su 12/1941   Sher	5,171,278 A	12/1992	Pisharodi	GB	2076657 A	12/1981
S.192,327 A	5,183,471 A	2/1993	Wilk	GB	2082754 A	3/1982
S.192,327 A	5,192,293 A	* 3/1993	Cartwright et al 606/172	GB	2126094 A	3/1984
S.197971 A 3/1993   Parodi et al.   IP   60-43984   21982						
S.219.355 A   6/1993   Rambin   IP   60-31706   21985						
S.242.443 A   9/1993   Graham   JP   60-43984   3/1985     S.248.031 A   1/1993   Graham   JP   61-122859   6/1986     S.258.033 A   1/1993   Bagby   SE   106 101   7/1939     S.269.772 A   12/1993   Wilk   SU   1063397   12/1983     S.269.772 A   12/1993   Wilk   SU   1063397   12/1983     S.282.856 A   2/1994   Edergerber   SU   1124960   11/1984     S.282.855 A   2/1994   Nickerson et al.   SU   117374   3/1986     S.304.132 A   4/1994   Wagner et al.   WO   8401298   4/1984     S.313.962 A   5/1994   Goble et al.   WO   WO   9000037   1/1990     S.314.27 A   S/1994   Goble et al.   WO   9000037   1/1990     S.352.29 A   10/1994   Mathews   128/898     S.364.399 A   11/1994   Mathews   128/898     S.363.90 A   1/1995   Michelson   S.393.036 A   3/1995   Kagan et al.   Cloward, Rajhb B.; Ruptured Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs & Springer Springer & Osteophytes & Anterior Cervical Intervertebral Discs & Springer & Osteophytes & Anterior Cervical Inte						
S.246,458 A   9/1993   Sailb et al.   IP   61-122859   61/986						
S.258,031 A   11/1993   Bagbry   SE   106 101   7/1939						
S.263,953 A						
S.209.772 A   12/1994 Baumann et al.   SU   1063397   12/1983						
S.279,292 A					106 101	7/1939
S282,856 A   2/1994   Cedergerber   SU   1124960   11/1984	5,269,772 A	12/1993	Wilk	su	1063397	12/1983
S.282,856 A   21/994   Gedegerber   SU   112/960   11/1984	5,279,292 A	1/1994	Baumann et al.	su	1107854	8/1984
Signature Serial   SU   1217374   31/986						
Signature   Sign	5 292 252 A	3/1994				
S,306,309 A 4/1994   Wagner et al.		4/1994				
S,313,962 A   5/1994 Goble et al.   WO   WO 90/00037   1/1990		4/1004	Wagner et al			
S,314,427 A 5/1994   Goble et al.   WO 91/06266   5/1991		5/1004				
S,324,295 A 6 (1994   Shapiro   WO 92/14423   9/1992						
Signature Signature Seried   Spine processed   Spine processed   Spine processed   Spine spine street   Spine street   Spine spine street   Spine street   Spine street   Spine spine street   Spine st						
5.357,983 A * 10/1994   Mathews						
S.364,399 A   11/1994   Lowery et al.   S.370,697 A   12/1994   Baumgartner   Cloward, Ralph B.; Surgical Techniques for Lumbar Disc Lesions; Codman; Signature Serial 3.   Cloward, Ralph B.; Ruptured Cervical Intervertebral Discs: Removal 5,395,317 A   3/1995   Kagan et al.   Cloward, Ralph B.; Ruptured Cervical Intervertebral Discs: Removal 6,397,364 A   3/1995   Kozak et al.   Gloward, Ralph B.; Ruptured Cervical Intervertebral Discs: Removal 6,1925   Kagan et al.   Gloward, Ralph B.; Ruptured Cervical Intervertebral Discs: Removal 6,1925   Kagan et al.   Gloward, Ralph B.; Ruptured Cervical Intervertebral Discs: Removal 6,1925   Kagan et al.   Gloward, Ralph B.; Recent Advances in Surgery of the Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Interbody Fusion with Threaded Cylindrical Bone; pp. 750-753; Journal of Neurosurgery, vov. 1985, vol. 63, No. 5.    15, 15, 15, 15, 15, 15, 15, 15, 15, 15,				WO	93/01771	2/1993
Signature   Sign	5,357,983 A	* 10/1994	Mathews 128/898		OTHER BUIL	OI ICATIONG
S.393,036 A   2/1995   Sheridan   Cloward, Raiph B.; Surgical Techniques for Lumbar Disc Lesions; S.395,317 A   3/1995   Kagan et al.   Cloward, Raiph B.; Surgical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Interbody Fusion (A.C. S.425,772 A   6/1995   Brantigan   Cloward, Ralph B.; Rueptured Cervical Interbody Fusion (A.C. S.425,772 A   6/1995   Brantigan   Cloward, Ralph B.; Recent Advances in Surgery of the Cervical S.415,157 A   8/1995   Steffee   Spine; pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical S.443,514 A   8/1995   Kuslich et al.   Spine; Operations; Excerpta Medica.   Otero-Vich, Jose M.; Anterior Cervical Interbody Fusion with S.489,307 A   2/1996   Kuslich et al.   Nov. 1985, vol. 63, No. 5.   Stensish and Posterior Lumbar S.505,732 A   4/1996   Michelson   Hutter, Charles George; Spinal Stenosis and Posterior Lumbar S.5797,909 A   8/1998   Michelson   Hutter, Charles George; Spinal Stenosis and Posterior Lumbar G.966,912 B2   11/2005   Michelson   Research; No. 193; The Association of Bone and Joint Surgeons.   Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122; Charles C. Thomas; Springfield, Illinois.   Charles C. Thomas; Springfield, Illinois.   Lin, Paul M.; Lumbar Interbody Fusion; pp. 81, 98, 120, 146, 7,326,214 B2   2/2008   Michelson   Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 7,387,365 B2   2/2011   Michelson   Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.   2006/0058793   Al. 3/2006   Michelson   Surgical Techniques of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques (Pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.   Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 83-91; Springe	5,364,399 A	11/1994	Lowery et al.		OTHER PU	BLICATIONS
S.393,036 A   2/1995   Sheridan   Cloward, Raiph B.; Surgical Techniques for Lumbar Disc Lesions; S.395,317 A   3/1995   Kagan et al.   Cloward, Raiph B.; Surgical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Intervertebral Discs: Removal of Disc & Osteophytes & Anterior Cervical Interbody Fusion (A.C. S.425,772 A   6/1995   Brantigan   Cloward, Ralph B.; Rueptured Cervical Interbody Fusion (A.C. S.425,772 A   6/1995   Brantigan   Cloward, Ralph B.; Recent Advances in Surgery of the Cervical S.415,157 A   8/1995   Steffee   Spine; pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical S.443,514 A   8/1995   Kuslich et al.   Spine; Operations; Excerpta Medica.   Otero-Vich, Jose M.; Anterior Cervical Interbody Fusion with S.489,307 A   2/1996   Kuslich et al.   Nov. 1985, vol. 63, No. 5.   Stensish and Posterior Lumbar S.505,732 A   4/1996   Michelson   Hutter, Charles George; Spinal Stenosis and Posterior Lumbar S.5797,909 A   8/1998   Michelson   Hutter, Charles George; Spinal Stenosis and Posterior Lumbar G.966,912 B2   11/2005   Michelson   Research; No. 193; The Association of Bone and Joint Surgeons.   Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122; Charles C. Thomas; Springfield, Illinois.   Charles C. Thomas; Springfield, Illinois.   Lin, Paul M.; Lumbar Interbody Fusion; pp. 81, 98, 120, 146, 7,326,214 B2   2/2008   Michelson   Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 7,387,365 B2   2/2011   Michelson   Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.   2006/0058793   Al. 3/2006   Michelson   Surgical Techniques of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques (Pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.   Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 83-91; Springe	5,370,697 A	12/1994	Baumgartner	C1 1	D 1 1 D G 2 1 T	1 ' C T 1 D' T '
S.395.317 A						nniques for Lumbar Disc Lesions;
5,396,880         A         3/1995         Kagan et al.         Cloward, Ralph B.; Ruptured Cervical Interbody Fusion (A.C. 5,397,364         A         3/1995         Kozak et al.         of Disc & Osteophytes & Anterior Cervical Interbody Fusion (A.C. 5,423,842         A         6/1995         Michelson         I.F.); Codman; Signature Series 4.           5,4425,772         A         6/1995         Brantigan         Cloward, Ralph B.; Recent Advances in Surgery of the Cervical Spine, pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine, opp. 391, 308.           5,485,638         A         10/1996         Michelson         Hutter, Old All All All All All All All All All A				Codman;	Signature Serial 3.	
5,397,364 A 3/1995 Kozak et al. 5,425,772 A 6/1995 Michelson 5,425,772 A 6/1995 Erickson et al. 5,441,527 A 8/1995 Erickson et al. 5,443,514 A 8/1995 Erickson et al. 5,458,638 A 10/1995 Kuslich et al. 5,488,437 A 1/1996 Michelson 5,489,307 A 2/1996 Kuslich et al. 5,489,308 A 2/1996 Kuslich et al. 5,505,732 A 4/1996 Michelson 5,797,909 A 8/1998 Michelson 6,436,098 B1 8/2000 Michelson 6,436,098 B1 8/2000 Michelson 6,436,098 B1 8/2000 Michelson 6,966,912 B2 1/12005 Michelson 6,966,912 B2 1/12005 Michelson 7,264,622 B2 9/2007 Michelson 7,264,622 B2 9/2007 Michelson 7,326,214 B2 2/2008 Michelson 7,330,301 B2 2/2008 Michelson 7,399,303 B2 7/2008 Michelson 7,387,565 B2 2/2011 Michelson 7,487,565 B2 2/2011 Michelson 7,887,565 B2 2/2011 Mic				Cloward.	Ralph B.: Ruptured Cer	vical Intervertebral Discs: Removal
5,423,842         A         6/1995         Michelson         I.F.); Codman; Signature Series 4.           5,425,772         A         6/1995         Brantigan         Cloward, Ralph B.; Recent Advances in Surgery of the Cervical Spine; pp. 285-293; German Society for Neurosurgery; vol. 2 Cervical Spine Operations; Excerpta Medica.           5,443,514         A         8/1995         Steffee         cal Spine Operations; Excerpta Medica.           5,483,638         A         1/1996         Michelson         Threaded Cylindrical Bone; pp. 750-753; Journal of Neurosurgery, Nov. 1985, vol. 63, No. 5.           5,489,307         A         2/1996         Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,505,732         A         4/1996         Michelson         Hutter, Charles George; Spinal Stenosis and Posterior Lumbar Interbody Fusion: pp. 103-114; Clinical Orthopaedics and Related Research; No. 193; The Association of Bone and Joint Surgeons.           6,436,098         B1         8/2000         Michelson         Lin, Paul M.; Posterior Lumbar Interbody Fusion: pp. 114-122;           6,875,213         B2         1/2005         Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in 7,264,622           B2         1/2005         Michelson         Spine Surgery; Techniques and Complications: pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,887,565         B						
5,425,772         A         6/1995         Brantigan         Cloward, Ralph B.; Recent Advances in Surgery of the Cervical Spine; pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Advances in Surgery of the Cervical Spine; pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine Operations; Excerpta Medica.           5,484,337 A         1/1995         Kuslich et al.         Otero-Vich, Jose M.; Anterior Cervical Interbody Fusion with Threaded Cylindrical Bone; pp. 750-753; Journal of Neurosurgery, Nov. 1985, vol. 63, No. 5.           5,489,308 A         2/1996         Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,505,732 A         4/1996         Michelson         Hutter, Charles George; Spinal Stenosis and Posterior Lumbar Interbody Fusion; pp. 103-114; Clinical Orthopaedics and Related 6,096,038 A         8/2000         Michelson         Research; No. 193; The Association of Bone and Joint Surgeons.           6,436,098 B1 Extraction States Stat						• ,
5,441,527 A         8/1995         Erickson et al.         Spine; pp. 285-293; German Society for Neurosurgery: vol. 2 Cervical Spine Operations; Excerpta Medica.           5,483,638 A         10/1995         Kuslich et al.         Otero-Vich, Jose M.; Anterior Cervical Interbody Fusion with Threaded Cylindrical Bone; pp. 750-753; Journal of Neurosurgery, Nov. 1985, vol. 63, No. 5.           5,489,307 A         2/1996         Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,505,732 A         4/1996         Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,797,909 A         8/1998         Michelson         Hutter, Charles George; Spinal Stenosis and Posterior Lumbar Interbody Fusion; pp. 103-114; Clinical Orthopaedics and Related Research; No. 193; The Association of Bone and Joint Surgeons.           6,436,098 B1         8/2002         Michelson         Research; No. 193; The Association of Bone and Joint Surgeons.           6,966,912 B2         11/2005         Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Pp. 114-122;           6,7326,214 B2         2/2007         Michelson         Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146           7,887,565 B2         2/2011         Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2008/0058793 A1         1/2008         Michelson         Michelson					•	
5,443,514 A 8/1995 Steffee 5,458,638 A 10/1996 Musich et al. 5,484,437 A 1/1996 Michelson 5,489,307 A 2/1996 Kuslich et al. 5,505,732 A 4/1996 Michelson 5,797,909 A 8/1998 Michelson 6,096,038 A 8/2000 Michelson 6,436,098 B1 8/2002 Michelson 6,966,912 B2 11/2005 Michelson 6,966,912 B2 11/2005 Michelson 6,966,912 B2 11/2005 Michelson 7,264,622 B2 9/2007 Michelson 7,326,214 B2 2/2008 Michelson 7,399,303 B2 7/2008 Michelson 7,387,565 B2 2/2011 Michelson 7,887,565 B2 2/2011 Michelson 2006/0036247 A1 2/2006 Michelson 2006/0036247 A1 2/2006 Michelson 2006/0036247 A1 2/2006 Michelson 2006/0036247 A1 2/2006 Michelson 2006/0038793 A1 3/2006 Michelson 2006/0038793 A1 3/2006 Michelson 2008/0257564 A1 10/2008 Michelson  FOREIGN PATENT DOCUMENTS  DE 2446 039 4/1975 DE 3101333 A1 12/1981 DE 3101233 A1 12/1981 DE 31012						
5,458,638 A 10/1995 Kuslich et al. 5,484,437 A 1/1996 Michelson 5,489,307 A 2/1996 Kuslich et al. 5,489,308 A 2/1996 Kuslich et al. 5,505,732 A 4/1996 Michelson 5,797,909 A 8/1998 Michelson 6,096,038 A 8/2000 Michelson 6,436,098 B1 8/2002 Michelson 6,436,098 B1 8/2002 Michelson 6,966,912 B2 11/2005 Michelson 7,264,622 B2 9/2007 Michelson 7,264,622 B2 9/2007 Michelson 7,326,214 B2 2/2008 Michelson 7,326,214 B2 2/2008 Michelson 7,399,303 B2 7/2008 Michelson 7,887,565 B2 2/2011 Michelson 8,180,184,204,224,225,231, Aspen Publishers, Inc.; 1989. 7,887,565 B2 2/2014 Michelson 8,180,184,204,224,225,231, Aspen Publishers, Inc.; 1989. 7,887,565 B2 2/2014 Michelson 8,180,184,204,224,225,231, Aspen Publishers, Inc.; 1989. 7,887,565 B2 2/2016 Michelson 8,180,184,204,224,225,231, Aspen Publishers, Inc.; 1989. 8,194,184,184,184,184,184,184,184,184,184,18				Spine; pp.	. 285-293; German Soci	ety for Neurosurgery: vol. 2 Cervi-
5,484,437 A         1/1996 Michelson         Oten-Victal, Jose Mr., Alteriol Cevitcal Interbody Pusion (Proposition)           5,489,307 A         2/1996 Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,489,308 A         2/1996 Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,505,732 A         4/1996 Michelson         Hutter, Charles George; Spinal Stenosis and Posterior Lumbar 1nterbody Fusion; pp. 103-114; Clinical Orthopaedics and Related Research; No. 193; The Association of Bone and Joint Surgeons.           6,436,098 B1 8/2002 Michelson         Michelson         Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122; Charles C. Thomas; Springfield, Illinois.           6,966,912 B2 1/2005 Michelson         Charles C. Thomas; Springfield, Illinois.           7,326,214 B2 2/2008 Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 7326, 214 B2 2/2008 Michelson           7,887,565 B2 2/2011 Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2006/0036247 A1 2/2006 Michelson         Michelson           2008/0287955 A1 11/2008 Michelson         Michelson           DE 310333 A1 1/2008 Michelson         Michelson           DE 310333 A1 1/2008 Michelson         Michelson           DE 310333 A1 1/2018 Michelson         Michelson				cal Spine	Operations; Excerpta N	Iedica.
5,484,437 A 1/1996 Kuslich et al. 5,489,307 A 2/1996 Kuslich et al. 5,489,308 A 2/1996 Kuslich et al. 5,505,732 A 4/1996 Michelson Hutter, Charles George; Spinal Stenosis and Posterior Lumbar Interbody Fusion; pp. 103-114; Clinical Orthopaedics and Related Research; No. 193; The Association of Bone and Joint Surgeons. 6,436,098 B1 8/2002 Michelson Research; No. 193; The Association of Bone and Joint Surgeons. 6,436,098 B1 8/2002 Michelson Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122; 6,875,213 B2 4/2005 Michelson Charles C. Thomas; Springfield, Illinois. 6,966,912 B2 1/2005 Michelson Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 7,3264,622 B2 9/2007 Michelson Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 7,326,214 B2 2/2008 Michelson Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989. 7,887,565 B2 2/2011 Michelson Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.  Michelson Michelson Surgery: Techniques, vol. 5, No. 3, 1990.  Michelson Michelson Michelson Surgical Techniques, vol. 5, No. 3, 1990.  Michelson Michelson Michelson Michelson Surgical Techniques, vol. 5, No. 3, 1990.  FOREIGN PATENT DOCUMENTS Por 3101333 A1 12/1981 Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1979.  Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1979.  Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.  BE 3101333 A1 12/1981 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Disk Syndromes of the Low				Otero-Vic	th. Jose M.: Anterior	Cervical Interbody Fusion with
5,489,308 A         2/1996         Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,489,308 A         2/1996         Kuslich et al.         Nov. 1985, vol. 63, No. 5.           5,505,732 A         4/1996         Michelson         Hutter, Charles George; Spinal Stenosis and Posterior Lumbar Interbody Fusion; pp. 103-114; Clinical Orthopaedics and Related Research; No. 193; The Association of Bone and Joint Surgeons.           6,436,098 B1         8/2002         Michelson         Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122; Charles C. Thomas; Springfield, Illinois.           6,966,912 B2         11/2005         Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,326,214 B2         2/2008         Michelson         Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,899,303 B2         7/2008         Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2006/0036247 A1         2/2006         Michelson         Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           DE         3101333 A1						,
Sistor   S	5,489,307 A					150 155, Journal of Neurosurgery,
5,505,732 A         4/1996 Michelson         Hutter, Charles George; Spinal Stenosis and Posterior Lumbar Interbody Fusion; pp. 103-114; Clinical Orthopaedics and Related Research; No. 193; The Association of Bone and Joint Surgeons.           6,966,038 A         8/2000 Michelson         Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122; Charles C. Thomas; Springfield, Illinois.           6,875,213 B2         4/2005 Michelson         Lin, Paul M.; Lumbar Interbody Fusion; pp. 114-122; Charles C. Thomas; Springfield, Illinois.           6,966,912 B2         11/2005 Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,399,303 B2         7/2008 Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2004/0073217 A1         4/2004 Michelson         Michelson           2006/0036247 A1         2/2006 Michelson         Surgical Techniques, vol. 5, No. 3, 1990.           Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           DE         34 46 039 (4)1975         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3101333 A1         12/1981	5,489,308 A	2/1996	Kuslich et al.			C I B
6,096,038 A 8/2000 Michelson 6,436,098 B1 8/2002 Michelson 6,436,098 B1 8/2002 Michelson 6,966,912 B2 11/2005 Michelson 7,264,622 B2 9/2007 Michelson 7,326,214 B2 2/2008 Michelson 7,399,303 B2 7/2008 Michelson 7,887,565 B2 2/2011 Michelson 7,889,4 Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990. Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979. Mierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1975. Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974. Mierholzer, G.; Manual on	5,505,732 A	4/1996	Michelson			
6,096,038 A         8/2000 Michelson         Research; No. 193; The Association of Bone and Joint Surgeons.           6,436,098 B1         8/2002 Michelson         Michelson         Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122;           6,875,213 B2         4/2005 Michelson         Charles C. Thomas; Springfield, Illinois.         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,393,303 B2         7/2008 Michelson         Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2004/0073217 A1         4/2004 Michelson         Michelson         Surgical Techniques, vol. 5, No. 3, 1990.         Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           DE         24 46 039 Michelson         Michelson         Michelson         Michelson         Michelson           DE         3 101333 A1         12/1981         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3 46 08 163 A1         6/1982         Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to	5,797,909 A					•
Lin, Paul M.; Posterior Lumbar Interbody Fusion; pp. 114-122;	· · · · · · · · · · · · · · · · · · ·			Research;	No. 193; The Associat	ion of Bone and Joint Surgeons.
6,875,213         B2         4/2005         Michelson         Charles C. Thomas; Springfield, Illinois.           6,966,912         B2         11/2005         Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,399,303         B2         7/2008         Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2006/0036247         A1         2/2006         Michelson         Michelson           2008/0255564         A1         10/2008         Michelson         Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           DE         24 46 039         4/1975         Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.           DE         3132520         A1         6/1982         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         36 08 163         A1         9/1987         Hermon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127;				Lin. Paul	M.; Posterior Lumba	r Interbody Fusion; pp. 114-122;
6,966,912         B2         11/2005         Michelson         Lin, Paul M.; Lumbar Interbody Fusion: Principles and Techniques in Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 7,326,214         B2         9/2007         Michelson         Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 180-184, 204, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,399,303         B2         7/2008         Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2004/0073217         A1         4/2004         Michelson         Surgical Techniques, vol. 5, No. 3, 1990.           2006/0036247         A1         2/2006         Michelson         Michelson         Michelson           2008/0255564         A1         10/2008         Michelson         Michelson         Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           DE         24 46 039         4/1975         Heirholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.           DE         3101333         A1         12/1981         Hyrical Allograft; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3105567         A1         6/1982         Harmon, Paul H.; Anterior						
7,264,622         B2         9/2007         Michelson         Spine Surgery; Techniques and Complications; pp. 81, 98, 120, 146, 173, 280, 224, 225, 231; Aspen Publishers, Inc.; 1989.           7,399,303         B2         7/2008         Michelson         Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.           2004/0073217         A1         4/2004         Michelson         Surgical Techniques, vol. 5, No. 3, 1990.           2006/0036247         A1         2/2006         Michelson         Michelson           2008/0255564         A1         10/2008         Michelson         Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           DE         24 46 039         4/1975         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3101333         A1         12/1981         Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini-20 Company Call Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-						
7,326,214 B2						
7,399,303 B2 7/2008 Michelson 7,887,565 B2 2/2011 Michelson 2004/0073217 A1 4/2004 Michelson 2006/0036247 A1 2/2006 Michelson 2006/0058793 A1 3/2006 Michelson 2008/0255564 A1 10/2008 Michelson 2008/02575564 A1 11/2008 Michelson  FOREIGN PATENT DOCUMENTS  DE 24 46 039 4/1975 DE 3101333 A1 12/1981 DE 31032520 A1 6/1982 DE 3505567 A1 6/1986 DE 36 08 163 A1 9/1987 DE 36 08 163 A1 9/1987 DE 36 08 163 A1 9/1987 DE 4 036 804 5/1992  Michelson Tan, S.B.; A Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990. Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979. Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985. Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974. Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical And Modified Technique of Anterior Lumbar Fusion with Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Action Surgical Actio				-	~ .	
7,887,565 B2 2/2011 Michelson 2004/0073217 A1 4/2004 Michelson 2006/0036247 A1 2/2006 Michelson 2006/0058793 A1 3/2006 Michelson 2008/0257564 A1 10/2008 Michelson 2008/0287955 A1 11/2008 Michelson EVALUATION FOREIGN PATENT DOCUMENTS  DE 24 46 039 4/1975 DE 3101333 A1 12/1981 DE 3132520 A1 6/1982 DE 3505567 A1 6/1986 DE 36 08 163 A1 9/1987 DE 36 08 163 A1 9/1987 DE 36 08 163 A1 9/1987 DE 4 036 804 5/1992  Tanh, S.B., A Modified reclinique of Anterior Limbal Putson With Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic Surgical Techniques, vol. 5, No. 3, 1990.  Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979. Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.  Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974. Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Company of the Company of th				173, 180	184, 204, 224, 225, 231	; Aspen Publishers, Inc.; 1989.
Femoral Cortical Allograft; pp. 83-93; The Journal of Orthopaedic				Tan, S.B.:	A Modified Technique	e of Anterior Lumbar Fusion with
2004/00/36247         A1         4/2004         Michelson         Surgical Techniques, vol. 5, No. 3, 1990.           2006/0058793         A1         3/2006         Michelson         Muller, M.E.; Manual of Internal Fixation: Techniques Recommended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           EVALUATION OF ALTERNY DOCUMENTS         POREIGN PATENT DOCUMENTS         Phierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.           DE         24 46 039         4/1975         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3101333         A1         12/1981         Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Phieracters of the August of Intervent Company of Interve					•	
2006/0058793 A1 3/2006 Michelson 2008/0255564 A1 10/2008 Michelson 2008/02575554 A1 11/2008 Michelson 2008/02575554 A1 11/2008 Michelson EVALUATION FOREIGN PATENT DOCUMENTS  DE 24 46 039 4/1975 DE 3101333 A1 12/1981 Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974. DE 3132520 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Company of the Asif-Group of the Asif-Group of the Asign o	2004/0073217 A	1 4/2004	Michelson			
2008/0257954         A1         3/2006         Michelson         mended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979.           FOREIGN PATENT DOCUMENTS           DE         24 46 039         4/1975         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3101333         A1         12/1981         Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion DE           DE         3505567         A1         6/1986         Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE           DE         4 036 804         5/1992         Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-	2006/0036247 A	1 2/2006	Michelson	-	•	
2008/0255564         A1         10/2008         Michelson         mended by the AO Group; Second Edition, Expanded and Revised; pp. 3-20, 27-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979. Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.           DE         24 46 039         4/1975         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3101333         A1         12/1981         Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-	2006/0058793 A	1 3/2006	Michelson			
2008/0287955         A1         11/2008         Michelson         pp. 3-20, 2/-41, 53-58, 71-78, 94, 311, 320; Springer-Verlag; 1979. Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator; pp. 85-91; Springer-Verlag; 1985.           DE         24 46 039         4/1975         Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.           DE         3132520         A1         6/1982         Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Policy Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clinical Policy Spine: Cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-					•	
Hierholzer, G.; Manual on the AO/ASIF Tubular External Fixator;  FOREIGN PATENT DOCUMENTS  DE 24 46 039 4/1975 Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.  DE 3101333 A1 12/1981 by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.  DE 3505567 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini-DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-				pp. 3-20, 2	27-41, 53-58, 71-78, 94	l, 311, 320; Springer-Verlag; 1979.
FOREIGN PATENT DOCUMENTS  DE 24 46 039 4/1975 Heim, Urs; Small Fragment Set Manual: Technique Recommended DE 3101333 A1 12/1981 by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974. DE 3132520 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion DE 3505567 A1 6/1986 Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini- DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-	2000/020/933 A	11/2008	MICHOISON			
DE 24 46 039 4/1975 Heim, Urs; Small Fragment Set Manual: Technique Recommended by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974. DE 3132520 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion DE 3505567 A1 6/1986 Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini-DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-	FOREIGN PATENT DOCUMENTS					
DE 3101333 A1 12/1981 by the ASIF-Group; pp. 5-7, 10, 20, 21, 30; Springer-Verlag; 1974.  DE 3132520 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion  DE 36 08 163 A1 9/1987 Operation for Intervertebral Disk Syndromes of the Lower Lumbar  DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini-  DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-				* *		
DE 3132520 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion DE 3505567 A1 6/1986 Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini- DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-		4 46 039	4/1975			
DE 3132520 A1 6/1982 Harmon, Paul H.; Anterior Excision and Vertebral Body Fusion DE 3505567 A1 6/1986 Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini- DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-	DE 3	3101333 A1	12/1981			
DE 3505567 A1 6/1986 Operation for Intervertebral Disk Syndromes of the Lower Lumbar DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini- DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-						
DE 36 08 163 A1 9/1987 Spine: Three- to Five-Year Results in 244 Cases; pp. 107-127; Clini- DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-				Operation	for Intervertebral Disk	Syndromes of the Lower Lumbar
DE 4 036 804 5/1992 cal Orthopaedics and Related Research, No. 26, J.B. Lippincott Com-						
				•		
рашу, 1905.				-		
	~~ T.			Pany, 190.	٥.	

Harmon, Paul H.; A Simplified Surgical Technic for Anterior Lumbar Diskectomy and Fusion; Avoidance of Complications; Anatomy of the Retroperitoneal Veins; pp. 130-143; Clinical Orthopaedics and Related Research, No. 37, J.B. Lippincott Company, 1964.

Bullough, Peter G.; Atlas of Spinal Diseases; Figure 5.7; J.B. Lippencott Company; 1988.

Butts, M.K.; Biomechanical Analysis of a New Method for Spinal Interbody Fixation; 1987 Symposium, American Society of Mechanical Engineers, "Advance in Bioengineering", Boston, MA (Dec. 13-18, 1987).

Canale, S. Terry; Campbell's Operative Orthopaedics; vol. 3, 9<sup>th</sup> Edition; pp. 2191, 2216, 2459; Mosby, 1998.

Bagby, George W.; Wobbler Syndrome in Horses (the Ataxic Horse); Spokane County Medical Society Bulletin; Spring 1979.

Albrektsson, T., et al.; Osseointegrated Titanium Implants; Acta. Orthop. Scand.; vol. 52:155-170 (1981).

Raveh, J., et al.; Neue Rekonstruktionsmoglichkeiten des Unterkiefers bei knochernen Defekten nach Tumorresektionen; Der Chirurg vol. 53:459-467 (1982).

DeBowes, R.M., et al.; Study of Bovine . . . Steel Baskets; Transactions of the 29th Annual Meeting; Orthopaedic Research Society, vol. 8, p. 407, Mar. 8-10 (1983).

Brandt, L., et al.; A Dowel Inserter for Anterior Cervical Interbody Fusion; J. Neurosurg. 61:793-794 (Oct. 1984).

Raveh, J., et al.; Use of the Titanium-coated Hollow Screw and Reconstruction Plate System in Bridging of Lower Jaw Defects; J. Oral Maxillofac Surg. 42:281-294 (1984).

Otero-Vich, Jose M.; Anterior Cervical Interbody Fusion with Threaded Cylindrical Bone; J. Neurosurg 63:750-753 (Nov. 1985). Morscher, E., et al.; Die vordere Verplattung der Halswirbelsäule mit dem Hohlschrauben-Plattensystem aus Titanium, *Der Chirurg*, vol.

57, pp. 702-707 (1986) with English Translation. Bagby, G.W.; Basket Implant Facilitates Spinal Fusion; Orthopedics Today, vol. 7, No. 10, (Oct. 1987).

Butts, M. K., et al.; Biomechanical Analysis of a New Method for Spinal Interbody Fixation; 1987 Symposium, American Society of Mechanical Engineers, "Advances in Bioengineering", Boston, MA (Dec. 13-18, 1987).

Crawley et al.; A Modified Cloward's Technique for Arthrodesis of the Normal Metacarpophalangeal Joint in the Horse; Veterinary Surgery, vol. 17, No. 3, pp. 117-127 (1988).

Raveh, J., et al.; Surgical Procedures for Reconstruction of the Lower Jaw Using the Titanium-Coated Hollow-Screw Reconstruction Plate System: Bridging of Defects; Otolaryngologic Clinics of North America; vol. 20, No. 3 (Aug. 1987).

Goldthwaite, N., et al.; Toward Percutaneous Spine Fusion; Ch. 45; Lumbar Spine Surgery; C.V. Mosby Company, pp. 512-522 (1987). Bagby, G.W.; Arthrodesis by the Distraction-Compression Method Using a Stainless Steel Implant; Orthopedics, vol. II, No. 6, pp. 931-934 (Jun. 1987).

Itoman, M., et al.; Banked Bone Grafting for Bone Defect Repair—Clinical Evaluation of Bone Union and Graft Incorporation; J. Jpn. Orthop. Assoc. 62:461-469 (1988).

Crock, H.V.; "Anterior Lumbar Interbody Fusion—Indications for its Use and Notes on Surgical Technique;" *Clinical Orthopedics*, vol. 165; 4 cover pages and pp. 157-163; May 1982.

Wiltberger, et al.; "Dowel Intervertebral Fusion as Used in Lumbar Disc Surgery;" *The Journal of Bone and Joint Surgery*, vol. 39A, cover page and pp. 284-292 and 331; Apr. 1957.

Harris, P., et al.; Spinal Deformity After Spinal Cord Injury; Paraplegia, vol. 6, No. 4, pp. 232-238 (Feb. 1969).

Gillingham, F.J., et al.; Head injuries; Proceedings of the 18<sup>th</sup> World Congress of the International College of Surgeons, Rome, pp. 68-71 (May 28-31, 1972).

Whatmore, W. J.; Sincipital Encephalomeningoceles; Brit. J. Surg., vol. 60, No. 4, pp. 261-270 (Apr. 1973).

Whatmore, W. J.; Meningioma Following Trauma; Brit. J. Surg., vol. 60, No. 6, pp. 496-498 (Jun. 1973).

Rathke, F.W., et al.; Surgery of the Spine; Atlas of Orthopaedic Operations, vol. 1, p. 137, W.B. Saunders Co., Philadelphia (1979). Crock, H.V.; Practice of Spinal Surgery; Springer-Verlag/Wien, New York (1983).

O'Neill, P., et al.; Spinal Meningoceles in Association with Neurofibromatosis; Neurosurgery, vol. 13, No. 1, pp. 82-84 (Jul. 1983)

Whatmore, W.J., et al.; The Coventry Cervical Spreader and Dowel Inserter; ACTA Neurochirurgica, vol. 70, FASC. 1-2 (1984).

Whatmore, W. J.; Proceedings of the Society of British Neurological Surgeons; Journal of Neurology, Neurosurgery, and Psychiatry, 50:1093-1100 (1987).

Itoman, M., et al.; Banked Bone Grafting for Bone Defect Repair—Clinical Evaluation of Bone Union and Graft Incorporation; J. Jpn. Orthop. Assoc. 62:461-469 (1988).

Kane, W.J.; Direct Current Electrical Bone Growth Stimulation for Spinal Fusion; Spine, vol. 13, No. 3, pp. 363-365 (Mar. 1988).

The SpF-T Spinal Fusion Stimulator: An Efficacious Adjunct that Meets the Diverse Needs of Spine Patients; EBI Medical Systems; (Aug. 1991).


Schmitz et al.; Performance of Alloplastic Materials and Design of an Artificial Disc; The Artificial Disc, Brock, Mayer, Weigel; pp. 23-34 (1991).

The Use of Direct Current for Electrically Induced Osteogenesis; The Positive Effect of an Electronegative charge on Bone Growth; EBI Medical Systems (Feb. 1993).


Mylonas, C., et al.; Anterior Cervical Decompression and Fusion Using the Coventry Cervical Spreader and Dowel Inserter; British Journal of Neurosurgery, 7:545-549 (1993).


Fusion of the Lumbar Spine; Anterior Monosegmental Fusion L5-S1, Atlas of Spinal Operations, Thieme, pp. 270-274 (1993). Spine Basics, Danek Group, Inc., Glossary (1993).


\* cited by examiner

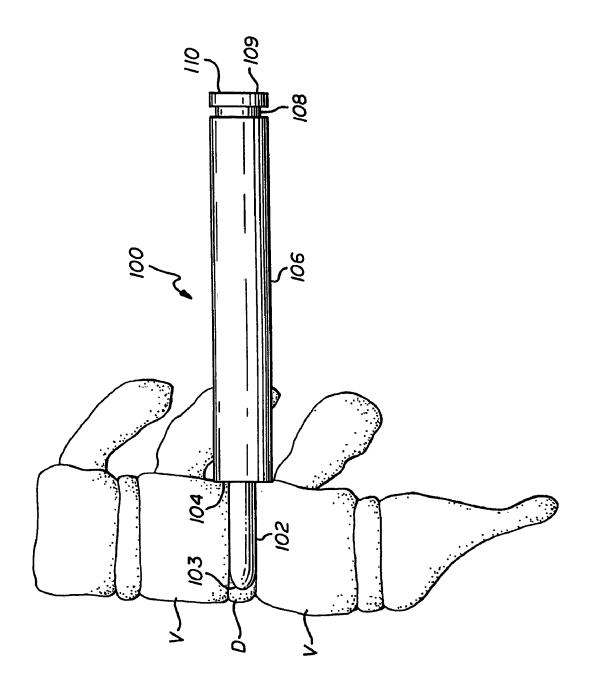

US 7,993,347 B1

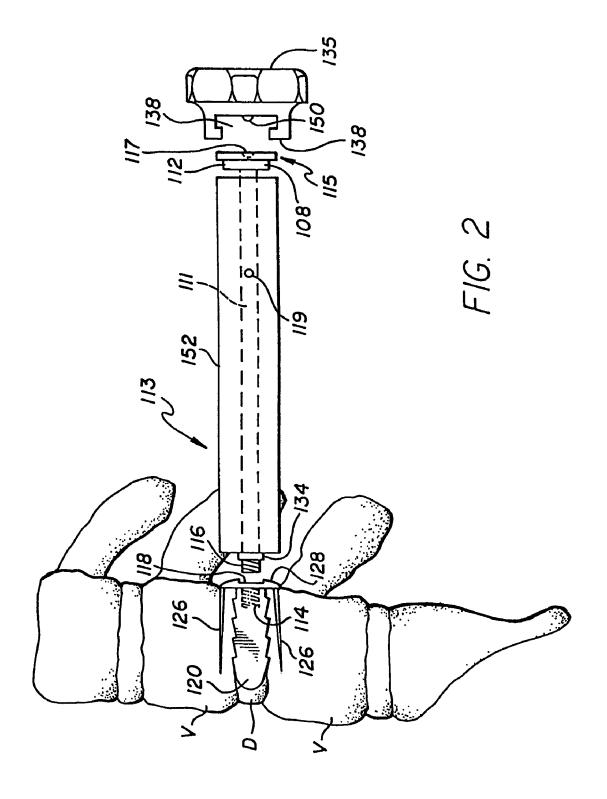

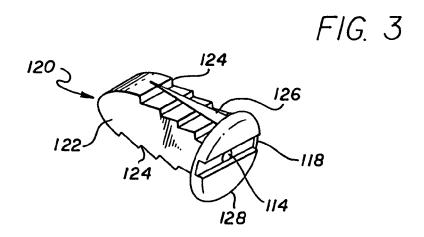

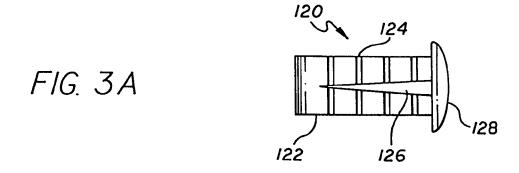

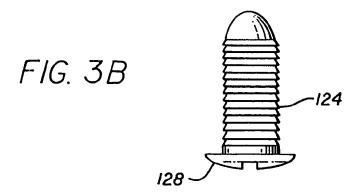


FIG. 3C



124 FIG. 3D

FIG. 3E



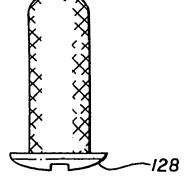
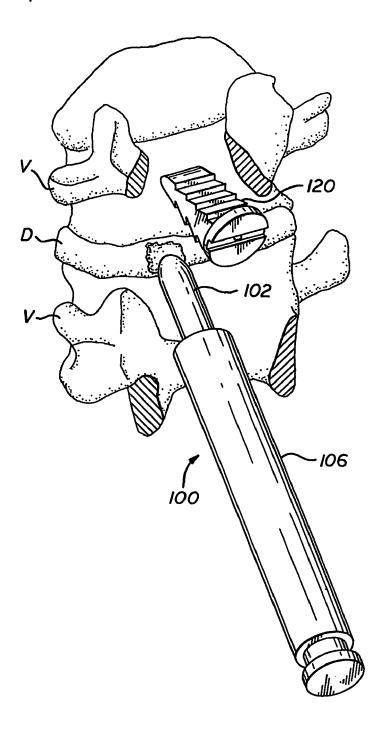
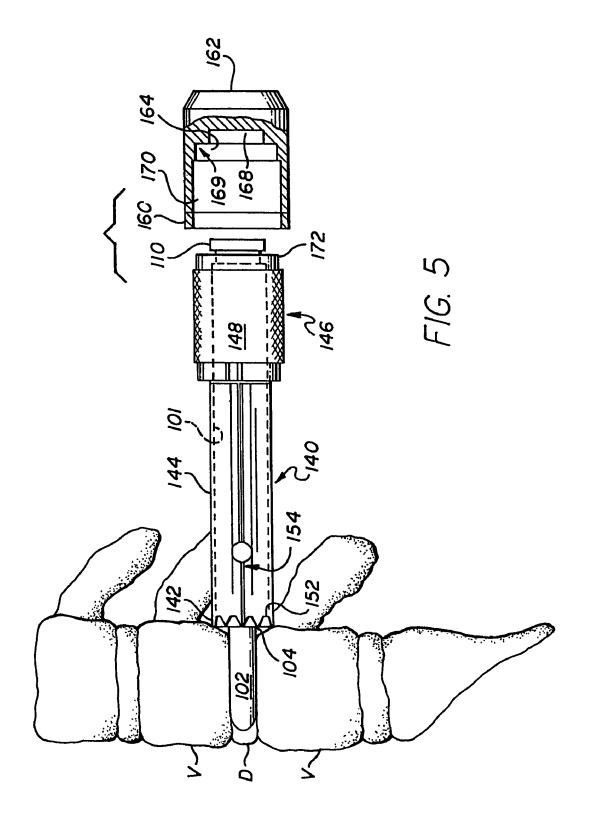
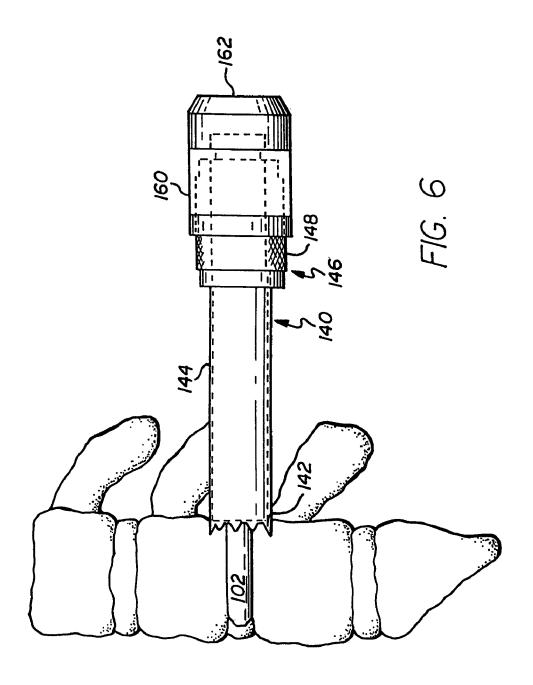



FIG. 3F


FIG. 4







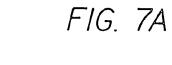

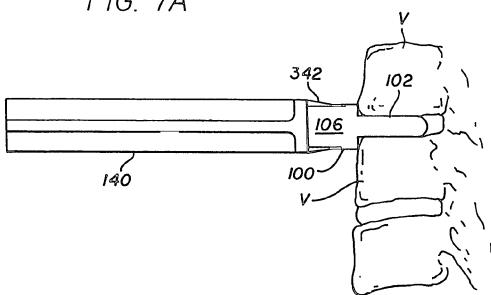

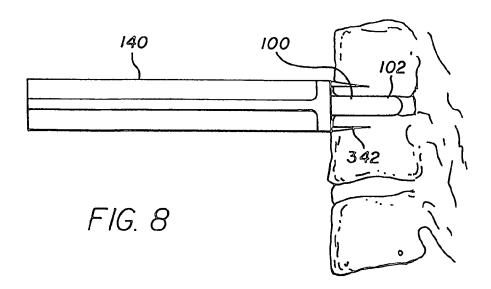


FIG. 7B

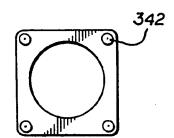


FIG. 7C

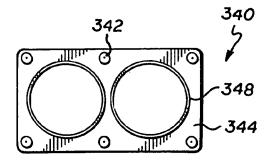


FIG. 7D

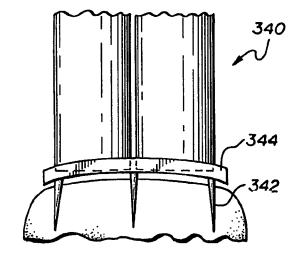
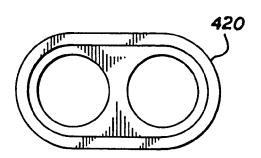
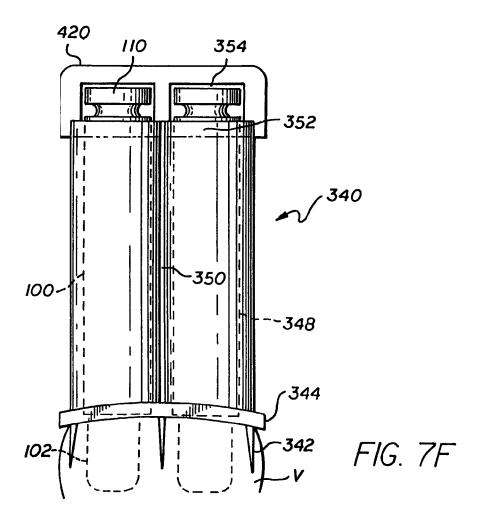
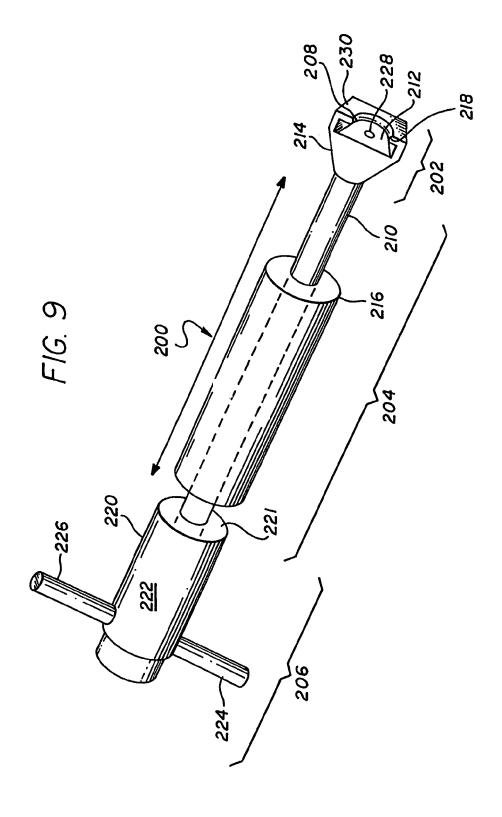
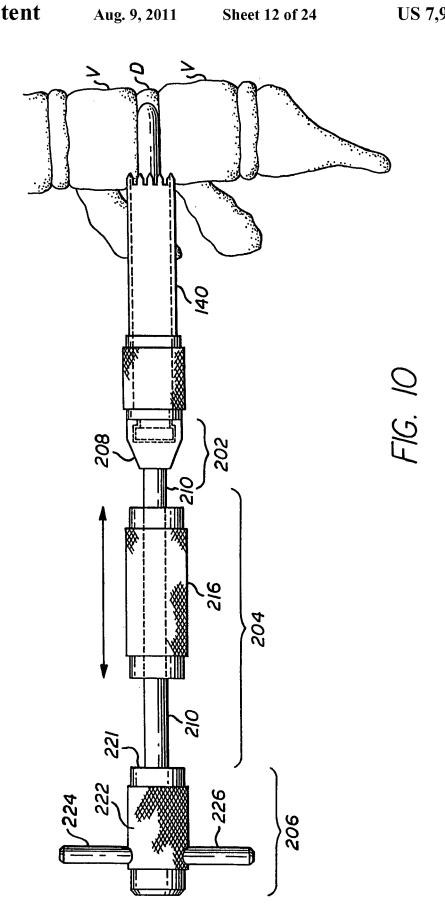



FIG. 7E









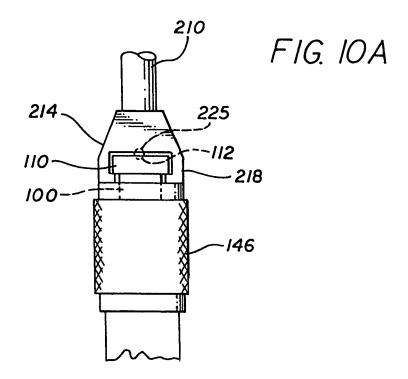

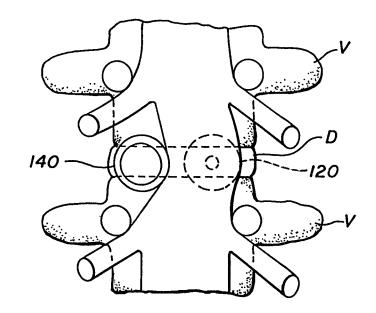


FIG. 10B

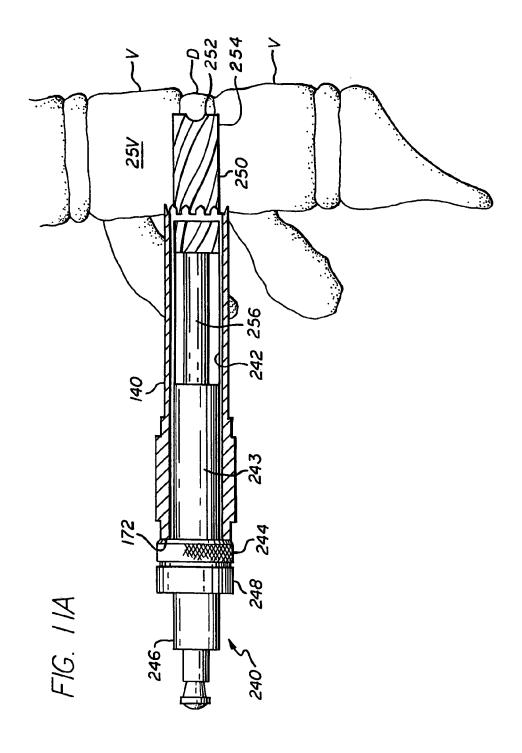
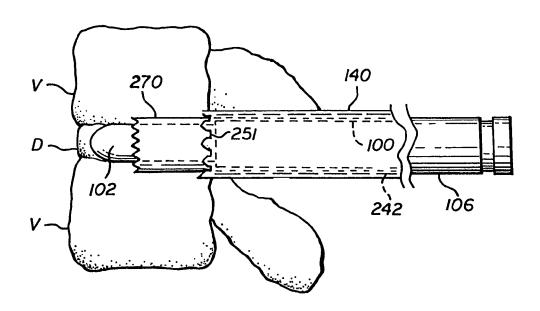
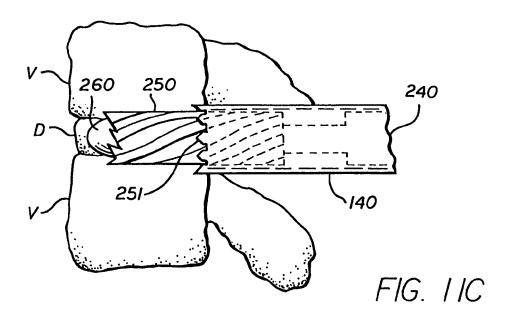
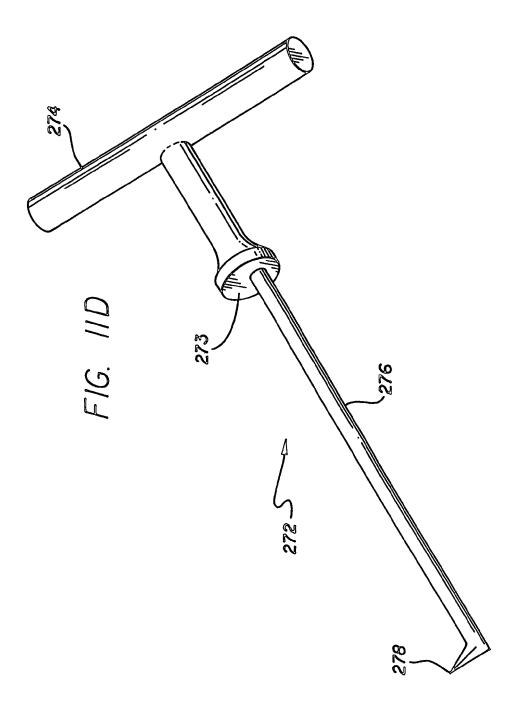
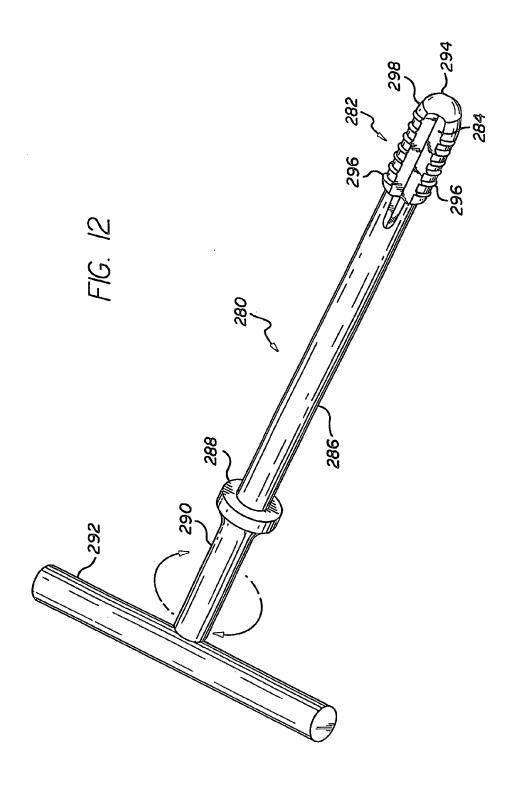


FIG. I IB









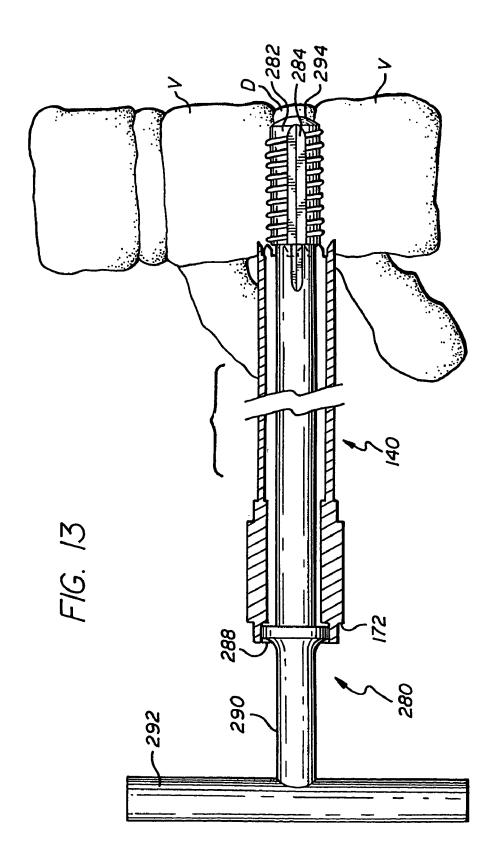
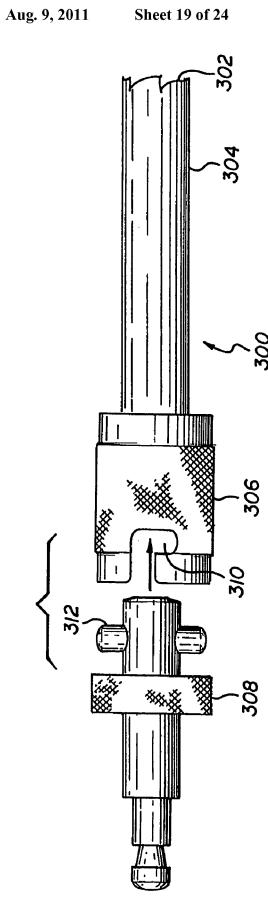
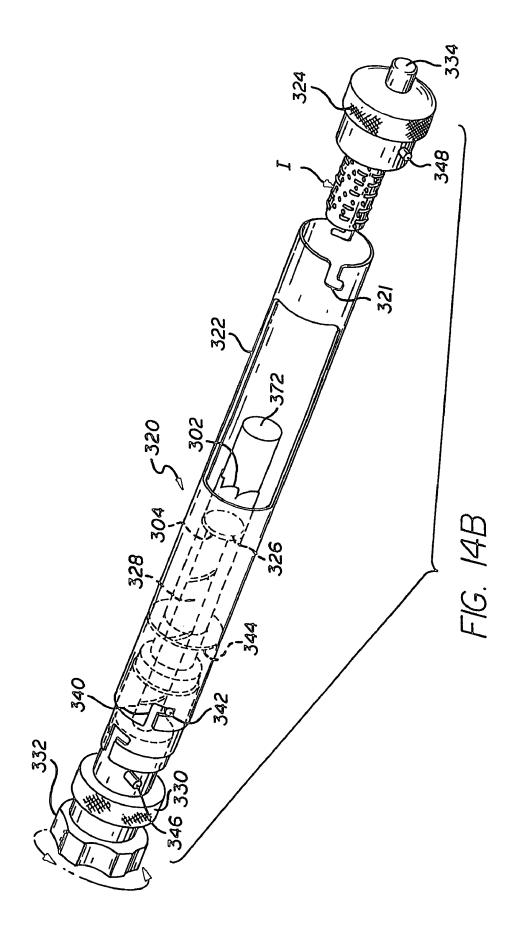
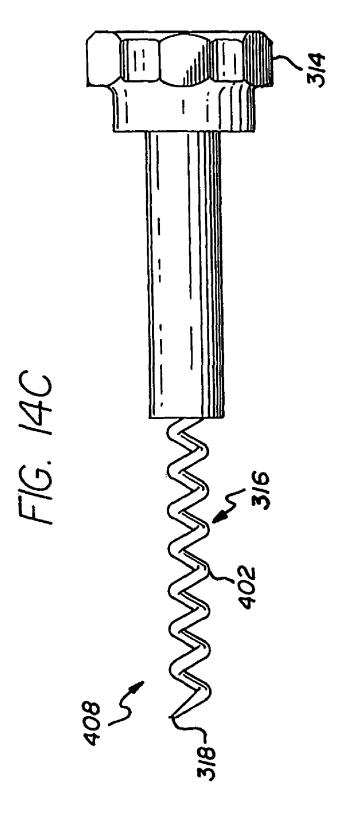
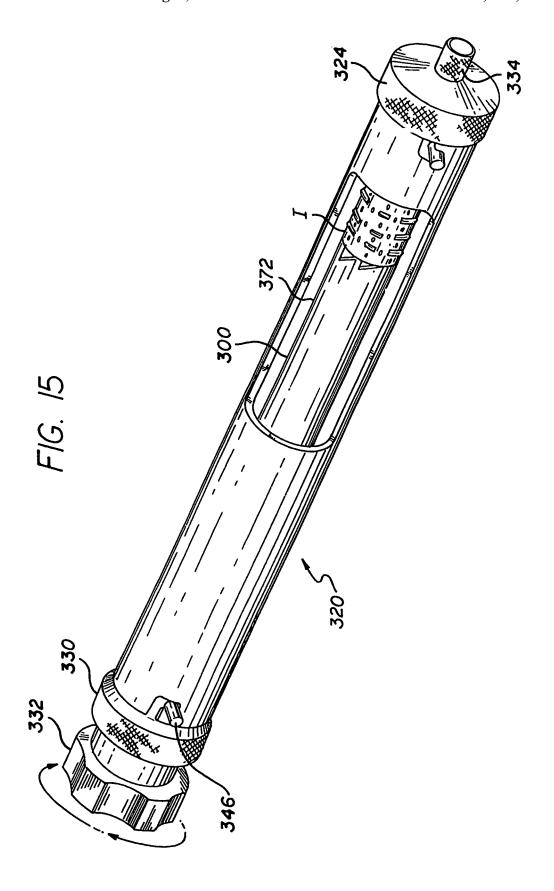



FIG. 14A









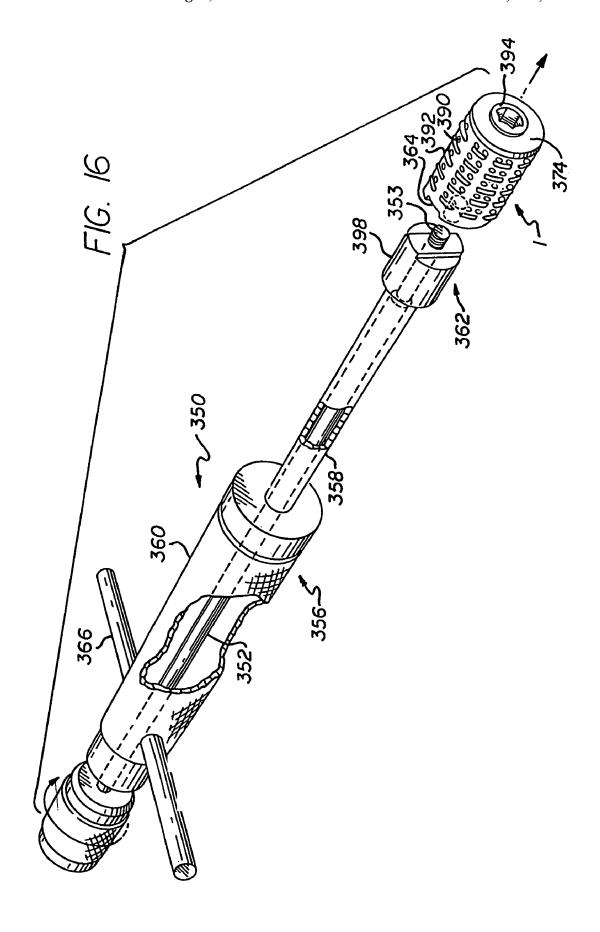


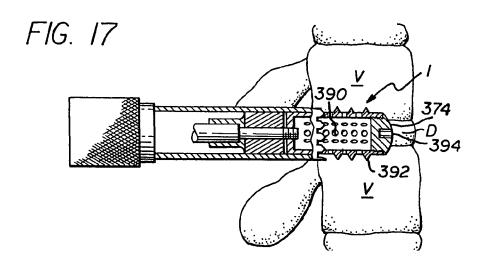


FIG. 18

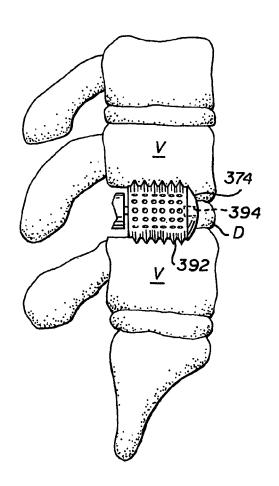

#### GUARD FOR USE IN PERFORMING HUMAN INTERBODY SPINAL SURGERY

#### RELATED APPLICATIONS

This is a continuation of application Ser. No. 08/484,927, filed Jun. 7, 1995, now U.S. Pat. No. 6,096,038 which is a divisional of application Ser. No. 08/074,781, filed Jun. 10, 1993, now U.S. Pat. No. 5,484,437, issued Jan. 16, 1996, all of which are incorporated herein by reference.

#### BACKGROUND OF THE INVENTION

#### 1. Field of the Invention

be placed into the intervertebral space left remaining after the removal of a damaged spinal disc and specifically to the apparatus for and method of, inserting the implants.

#### 2. Description of the Prior Art

For the purpose of achieving long term stability to a seg-20 ment of injured spine, a fusion (the joining together of two or more bones via a continuous bridge of incorporated bone) may be performed. Well-known to those skilled in such art is the interbody fusion wherein the disc is partially excised and bone placed within that space previously occupied by that 25 disc material (between adjacent vertebrae) for the purpose of restoring a more normal spatial relationship, and to provide for stability; short term by mechanical support, and long term by the permanent cross bonding of bone from vertebra to vertebra. For fusion to occur within the disc space, it is necessary to prepare the vertebrae to be fused by breaking through, or cutting into, the hardened outside plates of bone (the endplates) to allow the interposed bone graft to come into direct contact with the more vascular cancellous (spongy) bone, and to thereby trick the body into attempting to heal this 35 induced, but controlled, "fracturing" by both bone production and the healing of the grafts to both opposed vertebral surfaces such that they become one continuous segment of bone.

The purpose of the present invention is to provide an implant, and the apparatus and method of inserting the 40 implant within the intervertebral space left after the removal of the disc material and permanently eliminate all motion at that location. To do so, the device of the present invention is space occupying within the disc interspace, rigid, self-stabilizing to resist dislodgement, stabilizing to the adjacent spinal 45 vertebrae to eliminate local motion, and able to intrinsically participate in a vertebra to vertebra bony fusion so as to assure the permanency of the result.

At present, following the removal of a damaged disc, either bone or nothing is placed into the remaining space. Placing 50 nothing into this space allows the space to collapse which may result in damage to the nerves; or the space may fill with scar tissue and eventually lead to a reherniation. The use of bone to fill the space is less than optimal in that bone obtained from the patient requires additional surgery and is of limited 55 availability in its most useful form, and if obtained elsewhere, lacks living bone cells, carries a significant risk of infection, and is also limited in supply as it is usually obtained from accident victims. Furthermore, regardless of the source of the bone, it is only marginal structurally and lacks a means to 60 either stabilize itself against dislodgement, or to stabilize the adjacent vertebrae.

#### a. Prior Art Implants

There have been an extensive number of attempts to develop an acceptable disc prosthesis (an artificial disc). Such 65 devices by design would be used to replace a damaged disc and seek to restore the height of the interspace and to restore

2

the normal motion of that spinal joint. No such device has been found that is medically acceptable. This group of prosthetic or artificial disc replacements, seeking to preserve spinal motion and so are so different from the present invention, would include:

- U.S. Pat. No. 3,867,728 to STUBSTAD—describing a flexible disc implant.
- U.S. Pat. No. 4,349,921 to KUNTZ—describing a flexible disc replacement with file-like surface projections to discourage device dislocation.
- U.S. Pat. No. 4,309,777 to PATIL—describing a motion preserving implant with spiked outer surfaces to resist dislocation and containing a series of springs to urge the vertebrae away from each other.
- The present invention relates to artificial fusion implants to 15 U.S. Pat. No. 3,875,595 to FRONING—describing a motion preserving bladder-like disc replacement with two opposed stud-like projections to resist dislocation.
  - U.S. Pat. No. 2,372,622 to FASSIO (France)—describing a motion preserving implant comprising complimentary opposed convex and concave surfaces.

In summary, these devices resemble the present invention only in that they are placed within the intervertebral space following the removal of a damaged disc. In that they seek to preserve spinal motion, they are diametrically different from the present invention which seeks to permanently eliminate all motion at that spinal segment.

A second related area of prior art includes those devices utilized to replace essentially wholly removed vertebrae. Such removal is generally necessitated by extensive vertebral fractures, or tumors, and is not associated with the treatment of disc disease. While the present invention is to be placed within the disc space, these other vertebral devices cannot be placed within the disc space as at least one vertebra has already been removed such that there no longer remains a "disc space". Furthermore, these devices are limited in that they seek to perform as temporary structural members mechanically replacing the removed vertebrae (not a removed disc), and do not intrinsically participate in supplying osteogenic material to achieve cross vertebrae bony fusion. Therefore, unlike the present invention which provides for a source of osteogenesis, use of this group of devices must be accompanied by a further surgery consisting of a bone fusion procedure utilizing conventional technique. This group consisting of vertebral struts rather than disc replacements would include the following:

- U.S. Pat. No. 4,553,273 to WU—describing a turnbuckle-like vertebral strut.
- U.S. Pat. No. 4,401,112 to REZAIAN—describing a turnbuckle-like vertebral strut with the addition of a long stabilizing staple that spans the missing vertebral body.
- U.S. Pat. No. 4,554,914 to KAPP—describing a large distractible spike that elongates with a screw mechanism to span the gap left by the removal of an entire vertebra and to serve as an anchor for acrylic cement which is then used to replace the missing bone (vertebrae).
- U.S. Pat. No. 4,636,217 to OGILVIE—describing a vertebral strut mechanism that can be implanted after at least one vertebrae has been removed and consists of a mechanism for causing the engagement of screws into the vertebrae above and the vertebrae below the one removed.

In summary, this second group of devices differs from the present invention in that they are vertebral replacements struts, do not intrinsically participate in the bony fusion, can only be inserted in the limited circumstances where an entire vertebra has been removed from the anterior approach, and are not designed for, or intended to be used for the treatment of disc disease.

A third area of prior art related to the present invention includes all devices designed to be applied to one of the surfaces of the spine. Such devices include all types of plates, struts, and rods which are attached by hooks, wires and screws. These devices differ significantly from the present invention in that they are not inserted within the disc space and furthermore do not intrinsically participate in supplying osteogenic material for the fusion.

Therefore, where permanent spinal immobilization is desired, an additional surgery, consisting of a spinal fusion performed by conventional means or the use of supplemental methylmethacrylate cement is required. Such devices applied to the spine, but not within the disc space, would include the following:

- U.S. Pat. No. 4,604,995 to STEPHENS—describing a "U" shaped metal rod attached to the posterior elements of the spine with wires to stabilize the spine over a large number of segments.
- U.S. Pat. No. 2,677,369 to KNOWLES—describing a metal 20 column device to be placed posteriorly along the lumbar spine to be held in position by its shape alone and to block pressure across the posterior portions of the spinal column by locking the spine in full flexion thereby shifting the maximum weight back onto the patient's own disc.

Other devices are simply variations on the use of rods (e.g. Harrington, Luque, Cotrel-Dubosset, Zielke), wires or cables (Dwyer), plates and screws (Steffee), or struts (Dunn, Knowles)

In summary, none of these devices are designed to be nor 30 can be used within the disc space. Moreover, these devices do not replace a damaged disc, and do not intrinsically participate in the generation of a bony fusion.

Another area of related prior art to be considered is that of devices designed to be placed within the vertebral interspace 35 following the removal of a damaged disc, and seeking to eliminate further motion at that location.

Such a device is contained in U.S. Pat. No. 4,501,269 issued to BAGBY which describes an implantable device and limited instrumentation. The method employed is as follows: 40 a hole is bored transversely across the joint and a hollow metal basket of larger diameter than the hole is then pounded into the hole and then the hollow metal basket is filled with the bone debris generated by the drilling.

While the present invention (device, instrumentation, and 45 method) may appear to bear some superficial resemblance to the BAGBY invention, it is minimal, while the differences are many fold and highly significant. These differences include the following:

1. Safety—The present invention provides for a system of 50 completely guarded instrumentation so that all contiguous vital structures (e.g. large blood vessels, neural structures) are absolutely protected. The instrumentation of the present invention also makes overpenetration by the drill impossible. Such overpenetration in the cervical spine, for example, 55 would result in the total paralysis or death of the patient. In the thoracic spine, the result would be complete paraplegia. In the lumbar spine, the result would be paraplegia or a life-threatening perforation of the aorta, vena cave, or iliac vessels.

The present invention is atraumatically screwed into place 60 while the BAGBY device, in contradistinction, is pounded into position. BAGBY describes that its implant is significantly larger in size than the hole drilled and must be pounded in. This is extremely dangerous and the pounding occurs directly over the spinal cord which is precariously vulnerable 65 to percussive injury. Furthermore, while it is possible, for example in the lumbar spine, to insert the present invention

4

away from the spinal cord and nerves, the BAGBY device must always be pounded directly towards the spinal cord.

Furthermore, since the BAGBY device is pounded into a smooth hole under great resistance, and lacking any specific design features to secure it, the device is highly susceptible to forceful ejection which would result in great danger to the patient and clinical failure. The present invention, in contradistinction, is securely screwed into place, and possesses highly specialized locking threads to make accidental dislodgement impossible. Because of the proximity of the spinal cord, spinal nerves, and blood vessels, any implant dislodgement as might occur with the BAGBY device might have catastrophic consequences.

- 2. Broad applicability—The BAGBY device can only be inserted from the front of the vertebral column, however, in contrast, the present invention can be utilized in the cervical, thoracic, and lumbar spine, and can be inserted from behind (posteriorly) in the lumbar spine. This is of great importance in that the purpose of these devices is in the treatment of disc disease and probably greater than 99 percent of all lumbar operations for the treatment of disc disease are performed from behind where the present invention can easily be utilized, but the BAGBY device, as per BAGBY'S description, cannot.
- 3. Disc removal—The BAGBY invention requires the complete removal of the disc prior to the drilling step, whereas the present invention eliminates the laborious separate process of disc removal and efficiently removes the disc and prepares the vertebral end plates in a single step.
- 4. Time required—The present invention saves time over the BAGBY invention since time is not wasted laboring to remove the disc prior to initiating the fusion. Also, with the present invention the procedure is performed through a system of guarded instrumentation, time is not wasted constantly placing and replacing various soft tissue retractors throughout the procedure.
- 5. Implant stability—Dislodgement of the implant would be a major source of device failure (an unsuccessful clinical result), and might result in patient paralysis or even death. As discussed, the BAGBY device lacks any specific means of achieving stability and since it is pounded in against resistance to achieve vertebral distraction, and is susceptible to forceful dislodgement by the tendency of the two distracted vertebrae, to return to their original positions squeezing out the device. The present invention, however, is screwed into place. As there is no unscrewing force present between the vertebrae, compression alone cannot dislodge the implant. The implant is inherently stable by its design. Furthermore, the threads of the present invention are highly specialized in that they are periodically interrupted so that the tail ends of each of the tabs so formed are blunted and twisted so as to resist accidental unscrewing. The removal of an implant with such "locking threads" requires the use of a special extractor included within the instrumentation. The stability of the present invention is still further enhanced, again in contradistinction to the BAGBY device, by the presence of a "bone ingrowth" surface texturing, which both increases the friction of the fit and allows for the direct growth of the vertebral bone into the casing of the implant itself.
- 6. Spinal stability—The present invention is not only self-stabilizing, it also provides stability to the adjacent vertebrae in at least three ways that the BAGBY device cannot. First, the BAGBY device is placed transversely across the joint in the center, leaving both vertebrae free to rock back and forth over this round barrel shaped axis, much like a board over a barrel, being used for a seesaw.

Secondly, as the BAGBY device lacks any specific design features to resist sliding, it may actually behave as a third body allowing the translation of the vertebrae relative to the device and to each other.

Thirdly, any device can only provide stability if it remains 5 properly, seated. The present invention is inherently stable, and therefore assures that it will stabilize the adjacent vertebrae, rather than, as with the BAGBY, the instability of the spine to be treated may cause a dislocation of the BAGBY implant, with further loss of spinal stability.

7. The collapse of the interspace—While both the present invention and the BAGBY device can be fabricated to withstand the compression forces within the interspace, the interspace may nevertheless collapse under the superincumbent body weight as the implant settles into the vertebral bone. 15 This is related to the load per unit area. Again the present invention is superior to the BAGBY device in at least four ways.

First, the present invention offers considerably greater surface area to distribute the load. Secondly, while the BAGBY 20 device is placed centrally, the present device is placed bilaterally where the bone tends to be more cortical and much stronger out towards the rim. Thirdly, the present invention supports the load achieving an "I" beam effect, whereas the BAGBY implant does not. Fourthly, it is not pressure alone 25 that causes the collapse of the bone adjacent to the implant, but also bony erosion that is caused by the motion under pressure of the implant against the bone. As discussed in item 6 above, the present invention alone is highly resistant to such motion, again diminishing the likelihood of erosion and interspace collapse.

8. Bone ingrowth surface texturing—The present invention has a surface treatment of known and conventional technology to induce the growth of bone from the vertebrae directly into the casing material of the implant itself. The BAGBY 35 device has no similar feature. {L.A.—we may want to list examples of these bone growth factors}

9. Fusion mass—The BAGBY invention calls for removing the disc and then drilling a hole between the adjacent vertebrae. The bony debris so generated is then put into the 40 device. The present invention takes a core of pure bone producing marrow from the iliac crest, and then by use of a special press, forcibly injects the implant device with an extremely dense compressed core of that osteogenic material until the material itself virtually extrudes from every cell of 45 the implant.

10. The probability of achieving fusion—The fusion rate within the spine is known to be related directly to the amount of exposed vascular bone bed area, the quality and quantity of the fusion mass available, and the extent of the stabilization 50 obtained with all other factors being half constant. It would then be anticipated, that the fusion rate would be superior with the present invention as compared to the BAGBY device, because of optimal implant stability (#5), optimal spinal stability (#6), bone ingrowth surface treatment (#8), superior 55 fusion mass (#9), and the greater exposed vertebral bony surface area (#7).

The last area of prior art possibly related to the present invention and therefore, to be considered related to "bony ingrowth", are patents that either describe methods of producing materials and or materials or devices to achieve the same. Such patents would include:

U.S. Pat. No. 4,636,526 (DORMAN), U.S. Pat. No. 4,634, 720 (DORMAN), U.S. Pat. No. 4,542,539 (ROWE), U.S. Pat. No. 4,405,319 (COSENTINO), U.S. Pat. No. 4,439,152 65 (SMALL), U.S. Pat. No. 4,168,326 (BROEMER), U.S. Pat. No. 4,535,485 (ASHMAN), U.S. Pat. No. 3,987,499

6

(SCHARBACH), U.S. Pat. No. 3,605,123 (HAHN), U.S. Pat. No. 4,655,777 (DUNN), U.S. Pat. No. 4,645,503 (LIN), U.S. Pat. No. 4,547,390 (ASHMAN), U.S. Pat. No. 4,608,052 (VAN KAMPEN), U.S. Pat. No. 4,698,375 (DORMAN), U.S. Pat. No. 4,661,536 (DORMAN), U.S. Pat. No. 3,952, 334 (BOKROS), U.S. Pat. No. 3,905,047 (LONG), U.S. Pat. No. 4,693,721 (DUCHEYNE), U.S. Pat. No. 4,070,514 (ENTHERLY).

However, while the implant of the present invention would 10 utilize bone ingrowth technology, it would do so with conventional technology.

b. Prior Art Instrumentations and Methods

The following is a history of the prior art apparatus and methods of inserting spinal implants:

In 1956, Ralph Cloward developed a method and instruments which he later described for preparing the anterior aspect (front) of the cervical spine, and then fusing it. Cloward surgically removed the disc to be fused across and then placed a rigid drill guide with a large foot plate and prongs down over an aligner rod and embedded said prongs into the adjacent vertebrae to maintain the alignment so as to facilitate the reaming out of the bone adjacent the disc spaces. As the large foot plate sat against the front of the spine, it also served as a fixed reference point to control the depth of drilling. The reaming left two opposed resected arcs, one each, from the opposed vertebral surfaces. The tubular drill guide, which was placed only preliminary to the drilling, was thereafter completely removed. A cylindrical bony dowel, significantly larger in diameter than the hole formed, was then pounded into the hole already drilled. Cloward's method of instrumentation was designed for, and limited to, use on the anterior aspect and in the region of the cervical spine only. The hole was midline, which would preclude its use posteriorly where the spinal cord would be in the way.

As the bone graft to be inserted in Cloward's method was necessarily larger in diameter than the hole drilled, the graft could not be inserted through the drill guide. This mandated the removal of the drill guide and left the graft insertion phase completely unprotected. Thus Cloward's method and instrumentation was inappropriate for posterior application.

In addition, the failure to provide continuous protection to the delicate neural structures from the instruments, as well as the bony and cartilaginous debris generated during the procedure, made Cloward's method inappropriate for posterior application. Also, the drill guide described by Cloward could not be placed posteriorly within the spinal canal, as the foot plate would crush the nerves. Modifying Cloward's drill guide by removing the foot plate completely, would still leave the instrument unworkable as it would then lack stability, and would not be controllable for depth of seating.

Nevertheless, Wilterberger, (Wilterberger, B. R., Abbott, K. H., "Dowel Intervertebral Fusion as Used in Lumbar Disc Surgery," *The Journal of Bone and Joint Surgery*, Volume 39A, pg. 234-292, 1957) described the unprotected drilling of a hole from the posterior into the lumbar spine between the nerve roots and across the disc space, and then inserting a stack of button-like dowels into that space. While Wilterberger had taken the Cloward concept of circular drilling and dowel fusion and applied it to the lumbar spine from a posterior approach, he had not provided for an improved method, nor had he advanced the instrumentation so as to make that procedure sufficiently safe, and it rapidly fell into disrepute.

Crock (Crock, H. V., "Anterior Lumbar Interbody Fusion—Indications for its Use and notes on Surgical Technique," *Clinical Orthopedics*, Volume 165, pg. 157-163, 1981) described his technique and instrumentation for Anterior Interbody Fusion of the lumbar spine, wherein he drilled

two large holes side by side across the disc space from anterior to posterior essentially unprotected and then pounded in two at least partially cylindrical grafts larger than the holes

A review of the prior art is instructive as to a number of 5 significant deficiencies in regard to the method and instrumentation for the performance of Interbody Spinal Fusion utilizing drilling to prepare the endplates.

As the great majority of spinal surgery is performed in the lumbar spine and from posteriorly, a review of the prior art 10 reveals a number of deficiencies in regard to the spine in general, and to the posterior approach to the lumbar spine specifically. These deficiencies include the:

- 1. Failure to protect the surrounding tissues throughout the procedure, specifically, prior to drilling and until after the 15 fixation. insertion of the graft;
- 2. Failure to contain the debris, bony and cartilaginous, generated during the procedure;
- 3. Failure to optimize the contact of the cylindrical drill hole and bone graft, the mismatch in their diameters resulting 20 the following advantages: in incongruence of fit;
- 4. Failure to determine the optimal drill size prior to drill-
- 5. Failure to determine the optimal amount of distraction prior to drilling;
- 6. Inability to optimize the amount of distraction so as to restore the normal spatial relationships between adjacent ver-
- 7. Inability to create sufficient working space within the spinal canal (between the nerve roots and the dural sac) to 30 make the procedure safe;
- 8. Absent a foot plate on the drill guide, as necessitated by the close tolerances posteriorly, the inability to reliably insure that the drilling is parallel to the vertebral endplates;
- opposed vertebral surfaces; and
- 10. The inability to determine within the spinal canal, the proper side by side positioning for dual drill holes.

#### BRIEF SUMMARY OF THE INVENTION

The present invention comprises a series of artificial implants, the purpose of which is to participate in, and directly cause bone fusion across an intervertebral space following the excision of a damaged disc. Such implants are 45 structurally load bearing devices, stronger than bone, capable of withstanding the substantial forces generated within the spinal interspace. The devices of the present invention have a plurality of macro sized cells and openings, which can be loaded with fusion promoting materials, such as autogenous 50 bone, for the purpose of materially influencing the adjacent vertebrae to perform a bony bond to the implants and to each other. The implant casing may be surface textured or otherwise treated by any of a number of known technologies to achieve a "bone ingrowth surface" to further enhance the 55 stability of the implant and to expedite the fusion.

The devices of the present invention are configured and designed so as to promote their own stability within the vertebral interspace and to resist being dislodged, and furthermore, to stabilize the adjacent spinal segments.

The apparatus and method of the present invention for preparing the vertebrae for insertion of the implant allows for the rapid and safe removal of the disc, preparation of the vertebrae, performance of the fusion, and internal stabilization of the spinal segment.

The present invention is a method for Interbody Spinal Fusion utilizing novel instrumentation, whereby a protective

tubular member is placed prior to the drilling part of the procedure and is left in place until the graft is fully seated.

In the preferred embodiment two distractors are used to separate two adjacent vertebrae to a preferred distance. A hollow Outer Sleeve having teeth at one end is driven into the adjacent vertebrae on one side to hold the vertebrae in position when the distractor is removed, a diameter reducing hollow Inner Sleeve is introduced into the Outer Sleeve, a drill having a drill stop is passed through the hollow Inner Sleeve to drill a hole to a desired depth, and an implant is inserted in the hole. The method is repeated on the other side of the disc.

In summary then, the present invention, instrumentation, and method, provides for a single surgery providing for an integrated discectomy, fusion, and interbody internal spinal

#### Discussion of the Instrumentation

The apparatus and method of the present invention provide

- 1. The present invention is safer by providing protection of the surrounding tissues. An Outer Sleeve places all of the delicate soft tissue structures, nerves, blood vessels, and organs outside of the path of the various sharp surgical instruments and the implant. Further, it is an improvement upon hand held retractors in that it occupies the least possible amount of area, avoids the stretching associated with manual retraction, provides for the retraction and shielding of the surrounding tissues in all directions circumferentially and simultaneously, and it does so exclusively with smooth, curved surfaces.
- 2. The present invention is safer by providing protection against the danger of instrument or implant overpenetration.
- 3. The present invention is safer as the surgical site and 9. The inability to insure equal bone removal from the 35 wound are protected from the debris generated during the procedure.
  - 4. The present invention is safer because the method provides for absolute protection to the soft tissues directly and from indirect injury by overpenetration. It makes safe the use 40 of power instrumentation which is both more effective and efficient.
    - 5. The present invention maintains the vertebrae to be fused rigid throughout the procedure.
    - 6. The present invention holds the vertebrae to be fused aligned throughout the procedure.
    - 7. The present invention holds the vertebrae to be fused distracted throughout the procedure.
    - 8. The present invention assures that all instruments introduced through the Outer Sleeve are coaxial and equally centered through the disc space and parallel the endplates.
    - 9. The present invention facilitates the implant insertion by countering the high compressive forces tending to collapse the interspace, which if left unchecked would resist the introduction and advancement of the implant and make stripping more likely.
    - 10. The present invention extends the range and use of the procedure and similarly the interbody spinal implant itself by making the procedure safe throughout the spine.
    - 11. The present invention increases the ability to use a specifically sized implant.
    - 12. In the present invention the end of all the penetrating instrumentation is blunt faced.
    - 13. In the present invention all of the instruments have been stopped at a predetermined depth to avoid overpenetration.
    - 14. The design of the Outer Sleeve in the present invention conforms to the spacial limitations of the specific surgical

- 15. The design and use of a second or Inner Sleeve in the present invention allows for the difference in size between the inside diameter of the Outer Sleeve, and the outside diameter of the drill itself. This difference being necessary to accommodate the sum of the distraction to be produced, and the depth of the circumferential threading present of the implant.
- 16. In the present invention a specially designed drill bit with a central shaft recess allows for the safe collection of the drilling products, which can then be removed without disturbing the Outer Sleeve by removing the drill bit and Inner Sleeve as a single unit.
- 17. In the present invention a specially designed trephine for removing a core of bone slightly smaller in diameter than the internal diameter of the implant cavity itself, however of 15 a greater length.
- 18. In the present invention a specially designed press for forcefully compressing and injecting the long core of autogenous bone into the implant, such that it extrudes through the implant itself.
- 19. In the present invention a specially designed driver extractor, which attaches to the implant and allows the implant to be either inserted or removed without itself dissociating from the implant, except by the deliberate disengagement of the operator.
- 20. In the present invention predistraction increases the working space.
- 21. The Distractor in the present invention is self-orienting acting as a directional finder.
- 22. The Distractor in the present invention is self-centralizing between the opposed vertebral surfaces acting as a centering post for the subsequent bone removal.
- 23. In the present invention predistraction assures the equal removal of bone from the adjacent vertebral surfaces.
- 24. In the present invention predistraction assures the exact 35 congruence between the hole drilled and the device
- 25. In the present invention predistraction assures that the drilling is parallel to the vertebral endplates.
- 26. In the present invention predistraction allows for the determination of the optimal distraction prior to drilling.
- 27. In the present invention predistraction allows for the verification of the correct prosthesis size prior to drilling.
- 28. In the present invention predistraction facilitates device insertion by relieving the compressive loads across the interspace which would resist implantation.
- 29. In the present invention predistraction decreases the likelihood of stripping the bone during insertion.
- 30. In the present invention predistraction provides for the side by side positioning, spacing, and parallelism required prior to the irrevocable event of drilling.
- 31. In the present invention predistraction provides for the rigid stabilization of the vertebrae opposed to the disc space throughout the surgical procedure.
- 32. In the present invention predistraction provides for an implant easier to insert as the compressive loads of the 55 opposed vertebrae are held in check so that the device itself need not drive the vertebrae apart to be inserted.
- 33. In the present invention predistraction allows for the insertion of a more effective implant as more of the implant can be dedicated to its intended purpose and be full diameter, 60 whereas without the benefit of predistraction and the ability to maintain the same, a significant portion of the forward end of the implant would need to be dedicated to the purpose of separating the opposing vertebrae.
- 34. The present invention allows for the use of an implant 65 with a sharper thread or surface projections as there is no danger to the surrounding tissues.

10

35. The present invention allows for the implant to be fully preloaded as provided to the surgeon, or for the surgeon to load it with the material of his choice at the time of surgery.

36. The present invention allows for the loading of a spinal implant outside of the spinal canal and prior to implantation.

#### Objects of the Present Invention

It is an object of the present invention to provide an improved method of performing a discectomy, a fusion, and an internal stabilization of the spine, and specifically, all three of the above simultaneously and as a single procedure.

It is another object of the present invention to provide an improved method of performing a discectomy, a fusion, and an internal stabilization of the spine, which is both quicker and safer than is possible by previous methods.

It is another object of the present invention to provide an improved method of performing a discectomy, a fusion and an 20 internal stabilization of the spine, to provide for improved surgical spinal implants.

It is another object of the present invention to provide an improved method of performing a discectomy, a fusion, and an internal stabilization of the spine, which provides for an improved system of surgical instrumentation to facilitate the performance of the combined discectomy, fusion, and internal spinal stabilization.

It is another object of the present invention to provide an improved method of performing a discectomy, a fusion, and an internal stabilization of the spine procedures.

It is an object of the present invention to provide instrumentation and a method of spinal interbody arthrodesis that is faster, safer, and more efficacious than prior methods, and can effectively be performed in the cervical, thoracic, and lumbar spine anteriorly, as well as in the lower lumbar spine poste-

It is a further object of the present invention to provide a means for inserting a spinal implant between adjacent vertebrae while maintaining their optimal spacing, positioning, and alignment.

These and other objects of the present invention will be apparent from review of the following specification and the accompanying drawings.

#### BRIEF DESCRIPTION OF THE DRAWINGS

- FIG. 1 is a side view of the Long Distractor, of the present invention inserted into the intervertebral space.
- FIG. 2 is a side view of a Convertible Distractor assembly 50 in relation to the spine.
  - FIG. 3 is a perspective view of a high retention Short Distractor of FIG. 2.
  - FIG. 3A is a side view of the high retention Short Distractor of FIG. 2.
  - FIG. 3B is a side view of an alternative Short Distractor with circumferential forward facing ratcheting.
  - FIG. 3C is a top view of the alternative Short Distractor of FIG. 3B.
  - FIG. 3D is a perspective view of an alternative embodiment of a Short Distractor.
  - FIG. 3E is a top view of the alternative distractor of FIG.
  - FIG. 3F is a side view of a further alternative rectangularized Short Distractor with knurled surfaces.
  - FIG. 4 is a perspective view of a spinal segment (two vertebrae and an interposed disc) with a Short Distractor in place, with a portion of the upper vertebrae and disc cut away

11

to show the Short Distractor on one side of the spine and the Long Distractor about to be placed contralaterally.

FIG. 5 shows a side view of the Outer Sleeve in place over the Long Distractor, and about to receive the Driver Cap in preparation for being seated.

FIG. 6 shows the Long Distractor, Outer Sleeve, and Driver Cap following the proper seating of the Outer Sleeve into the two adjacent vertebrae.

FIG. 7A is a side view of the cervical Outer Sleeve being placed over a Long Distractor which is in place within the disc space anteriorly.

FIG. 7B is a bottom view of the single Outer Sleeve of FIG. 7A.

FIG. 7C is a bottom view of a Dual Outer Sleeve.

FIG. 7D is an enlarged side view of the proximal portion of FIG. 7C.

FIG. 7E is a bottom view of a Dual Driver Cap for driving two distractors.

FIG. 7F is a side sectional view showing the Dual Outer 20 Sleeve of FIGS. 7C and 7D, Distractors and Dual Cap of FIG. 7E seated.

FIG. 8 is a side view of the Outer Sleeve of FIG. 7A centered on the Long Distractor and fully seated on the anterior aspect of the cervical spine.

FIG. 9 is a perspective view of the Distractor Puller.

FIG. 10 is a cutaway partial side view of the Proximal Puller engaging the extraction ring of the Long Distractor over the end of the Outer Sleeve.

FIG. 10A is a side view of the Puller coupled to the Long Distractor just prior to its extraction.

FIG. 10B is a posterior view of the proximal Outer Sleeve and a Short Distractor in place in regard to the vertebrae, disc and nerves

FIG. 11A is a side sectional view of the Drill and Inner Sleeve within the Outer Sleeve and drilling across the intervertebral space and cutting partially cylindrical arcs from the adjacent vertebrae.

FIG. 11B is a sectional side view of preparation of the 40 intervertebral space by the alternative "Trephine Method" showing the Distractor, Trephine, Inner Sleeve, and Outer Sleeve in place.

FIG. 11C is a sectional side view as in FIG. 11A, but showing the use of an alternative drilling conformation 45 wherein the extended proximal portion is both distracting and self-centering.

FIG. 11D is a side view of an instrument for removing arcs of bone from vertebrae following drilling.

FIG. 12 is a perspective view of the surgical Tap.

FIG. 13 is a side view of the Outer Sleeve and the surgical Tap fully threaded within the interspace.

FIG. **14**A is a side view of the bone harvesting Trephine and motor adapter.

FIG. 14B is a perspective view of the implant Bone Load- 55 ing Device.

FIG. **14**C is a perspective view of the Corkscrew bone freeing and extracting instrument.

FIG. 15 is a partial perspective view of the Bone Loading Device in operation.

FIG. 16 is a perspective view of the Implant Driver about to engage the spinal implant.

FIG. 17 is a side view of the spinal implant being fully seated within the intervertebral space by means of the Driver apparatus in place within the Outer Sleeve.

FIG. 18 is a side view of the lumbar spine showing the end result of the device implantation via the posterior route.

12

#### DETAILED DESCRIPTION OF THE DRAWINGS AND DETAILED DESCRIPTION OF METHOD OF INSERTION

The following discussion will be in regard to application in the lumbar spine via the posterior approach. In its simplest form, the method of the present invention involves the following steps. The patient is placed on a spinal surgery frame, which allows for the distraction and alignment of the disc space to be fused. A bilateral posterior exposure of the interspace, with or without partial discectomy is then performed. Utilizing distractors the disc space is distracted, and a hollow Outer Sleeve is fitted over one of the distractors. The end of the Outer Sleeve has teeth for engaging the two adjacent vertebrae. The Outer Sleeve is driven into the vertebrae and the distractor is then removed. A hollow Inner Sleeve is then inserted into the Outer Sleeve and a stopped Drill is utilized to prepare the opposed vertebral surfaces. The Drill and the Inner Sleeve are removed as a single unit. The space is tapped if so required. The prepared spinal implant is then inserted via the Outer Sleeve utilizing a stopped inserter. The instruments are then removed and the procedure repeated on the contralateral side of the spine.

Detailed Description of the Preferred Embodiment

Step 1a. Prior to surgery, translucent implant templates appropriately adjusted for scale are superimposed on AP, lateral, and axial images of the interspace to be fused, for the purpose of selecting the optimal implant size and to determine the desired distraction.

Step 1b. The patient is preferably placed onto a spinal surgery frame capable of inducing both distraction and vertebral alignment.

Step 2. In the preferred embodiment, a standard bilateral (partial) discectomy is performed and any posterior lipping of the vertebral bodies adjacent the interspace is removed. Alternatively, no disc material need be removed. In the preferred embodiment, the interspace is exposed by performing bilateral paired semihemilaminotomies and resecting the inner aspects of the facet joints adjacent the spinal canal while preserving the supra and interspinous, ligaments.

Step 3. Beginning on the first side, the dural sac and traversing nerve root at that level are retracted medially and a Long Distractor then inserted and impacted flush to the posterior vertebral bodies adjacent that interspace. Long Distractors with working ends of increasing diameter are then sequentially inserted until the optimal distraction is obtained. This optimal distraction not only restores the normal height of the interspace, but further achieves a balance wherein the tendency for the space to collapse is resisted, which in urging the vertebral bodies apart is being equally resisted by the powerful soft tissue structures about the spinal segment including the outer casing of the disc (the annulus fibrosus), various ligaments, capsular structures, as well as the muscles and other soft tissue structures. This balanced distraction not only provides for the spatial restoration of the height of the interspace, but for considerable stability as the space now resists further distraction or collapse.

In the preferred embodiment, as the desired distraction is approached, the use of the solid bodied Long Distractors is terminated and a disassemblable Convertible Distractor is placed with tactile and/or radiographic confirmation of ideal distraction. The Convertible Distractor is then disassembled such that the Short Distractor portion is left in place and the ultra-low profile head portion being positioned adjacent to the canal floor and safely away from the neural structures. To

insure that the Short Distractor remains in place until its removal is desired, various embodiments of the Short Distractor are available with varying degrees of resistance to dislodgment. In the preferred embodiment of the procedure, attention is then directed to the contralateral side of the spine. 5

Step 4. On the contralateral side of the same interspace the Long Distractor having at its working end the diameter Matching the Short Distractor already in place, is then inserted. If however, due to an asymmetrical collapse of the interspace it is then determined that greater distraction is required on the second side to achieve the optimal stability, then the appropriate Short Distractor would be placed on the second side. Then the Short Distractor would be removed from the first side and replaced with a larger Long Distractor so as to bring the interspace into balance.

In an alternative embodiment, the entire procedure is performed on the one side of the spine utilizing only the Long Distractor prior to repeating the procedure on the contralateral side of the spine. While this method can be performed in ferred embodiment, when utilized it is best performed using a Trephine which allows the Long Distractor to remain in place, thereby allowing for interspace distraction otherwise provided in the first method by the Short Distractor. This alternative method then requires the use of a Trephine over the 25 Long Distractor in lieu of a reamer and is therefore called the "Trephine Method", which will be discussed in detail later.

Step 5. With the Short Distractor in place on the first side of the spine, and the matching Long Distractor in place on the second side of the spine, and with the dural sac and traversing 30 nerve root safely retracted, the Outer Sleeve is placed over the Long Distractor and firmly impacted to its optimal depth using the Impaction Cap and a mallet. The Long Distractor is then removed.

Step 6. An Inner Sleeve is then placed within the Outer 35 Sleeve, and the interspace is then prepared on that side by utilizing a Drill, Endmill, Reamer, or Trephine to drill, ream, or cut out the bone to be removed to either side, as well as any remaining interposed discal material. In the preferred method, utilizing a specially designed Endmill-Drill, it and 40 the Inner Sleeve are removed as a unit, safely carrying away the bone and disc debris trapped within them from the spinal

Step 7. If required, a thread forming Tap with penetration limiting means to control the depth of insertion, is then 45 inserted through the Outer Sleeve.

Step 8. The prepared implant is then inserted utilizing the specialized Driver unit. It should be noted that the implant may be coated with, made of, and/or loaded with substances consistent with bony fusion. However, in the preferred 50 embodiment, the implant is treated with bone promoting and inducing substances, but is loaded with materials suitable for participating in a fusion.

While substances both natural and artificial are covered by the present invention, the preferred embodiment is in regard 55 to the use of the patient's own bone by the following method. A hollow Trephine is utilized to harvest a core of bone from the posterior superior aspect of the iliac crest adjacent the sacroiliac joint. This core of bone is at its outside diameter, slightly smaller than the inside diameter of the spinal implant 60 to be loaded, but longer than the spinal implant. Utilizing an instrument designed for that purpose, the core of bone is then injected from within the Trephine into the central cavity of the implant causing a superabundance of the bone material within the implant such that the bone material tends to press out through the openings communicating with the outside surface of the implant.

14

Step 9. Using the Driver Extractor instrument, the prepared implant is threaded into the prepared interspace. The instrumentation is removed from that side of the spine and attention is then redirected to the first side of the spine. A small retractor is utilized to move the dural sac and traversing nerve root medially and to protect them and allowing the direct visualization of the retained Short Distractor unit. Without removing the Short Distractor, it is reassembled to its shaft portion, essentially reconstituting itself into a Long Distractor. With the inserted implant now acting as the distractor on the opposite side, the Long Distractor is utilized to guide the Outer Sleeve down where it is impacted as described in Step 5.

Steps 6 & 7 are then repeated, completing the procedure at that level. The wound is then irrigated and closed in the 15 routine manner.

#### Representative Example of the Preferred Method

Through preoperative templating of the patient's anterior accordance with the remaining steps as described in the pre- 20 posterior, lateral, and axially imaged MRI scan in conjunction with translucent overlays of the various sized implants, the correct implant diameter and length are accurately assessed, as well as the correct amount of distraction needed to restore the interspace to its premorbid height. The patient is then properly positioned and a bilateral partial discectomy performed via paired semihemilaminotomies.

> For the purpose of this example, it will be assumed that by preoperative assessment it was determined that the correct implant would have an external diameter of 18 mm and be 26 mm long. Further, the distraction necessary to restore the height of the interspace would be approximately 10 mm. The dural sac and traversing nerve root would then be retracted medially and protected, while a Long Distractor having an outside diameter to the barrel portion corresponding to the implant to be inserted, that is 18 mm, and having a diameter at the working end of perhaps 8 mm, would be inserted. This then being found to be slightly smaller than optimal by direct observation, a Convertible Distractor having in its barrel portion an 18 mm outside diameter, but having in its working portion a 10 mm diameter would then be inserted. Direct observation and/or x-ray then confirming the ideal distraction, the Convertible Distractor would then be disassembled, the barrel and head portion removed, and the Short Distractor portion left deeply embedded and with its flanged head flat against the canal floor and deep to the neural structures. It would then be safe to allow the dural sac and nerve root to return is to their normal positions, which would be superficial to the flanged portion of the Short Distractor.

> Attention would then be directed to the contralateral side. The dural sac and nerve root would then be retracted medially on this second side, and a Long Distractor with an 18 mm diameter barrel portion and a 10 mm working portion would then be inserted into the interspace and driven flush to the bone if necessary, such impaction imploding any osteophytes not already removed, and assuring that the shoulder portion of the barrel comes to lie flat against the posterior aspects of the adjacent bodies. With the dural sac and nerve root still safely retracted, the Outer Sleeve would then be placed over the Long Distractor and utilizing the Driver Cap and a mallet, seated to the optimal depth.

> In the preferred embodiment, the Long Distractor is then removed and the Inner Sleeve is inserted into the Outer Sleeve. Since the purpose of the Inner Sleeve is to support the drill and allow for the increased size of the implant over the size of the drill, thus making it possible for the insertion of the implant to occur through the Outer Sleeve, the Inner Sleeve therefore measures 18 mm in its outside diameter, and 16.6

mm in its inside diameter. This allows it to fit within the Outer Sleeve, the diameter of which is 18.1 mm and to admit the drill bit which is 16.5 mm in diameter.

Following the drilling procedure, the Drill and Inner Sleeve are removed as a single unit with the trapped interposed cartilaginous and bony debris. The depth of drill penetration is preset and limited by the fixed rigid column of the Outer Sleeve. In this example, the space will be prepared to a depth of 28 mm in anticipation of countersinking a 26 mm long implant at least 2 mm. If a Tap were to be utilized, it would be inserted at this time and be appropriate to the minor and major diameters of the implant to be inserted and as with the Drill, controlled for its depth of penetration. The spinal implant would then be prepared for implantation by utilizing a Trephine to harvest a core of posterior iliac bone greater than 30 mm long and approximately 14.5 mm in diameter.

Using the Bone Loading Device, this core of bone would be forcefully injected into the internal chamber of the spinal implant which would then be capped. Cap end forward, the 20 fully loaded implant would then be attached to the Insertion Driver, down the Outer Sleeve and screwed into place with the depth of penetration limited by the Insertion instrument. The Insertion Driver is then unscrewed from the implant and removed from the Outer Sleeve. With the dural sac and nerve 25 root retracted and protected, the Outer Sleeve would then be removed. This would complete the fusion procedure on that side, and then as described, the procedure would be repeated on the other (first) side of the same interspace.

#### Alternative Methods

An alternative and extremely useful method is the "Trephine method". Its advantages include that it may be used in conjunction with the preferred embodiment substituting the 35 use of a hollow, tubular cutter, called a Trephine for the use of the Drill in Step 5 of the preferred embodiment. Additionally, it may be utilized so as to obviate the need for the placement of the Short Distractor and to allow the procedure to be effectively performed from start to finish on one side prior to 40 initiating the procedure on the opposite side, and while nevertheless maintaining distraction at the site of the bone removal.

The following is a description of the "Trephine Method". Having completed the exposure of the interspace on at least 45 one side, the dural sac and nerve root are retracted. A Long Distractor differing from the Long Solid Bodied Distractor of the preferred embodiment only in that the barrel portion is of a precisely lesser diameter than the spinal implant. As in the preferred embodiment, the Outer Sleeve has an inner diam- 50 eter only slightly greater than the implant to be inserted. Therefore, at this time, a first Inner Sleeve is inserted into the Outer Sleeve to make up the difference between the outside diameter of the Long Distractor and the inside diameter of the Outer Sleeve. With the Outer Sleeve and first Inner Sleeve 55 thus assembled, they are placed over the Long Distractor and the Outer Sleeve is optimally seated using the Impaction Cap. The Cap and first Inner Sleeve are removed, but the Long Distractor and Outer Sleeve are left in place.

With the Long Distractor maintaining optimal distraction 60 and with the Outer Sleeve locking the vertebrae together so as to resist any movement of the vertebrae, a hollow, tubular cutter known as a Trephine is then inserted over the Long Distractor and its barrel portion and within the Outer Sleeve. The Trephine, which is stopped out to the appropriate depth, 65 can then be utilized to cut equal arcs of bone from the opposed vertebral endplates.

16

Alternatively, a second Inner Sleeve may be placed within the Outer Sleeve prior to placing the Trephine over the Long Distractor and within that second sleeve. This second Inner Sleeve would be just greater in its internal diameter than the Long Distractor and just smaller in its outside diameter than the inner diameter of the Outer Sleeve. While it would provide enhanced stability to the Trephine, provision would then need to be made in the way of large flutes passing longitudinally or obliquely along the outer surface of the Distractor to its barrel portion to accommodate the bony and cartilaginous debris generated during the cutting procedure.

Following the use of the Trephine to the appropriate depth by either of these methods, the Trephine, the Long Distractor, and the second Inner Sleeve, if utilized, are all removed. Since the Trephine cuts two arcs of bone but does not ream them out, a shafted instrument with a perpendicular cutting portion at its working end is then inserted parallel to the disc space and then rotated through an arc of motion cutting the bases of the two longitudinally cut arcs, thus freeing them for removal through the Outer Sleeve. The space may then be tapped if required, and the implant is inserted as per the preferred method. As already mentioned, the "Trephine Method" can be used with or without the use of the Short Distractor on the contralateral side.

#### Applications of Method in Other Areas of the Spine

The following method is the preferred embodiment for performing anterior interbody fusion in the thoracic and lum30 bar spines. It is also appropriate in the cervical spine when the width of the spine anteriorly is sufficient so that it is possible to place two implants side by side and such that each intrudes at least several millimeters into the substance of the opposed vertebrae and for the length of the implants.

The interspace to be fused is adequately exposed and the soft tissues and vital structures retracted and protected to either side. Visualization of the broad width of the interspace anteriorly is made possible by the absence of the neurological structures in relation to this aspect of the spine. The center line of the anterior aspect of the interspace is noted and marked. The disc is removed using first a knife and then curettes and rongeurs as needed. Alternatively, the disc may be left intact to be removed during the drilling stage of the procedure. However, as per the preferred embodiment of the procedure, having removed the great mass of the nucleus and the greater portion of the annulus anteriorly, Long Distractors with progressively increasing diameters to their working ends are inserted into the interspace at a point midway between the central marking line and the lateral extent of the anterior aspect of the spine as visualized.

The Dual Outer Sleeve with its common Foot Plate and Retention Prongs is then inserted over either a singly placed Long Distractor and then the second Distractor placed, or is placed over both Distractors if already placed. The Dual Outer Sleeve is then seated firmly against the anterior aspect of the spine. Any spurs which would interfere with the flush seating of the Foot Plate to the anterior aspect of the spine should be removed prior to inserting the Long Distractors. Once the Outer Sleeve has been optimally seated, one of the Long Distractors is removed and in its place is inserted an Inner Sleeve and drill bit. The drill bit has as its outside diameter the minor diameter of the implant to be inserted. The Inner Sleeve is essentially equal in thickness to the difference between the minor and major diameters of the threaded implant.

A Stopped Drill is then utilized to prepare the opposed vertebral surfaces and to remove any remaining disc material

interposed. If required, a Stopped Tap may be inserted through the Outer Sleeve and into the interspace to create a thread form. The properly prepared implant is then affixed to the Insertion Driver and passed through the Outer Sleeve down into the interspace and inserted until its depth of penetration is limited by the stop on the Insertion Driver. With the implant itself now in a position to act as a distractor, the Long Distractor is then removed from the contralateral side and the procedure repeated. When both implants are firmly in place, the outer sleeve may then be removed. The amount of countersinking of the implants may then be adjusted under direct vision.

Detailed Description of the Preferred Embodiment

#### Method and Instrumentation

In the preferred embodiment, the disc (D) between adjacent vertebrae (V) is approached via bilateral paired semi-hemilaminotomies of the adjacent vertebrae. In the preferred embodiment the supraspinous ligament, the interspinous ligament, the spinous process, portions of the lamina, and most of the facet joints are preserved. However, while less desirable, these structures may be removed.

In the preferred method, a bilateral partial nuclear discectomy is then performed through bilateral openings created through the posterior aspect of the annulus fibrosus. While considered less desirable, disc excision can be delayed and performed simultaneously with the vertebral bone resection 30 during the drilling procedure. Starting on the first side a dural nerve root retractor is placed such that the dural sac and lower nerve root are retracted medially allowing exposure to one side of a portion of two adjacent vertebral bodies and the interposed disc posteriorly.

Referring now to FIG. 1, preferably after removing some portion of nuclear disc material, a Long Distractor 100 is inserted under direct vision into the intervertebral space. The disc penetrating portion 102 is essentially cylindrical with a bullet-shaped front end 103 and a shoulder portion 104 where 40 the penetrating portion 102 extends from barrel 106. The penetrating portion 102 urges the vertebral bodies apart, facilitating the introduction of the instruments. Long Distractors with sequentially increasing diameter penetrating portions 102 are then introduced. As the optimal diameter of 45 penetrating portion 102 is achieved, the vertebral bodies to either side are forced into full congruence and thus become parallel, not only to the penetrating portion 102, but to each other. At this time, any remaining excrescences of bone of the posterior vertebral bodies adjacent the posterior disc which 50 have not already been removed are flattened flush to the vertebral body by the forced impaction, such as by hitting with a hammer flat surface 109 of crown 110, driving the shoulder 104 against the lipped portions of vertebrae V. Because of the forced opposition of the vertebral endplates to 55 portion 102 with optimal distraction, unit 100 will then come to lie absolutely perpendicular to the plane of the posterior bodies and absolutely parallel to the vertebral endplates, allowing optimal alignment for the procedure to be performed.

Penetrating portion 102 is available in various diameters, but all are of a constant length, which is less than the known depth of the interspace. This combined with the circumferential shoulder 104, which is too large to fit within the interspace, protects against the danger of overpenetration. Barrel 106 is of the same diameter as the external diameter of the device to be implanted. A recessed portion 108 below the

18

crown 110 allows for the Long Distractor 100 to be engaged by an extractor unit shown in FIG. 9.

In the preferred embodiment, a Convertible Long Distractor 113 is used on the first side of the spine. As shown in FIG. 2, the Convertible Long Distractor 113 has a barrel portion 152 separable from the Short Distractor portion 120. While the initial distraction may be performed with a solid Long Distractor, as the optimal distraction is approached the appropriate Convertible Long Distractor is utilized. The Convertible Long Distractor 113 consists of a Short Distractor portion 120 and a barrel 152 having a rectangular projection 134 at one end. The Short Distractor 120 has an increased diameter head 128, a rectangular slot 118 and an internal threaded opening 114. The barrel 152 is hollow and has an internal 15 shaft 111 terminating in a large diameter hexagonal crown 115 at one end and a reduced diameter portion 112. The crown has a detent portion 117 in its flat surface. The other end of the shaft 111 has a threaded small member 116 that corresponds to threaded opening 114. The shaft 111 is prevented from removal from the barrel 152 by set pin 119, passing through the wall of barrel 152 in a convenient manner. The Short Distractor portion 120 is removably attached to the barrel portion 152 via the mating of female rectangular slot 118 and the male mating member 134. The mating held together by utilizing knob 136 to drive the crown 110 connected to interior shaft 111 having a threaded working end screw 116 that threads into the female aperture 118 of the Short Distractor portion 120.

Cap 136 has an open socket 138 for fitting around Crown 115 and engages the reduced diameter hexagonal portion 112 so as to permit the rotation of shaft 111 and threaded male member 116. A detent ball 150 in the inside of the socket 138 engages detent 117 in the crown 115, holding them together.

The Short Distractor portion 120 of FIGS. 2, 3, and 3A-3F are designed to provide for high stability when temporarily situated so as to resist inadvertent migration while the surgeon is working on the second side. To that end, the embodiment of the Short Distractor 120 shown in FIGS. 3 and 3A has a pair of sharp pegs 126, to embed into the opposing vertebral bodies and forward facing ratchetings 124, that further resist backward movement. FIGS. 3B and 3C, which show the preferred embodiment, are side and top views of an alternative embodiment of the distractor portion such that the distractor portion to be interposed between the vertebrae is essentially cylindrical, but with circumferential forward facing ratchetings 124.

A further alternative embodiment is shown in FIGS. 3D and 3E. This is a more rectangularized design, with forward facing ratchetings, without the sharp prongs 126 of FIG. 3. FIG. 3F is a side view of a further embodiment of the Short Distractor 120 shown with knurling, to increase the interference with the bone surface so as to add stability to the unit and to resist dislodgment. To this end, it is apparent that the working ends of both the Long and Short Distractors can have a variety of configurations consistent with their purpose, and that surface irregularities as well as the shape of the ends themselves, with or without prongs 126, may be utilized to make the Short Distractor 120 more resistant to migration.

Once the ideal distraction has been achieved on the first side of the spine, the Convertible Distractor is dissociated, leaving Short Distractor 120 in place with its rounded external end 128, safely on the canal floor and deep to the dural sac and nerve root.

As shown in FIG. 4, the surgeon then moves to the other side of the spine at the same disc (D) level, and retracts the dural sac and nerve root medially, exposing the disc on that side. Long Distractors 100 are then sequentially inserted into

the disc space until the diameter of the distractor on the second side is at least as big, as that on the first side. If because of some asymmetry of the interspace a larger diameter distractor is required on the second side to achieve the ideal distraction as compared to the first side, then the second side 5 is fitted with a Short Distractor of the larger diameter, and the surgeon would then return back to the first side. In that event, the first side Short Distractor would then be removed and the Long Distractor 100 corresponding to the increased diameter of the already placed Short Distractor 120 would then be inserted. In either event, the operation is continued by working on the one side where the Long Distractor is in place. In this regard, it should be noted, that by the use of such a device as the Michelson Spinal Surgery Frame, it may be possible to obtain adequate distraction preoperatively such that the sur- 15 geon is either disinclined to use a distractor, or to simply place the correct Long Distractor on the first side and then proceed with the surgical procedure on that side before moving to the opposite side. These variations are within the scope of the present invention.

19

The Long Distractor now serves as both a centering post and an alignment rod for the hollow Outer Sleeve 140 shown in FIG. 5 which is fitted over the Long Distractor 100, shown by phantom lines 101 in FIG. 5. The Outer Sleeve 140 is metal and has a sharp toothed front end 142 that is capable of 25 penetrating into and holding fast the two adjacent vertebrae (V). Interrupting the circumferential sharp teeth of 142 are flat, planar areas 152 which serve to resist the further insertion of the sharp teeth into the vertebral bodies. The toothed front end 142 of the Outer Sleeve 140 is a continuation of the 30 tubular shaft 144, which in turn is connected to circumferentially enlarged tubular back end 146 having a knurled outer surface 148 for easier manipulation. An alternative embodiment of an Outer Sleeve incorporates an expansile key hole and slot configuration 154 to either side of shaft 144 along the 35 mid-plane of the interspace and parallel to it such that the end 142 resists the collapse of the vertebrae (V) to either side of the disc (D), but may nevertheless allow for their further distraction, in the event the only diameter or the root diameter of the implant is larger than the hole drilled.

A Driver Cap 160 in the form of an impaction cap has at is its far end a flat, closed-back surface 162 and at its other end a broad, circular opening. The Driver Cap 160 fits over both the Outer Sleeve 140 and the Long Distractor 100. As the Driver Cap 160 is seated, interior surface 170 circumferen- 45 tially engages portion 146 of the Outer Sleeve until the back end 172 engages the internal shoulder 164. As mallet blows are applied to surface 162, that force is transmitted via the internal shoulder 164 to the Outer Sleeve 140 via its far end 172, seating teeth 142 into the vertebral bodies adjacent the 50 disc space D and to the depth of the teeth 142 to the flat portions 152. As the Outer Sleeve 140 is advanced forward, crown portion 110 of the Long Distractor is allowed to protrude within the Driver Cap 160 unobstructed until it contacts the interior flat surface 168. Once crown 110 comes into 55 contact with the flat interior surface 168, then further taps of the mallet will not advance the Outer Sleeve, any further motion being resisted by the flat shoulder portion 104 of the Long Distractor abutting the hard surfaces of the posterior vertebral bodies. In this way, the Outer Sleeve 140 is safely 60 and assuredly inserted to its optimal depth and rigidly securing the two opposed vertebrae as shown in FIG. 6.

The Cap 160 is then removed and the Distractor Puller 200 of FIG. 9 utilized to remove the Long Distractor 100 from the spine leaving the Outer Sleeve 140 in place. The Distractor 65 Puller 200 has front portion 202, a mid portion 204, and a back handle portion 206. At the front portion 202 of the

20

Distractor Puller 200, a socket 208 is connected to one end of shaft 210 which at its far end is connected to back handle portion 206. The socket 208 has defined within it a cavity 212 that is open at its front end and funnelized on the interior aspect of its sides. The cavity 212 is constructed so that the head of the Distractor Puller 200 and the partially circumferential flange 218 engages the circumferential recess 108 of the Distractor 100. The entrance to cavity 212 is slightly funnelized, and the leading edges of flange 218 slightly rounded to facilitate the engagement of recess 108 and head 110 of Distractor 100, which is further facilitated in that the Driver Cap 160 leaves portion 108 of Distractor 100 precisely flush with the back surface 172 of the Outer Sleeve 140. This provides a large, flat surface 172 to precisely guide surface 230 of socket 208, and open portion 212 around head 110 while flange 218 engages recess 108. The springloaded detent ball 228 engages hemispherical depression 112 in the crown 110, shown in FIG. 2. This springloaded detent 228 in engagement with complimentary indent 218 protects against 20 the inadvertent dissociation of the Long Distractor from the Puller 200 after the Distractor has been removed from within the Outer Sleeve 140 and prior to its removal from the wound. Once out of the body, the two instruments are easily disassociated by freeing the crown portion 110 from cavity 212 by a manual force applied perpendicular to their relative long axes at this location.

A cylindrical and free removable weight 216 is fitted around shaft 210 between the front portion 202 and the rear handle portion 206. Gently, but repeatedly sliding the weight 216 along shaft 210 and driven rearwardly against flat surface 228, transmits a rearward vector to proximal end 202 and thereby to the Long Distractor 100 to which it is engaged.

Paired extended handle 224 and 226, allow the surgeon to resist any excessive rearward motion as the instrument is used to liberate the Long Distractor 100. Paired handles 224 and 226 are also useful in that they allow a rotational directing of portion 208, via the shaft 210. This allows the surgeon to control and manipulate rotationally the orientation of the opening of cavity 212 to facilitate its application, to the head 110 of the distractor 100.

The Distractor Puller 200 is a significant improvement over the alternatives of striking a remover instrument with an independent hammer over the exposed surgical wound, or manually extracting the distractor by forcefully pulling. The use of a free hammer over the open wound is dangerous because the neural structures can be impacted on the back swing which is made even more likely by the effects of gravity on the mallet head. Manual extraction by pulling is dangerous because of the significant interference fit of portion 102 within the spine such that significant force would be required to remove the Distractor 100, and if force were not coaxial then the Outer Sleeve might be dislodged or misaligned. Further, once the flat portion 102 became free of the interspace, all resistance to withdrawal would be lost and in the face of the considerable force necessary to free it, the Distractor 100 might easily become projectile imparting injury to the patient and/or the surgeon.

Once the Long Distractor 100 has been fully removed from the Outer Sleeve 140, the toothed end 142 of the Outer Sleeve 140, working in conjunction with the Short Distractor 120 on the contralateral side rigidly maintains the relative position of the adjacent vertebrae V. Further, since the remainder of the procedure on that side of the spine occurs entirely through the protective Outer Sleeve 140, and as the nerves and dural sac are external to that Outer Sleeve and superficial to the toothed end 142 of the Outer Sleeve 140, which is firmly embedded into the adjacent vertebrae V, the Outer Sleeve 140 serves to

insure the safety of these delicate neural structures. Further, since the Outer Sleeve 140 is of a fixed length and rigid, its flat rearward surface 172 may be used as a stop to the advancement of all instruments placed through the Outer Sleeve 140, thus protecting against accidental overpenetration. Further, 5 the Outer Sleeve 140 assures that the further procedure to be performed will occur coaxial to the disc space D and further, be symmetrical in regard to each of the opposed vertebral surfaces.

FIG. 10B is a posterior view of the spine at this stage of the 10 procedure, showing a Short Distractor 120 in place on one side of the spine and the bottom portion of Outer Sleeve 140 in place on the opposite side of the spine.

Referring to FIG. 11A, an Inner Sleeve 242 is inserted from the rear within the Outer Sleeve 140. This Inner Sleeve has a 15 collar portion 244 of a known thickness which seats against the top edge surface 172 of Outer Sleeve 140. The cylindrical barrel portion of Inner Sleeve 242 comes to approximate the posterior aspect of the vertebral bodies interior the Outer Sleeve when fully seated. A Drill **240**, having a known 20 selected length is then introduced through the rearward aperture of the Inner Sleeve 242 and utilized to ream out the arcs of bone which it engages from the opposed vertebral endplates as well as any discal material within its path down to its predetermined and limited depth. The Drill 240, has a narrow 25 engagement portion 246, which allows it to be affixed to a drill mechanism which may be either a manual or a power unit. A circumferential collar 248 of an increased diameter serves to limit the depth of penetration of the drill 240 and may be fixed, or lockably adjustable.

Not shown here, but well known to those skilled in the art, are various mechanisms to lockably adjust such instruments as drills. Such mechanisms include, but are not limited to, the use of collets, threaded shafts with lock nuts, and flanges engaging grooves forced therein by either a cap pulled over 35 the flanges or screwed down upon them.

In the preferred embodiment, the forward cutting edge 252 of Drill 240 is a modification of a large fluted drill design such to that the end resembles an end cutting mill which may contain any workable number of cutting surfaces, but preferably four or more, and such cutting surfaces being relatively shallow such that the advancement of the instrument occurs more slowly. The outside diameter of the Drill 240 corresponds to the minor diameter of the threaded spinal implant. The Inner Sleeve 242 has an inner diameter slightly greater 45 than that dimension and its outer diameter is slightly smaller than the inside diameter of the Outer Sleeve 140 which has the same outer diameter as the major diameter of the threaded implant.

The drill shaft of drill 240 comprises an upper portion 243, 50 a central recessed portion 256 of a smaller diameter and a lower cutting drill portion 250. The upper portion 243 and lower portion 256 of the drill 240 have the same outside diameter.

The Inner Sleeve **242** serves many functions. First, it provides a more intimate drill guide for drill **240** in the event a smaller diameter hole is to be drilled than that of the inside diameter of the Outer Sleeve **140**. Second, since it now guides the Drill, it allows for the Outer Sleeve **140** to have an internal diameter large enough to admit the threaded spinal implant, 60 which is indeed considerably larger in diameter than the Drill **240** itself.

If a larger Outer Sleeve **140** were utilized absent the Inner Sleeve **242**, then the Drill **240** would be free to wander within the confines of that greater space and would not reliably make 65 parallel cuts removing equal portions of bone from the adjacent vertebrae V. Further, the bone removal not only needs to

22

be equal, but must be correctly oriented in three dimensions. That is, the path of the Drill **240** must be equally centered within the disc space, parallel the endplates, and parallel to the sagittal axis dissecting the interspace.

A further purpose of the Inner Sleeve 242 is that it may be removed simultaneously with the Drill 240, thereby trapping the debris, both cartilaginous and bony generated during the drilling procedure, which are guided rearward by the large flutes 251 of Drill portion 250, where they are collected around recessed portion 256 between the recessed portion 256 and the inner wall of the Inner Sleeve 242 are there contained therein. Thus, by removing the Drill 240 in conjunction with the Inner Sleeve 242, all of the debris generated by the reaming procedure is safely removed from the spinal canal and wound area.

Further, if the disc tissue in the area to be reamed has been removed previously, as per the preferred method, then the patient's own bone of good quality and useful within the operation will then be contained between the Inner Sleeve 242 and the shaft portion 256. Once away from the surgical wound, this material may be used to load the spinal implant or placed deep within the interspace to participate in the fusion.

The method of actually producing the surgical hole within the spine is variable. As shown in FIG. 11C, in an alternative embodiment Drill end 250 has a forward projecting nipple 260, which itself is bullet-shaped in its leading aspect so as to ease its entrance into the disc space and to urge the vertebrae apart. Nipple 260 is distracting, stabilizing as it resists any tendency of the vertebrae to move together, is self-centering to the Drill portion 250 when working in conjunction with Sleeves 140 and 242, and virtually assures the symmetrical resection of bone from the opposed vertebral surfaces.

The alternative "Trephine Method" referred to earlier in this application, is shown in FIG. 11B. In this alternative, a Long Distractor 100 is left in place after the Outer Sleeve 140 is seated. The Long Distractor 100 in this case differs from the Long Distractor of the preferred embodiment in that its outside diameter Of the barrel 106 is of a smaller diameter than in the prior version. This is made necessary because regardless of the method, the hole to be formed corresponds to the minor diameter of the spinal implant. Trephine 270, a hollow, tubular member with sharp cutting teeth 251 at its proximal end, has a wall thickness and since the outside diameter of that trephine 270 must correspond to the root diameter of the implant, then the wall thickness of the trephine 270 must be allowed for by a corresponding reduction in the diameter of the Long Distractor 100.

A further modification of the Long Distractor 100 to the "Trephine Method" would use longitudinal grooves (not shown) along the barrel surface 106 for the purpose of transmitting any debris generated during the cutting procedure, rearward. Since the cutting element is both centered and aligned by the Long Distractor, the use of the Inner Sleeve 242 is not mandatory, but may once again be useful in controlling the path of the debris. To that end, little debris is generated in the "Trephine Method" as the bony arcs are not so much being reamed out and removed as they are simply being cut into the bone where these arcs of bone are left connected at their far ends. Thus, when the Trephining Method has been completed and the Trephine 270 and Inner Sleeve 242 removed, unlike in the preferred embodiment where the hole is drilled out, it remains necessary to remove both the two arcs of bone, and any interposed material. Nevertheless, this is very easily performed by various means, one of which is depicted in FIG.

Instrument 272 consisting of a shaft 276 attached off center to the lower surface 273 handle 274. The shaft 274 terminates

in a cutting arm 278. The instrument 272 is inserted through Outer Sleeve 140 where the lower surface 273 of handle 274 abuts the top 172 of the Outer Sleeve 140, both stopping downward motion of instrument 272 and precisely placing the perpendicularly cutting arm 278 of instrument 272 so that as handle portion 274 is rotated, the cutting arm 278 is also rotated, cutting the arcs of bone and liberating them from their last attachments. These portions of bone are then removed utilizing this instrument or a long forceps, and then placed within the implants or otherwise used to participate in the 10 fusion.

While in the preferred embodiment of the present invention the spinal implant I, is essentially self-tapping, if the bone is unusually hard it may be desirable to form the thread pattern within the interspace prior to the insertion of the implant I. To that end, as shown in FIG. 12, Tap 280 has a threadcutting portion 282 connected by a shaft 286 to a handle portion 292, which has been designed to give mechanical advantage to the rotation of the instrument for the purpose of cutting threads. The lower portion of handle 290 has a forward facing flat 20 surface 288 too large to fit through the opening of Outer Sleeve 140 which thus safely limits the depth of penetration of the cutting element 282. This tap 280 is further made safe by blunt end 294 which will engage the uncut portions of the vertebral bone just prior to the engagement of shoulder 288 25 against surface 172. This feature allows the surgeon to appreciate a less harsh resistance as the blunt nose 294 encounters the remaining unresected bone for the drill hole and prior to the sudden increase in resistance caused by the seating of shoulder 288 against top edge 172, which first resistance 30 serves as a warning to the surgeon to discontinue the tapping procedure. Thus, the surgeon has both visual (as shoulder 288 approaches top edge 172) and tactile warnings to avoid stripping the thread form. Tap end 282 is highly specialized for its specific purpose. Rearward to the specialized blunt tip 294 is 35 a truncated bullet-shaped area 298 which ramps up to the constant diameter intermediate the cutting ridges 296. Ramp portion 298 urges the opposed vertebral bodies apart, which motion is resisted by Outer Sleeve 140, thus progressively driving, the sharp leading edges of thread forms 296 into the 40 vertebral bodies. The periodic longitudinal grooves 284 interrupting the thread forms, which may number 1 to 8, but preferably 4, function to accumulate the bony material which is removed during the thread cutting process. In that regard, in the ideal embodiment, the thread cutting form is designed to 45 compress the bone to be formed rather than to trough through it. Further, while both the major and minor diameters of the Tap 280 may be varied, in the preferred embodiment, the minor diameter corresponds to the minor diameter of the implant I, but the major diameter is slightly less than the 50 major diameter of the implant.

With Tap 280 now removed, and Sleeve 140 still in place, the surgical site is now fully prepared to receive the spinal implant I. In the preferred embodiment of the spinal implant, the implant has been enhanced by the use of, application to, 55 and filling with fusion promoting, enhancing, and participating substances and factors. Thus, the implant may be fully prepared for insertion as provided to the operating surgeon. However, at the present time, human bone is most commonly used as the graft material of choice, with the patient's own 60 bone being considered the best source.

FIG. 14a shows a trephine 300 with an exceedingly sharp front cutting edge 302 for quickly and cleanly coring into the patient's posterior iliac crest, or any other bony tissue, and for the purpose of producing a core of bone then contained within 65 the hollow 304 of the trephine 300. Trephine 300 has a rear portion 306 with a pair of diametrically opposed slots 310,

24

and disposed clockwise from their longitudinally oriented rearward facing openings so as to engage diametrically and opposing members 312 of Drive unit 308, by which trephine 300 may be attached to either a hand or power drill. It can be appreciated that engagement mechanism 312 is stable during the clockwise cutting procedure, and yet allows for the rapid disconnection of the two components once the cutting is completed.

Because of the high interference between the graft and the inner wall of hollow portion 304, and the relative weakness of the cancellous bone being harvested, it is possible to remove the Trephine 300 while still drilling, and to have it extract the core of bone with it. However, in the highly unlikely event that the core of bone would remain fixed at its base, then with the drive mechanism 308 removed, a corkscrew 408 shown in FIG. 14C is introduced though the central opening of rear portion 306 and threaded down and through the core of bone within 304 and to the depth of teeth 302. The tip 318 of the corkscrew 408, which extends substantially on line with the outer envelope of the corkscrew, then cuts radially through the base of the bone core. As the handle portion 314 of the corkscrew 408 abuts the flat, rearward surface of portion 306 and it can no longer advance. As corkscrew 408 is continued to be turned further, it will cause the core of bone to be pulled rearward, as in removing a cork from a wine bottle. Trephine 300 has a barrel portion 304 continuous with sharp toothed portion 302 having an inner diameter just less than the inner diameter of the spinal implant I to be loaded.

The Trephine 300 with its core of harvested bone is then placed as shown in FIG. 14B, through opening 340 of Implant Bone Loading device 320, where the barrel portion 304 then passes through and is stopped by circular flange 344. The plunger shaft 326 of instrument 320 is then prepared for attachment by rotating knob 332 counterclockwise such that the plunger 372 is pulled via the long threaded shaft portion 328 back to the base of collar 330 at its proximal end. In this position, knob 332 is considerably extended rearward from collar 330. With plunger shaft 326 in this position, the plunger head 372 is inserted into the central hollow of portion 306 of Trephine 300 as the proximal cylindrical portion of collar 330 then follows it, such that the plunger 372 then occupies the rearward portion of barrel 304 and the proximal cylindrical portion of collar 330 occupies the central hollow of portion 306. A pair of diametrically opposed radially projecting arms 346 on collar 330 are then advanced longitudinally into diametrically opposed paired L slots 340 and then rotated clockwise to complete this assembly.

At the other end of instrument 320, a spinal implant I is engaged through its female rectangular slot 364 by a rectangular protruding bar extending from rearward facing surface of end plug 324, (not shown) and secured there by knob 334 which extends as a rod through a central aperture within end plug 324 to extend at the far end as a small bolt which threads to a female aperture centered within the female slot 364 of the spinal implant. With the spinal implant I secured to end plug 324 and the opposite end of the implant I presenting as a hollow, tubular opening, end plug 324 is advanced into device 320 where it is secured by rotationally engaging diametrically opposed L-shaped slots 321. With device 320 fully assembled, end 302 of trephine 300 lies coaxial and opposed to the open end of implant I.

As shown in FIG. 15, as knob 332 is then rotated clockwise, the plunger 372 proximal the threaded shaft 328 is then forcibly, but controllably driven forward down the barrel 304 ejecting the bone graft directly into the spinal implant I. As the bone graft is greater in length than the interior of the spinal implant, with further compression the bone is forced into the

radially disposed apertures through the wall of the device communicating from the central cavity to the exterior.

End plug **324** is then removed from apparatus **320**. Using end plug **324** as a handle, end cap **374** shown in FIG. **16** is secured to the open end of the spinal implant I. The implant is then disassociated from end plug **324** by rotating knob **334** counterclockwise.

FIG. 16 shows an Implant Driver instrument which may be used to either insert or to remove said implant I. Driver 350 has at its far end 362, a rectangular protrusion 398, which protrusion intimately engages the complimentary rectangular slot 364 of implant I. Protruding from slot 393 of end 362 is threaded portion 353, which extends as a rod through hollow shaft 358 and hollow hand barrel 360 to knob 354 where it can be rotationally controlled. Threaded portion 353 screws into a female aperture central slot 364, urging 353 into 364, and binding them together such that instrument 350 can be rotated via paired and diametrically opposed extending arms 366 and in either direction while maintaining contact with the implant.

Affixed to the Driver **350**, the implant is then introduced through the Outer Sleeve **140** and screwed into the interspace opposed between the two prepared vertebrae V until such time as the leading edge of the Implant Cap **374** reaches the depth of the prepared hole at which time its forward motion is 25 impeded by the bone lying before it which had not been drilled out. This allows for a progressive feel to the surgeon as the implant is screwed home.

As described previously, with the use of the Tap 280, this terminal resistance to further seating provides significant tactile feedback to the surgeon. Again, as with the Tap 280, visual monitoring of the depth of insertion of the implant is provided to the surgeon by observing the progressive approximation of the forward surface 370, of barrel portion 360, as it approaches the rearward facing surface 172 of Outer Sleeve 140. Nevertheless, a final safety mechanism, when the full depth of insertion has been achieved, surface 370 of instrument 350 will abut surface 172 of the Outer Sleeve 140, prohibiting any further installation of the spinal implant.

Once the implant has been fully installed, the Driver **350** is dissociated from the implant by turning knob **354** in a counterclockwise direction. The Driver **350** is then withdrawn from the outer sheath, then the Outer Sleeve **140** is removed. This leaves the implant fully installed and inset to the determined depth as shown in FIG. **18**.

Attention is then redirected to the other, or first, side of the spine. A dural nerve root retractor is used to retract the neural structures medially, bringing into full view the head 128 of the Short Distractor 120, lying flush on the canal floor. Utilizing apparatus 152, extended screw portion 116 is inserted into the female threaded portion 114 of the Short Distractor 120 as the extended rectangular portion 134 of apparatus 152 is engaged to the female rectangular portion 118 of the Short Distractor 120. Then turning rearward facing portions 108 and 110, utilizing the knob 136 of FIG. 2, the Long Distractor configuration is restored.

With the dural sac and nerve roots still retracted and protected, the Outer Sleeve **140** is slipped over the reconstituted 60 Long Distractor and seated using the Driver Cap **162**. The entire sequence of events as described for the implantation of the spinal implant I as already placed, is then repeated such that both spinal implants come to lie side by side within the interspace. Though not necessary, circlage or other internal 65 fixation of the levels to be fused may additionally be performed, and then the wound is closed in the routine manner.

26

Brief Discussion with Reference to the Drawings of the Preferred Method and Instrumentation for Anterior Interbody Fusion Incorporating Intercorporeal Predistraction and Utilizing a Guarded Sleeve System is Disclosed

Because of the absence of the spinal cord and nerve roots, it is generally possible to visualize in one instance the entire width of the disc space from side to side throughout the cervical, thoracic, or lumbar spine. In the preferred embodiment of the anterior interbody fusion, implants are placed side by side from anterior to posterior parallel to the interspace and extending through into the adjacent vertebral bodies. Where the transverse width of the disc space is insufficient to allow for the use of two implants, each of which would be large enough to protrude to the required depth into the adjacent vertebrae, then a singular and significantly larger implant may be placed centrally. With this in mind, and in light of the very detailed description of the technique and instrumentation already provided in regard to the method of posterior lumbar interbody fusion, a brief discussion of anterior spinal interbody fusion with dual implant installation will suffice, and the method for installation of a large, singular midline graft will become obvious.

The interspace to be fused is exposed anteriorly. The soft tissues are withdrawn and protected to either side, and if necessary, above and below as well. It is then possible to visualize the entire width of the vertebrae anteriorly adjacent that interspace. As discussed above, the surgeon has already templated the appropriate patient radiographs to determine the requisite distraction and optimal implant size. In the preferred method, the surgeon then broadly excises the great bulk of the nuclear disc portion. (Alternatively, the disc can be left to be removed via the drill later.) The surgeon then notes and marks a point midway from side to side anteriorly. He then inserts Long Distractor 100 centering it on a point midway between the point just noted and the lateral extent of the intervertebral space visualized anteriorly. The outer barrel portion 106 of the Distractor 100 utilized, will correspond to 40 the outside diameter of the implants to be installed. The Distractor tips 102 inserted are sequentially larger in diameter until the optimal distraction is achieved. This optimal distraction, although suggested by the initial templating, may be visually and tactilely confirmed as performed. When the optimal distraction is achieved, the vertebral endplates will come into full congruence and parallel to the forward shaft portion 102 of the Distractor 100, causing an alteration in the alignment of the vertebrae and a significant increase in the interference fit and pressurization at the tip, such that the instrument becomes exceedingly stable.

There is a sensation imparted to the surgeon of the tissues having moved through their elastic range to the point where the two adjacent, vertebrae V begin to feel and move as if a single solid. These changes are easily appreciated visually as the vertebrae realign to become congruent to tip 102, and can also easily be appreciated via lateral Roentgenography. However, should the surgeon fail to appreciate that optimal distraction has been achieved and attempt to further distract the interspace, he would find that extremely difficult to do because of the increased resistance as the tissues are moved beyond their range of elastic deformation. Further, there would be no elasticity left to allow the vertebrae to move further apart and the sensation to the surgeon should he attempt to gently tap the oversized Distractor forward with a mallet, would be one of great brittleness.

Returning now to the procedure, when the correct intercorporeal Distractor 100 producing the ideal interspace distrac-

tion having its barrel portion 106 corresponding to the implant to be installed has been inserted, then its exact duplicate is inserted anteriorly equidistant to the other side of the spine. As the barrel portion 106 of Long Distractor 100 is exactly of the same major diameter as the spinal implant I looking coaxially on end, the surgeon can then assess the anticipated side by side relationship of the dual implants when implanted.

As shown in FIGS. 7C and 7D, a Dual Outer Sleeve 340 consisting of a pair of hollow tubes is then introduced over the side by side Long Distractors protruding anteriorly from the spine. The Dual Outer Sleeve 340 is comprised of two hollow tubular members identical in size displaced from each other ideally the sum of the difference between the minor and major 15 diameters of both implants combined, but not less than that difference for one implant, as it is possible to have the threads of one implant nest interposed to the threads of the other, such that they both occupy a common area between them. However, while the preferred embodiment is slightly greater than 20 two times the difference between the major and minor diameters of the implant (the sum of both) the distance may be considerably greater. Whereas in the preferred embodiment extending tubular portions 348 of instrument 340 are parallel, when the area between them 350, is sufficiently great, these 25 elements may be inclined or declined relative to each other such that they either converge or diverge at their proximal ends. Paired tubular structures 348, may be bridged in part or wholly throughout their length, but are rigidly fixed by Foot Plate **344**. In its preferred embodiment, a top view shows the Foot Plate to be essentially rectangular, but without sharp corners.

Other shapes can be utilized. In side view 7D it can be appreciated, that Foot Plate 344 is contoured so as to approximate the shape of the vertebrae anteriorly. Extending forward from Foot Plate 344 are multiple sharp prongs 342 sufficiently long to affix them to the vertebrae. The prongs 342 are limited in length so as to not penetrate too far posteriorly and number from 2 to 10, but preferably 6. As the Dual Outer Sleeve 340 is driven forward utilizing Dual Driver Cap 420, of FIG. 7E, engaging the rearward end 352, the prongs 342 extending from Foot Plate 344 are embedded into the opposed vertebral bodies until their forward motion is inhibited by the curved Foot Plate 344 becoming congruent to and being 45 stopped by, the anterior aspect of the vertebral bodies.

As already taught in FIG. 5, the Dual Driver Cap 420 is of the same design as Single Driver Cap 160, in that there is a recess 354 as per 168, allowing the Outer Sleeve to be fully seated without impeding the rearward projection of the Long Distractor unit. However, unlike in Cap 160, area 354 is more relieved as it is unnecessary for the Dual Cap 420 to contact the Long Distractor through portion 110 to inhibit its forward motion, as the Foot Plate 344 functions to that effect. Further, the Dual Cap 420 for the Dual Outer Sleeve 340 is correspondingly dual itself and engages the rearward facing dual tubular portion 352. Once the Dual Outer Sleeve has been fully seated, the vertebrae adjacent the interspace to be fused are rigidly held via Foot Plate 344 and the prongs 342. Thus, it is possible to remove either one, or if desired, both of the Long Distractor rods utilizing Long Distractor puller 200, as per the method already described. It is then the surgeon's choice to work on one or both sides of the spine. As per previous discussion, the surgeon may drill the interspace utilizing the Inner Sleeve 242 or leave the Long Distractors in place as per the "Trephine Method".

28

Tapping, if necessary, and the insertion of the implants then occurs through the protective Outer Sleeve **340**. Once the implants have been fully inserted, the Outer Sleeve is removed.

Having utilized the Drill method, or "Trephine Method", with or without an Inner Sleeve to prepare the fusion site, it is the preferred embodiment to leave the Outer Sleeve 340 in place as it provides for the ideal placement and alignment of the Tap 280 and implant I.

It is anticipated that the surgeon wishing to work deep within the interspace, or preferring the ability to directly visualize the tap being used, or the implant being inserted, may choose to remove the Outer Sleeve after the insertion of the first prosthesis to maintain stability, or prior to that, which while not the preferred embodiments, are nevertheless within the scope of the present invention.

Alternative Methods to the Preferred Embodiment for Method of Anterior Interbody Fusion

As previously described for the posterior lumbar spine, alternatively, one can employ the "Trephine Method" as has been described in detail.

As a further alternative, it should be noted that the key element in the anterior method is the use of the predistraction principle, where such distraction is maintained by the Outer Sleeve with or without the Long Distractor. Therefore, once the preparation of the interspace has been completed, while not the preferred embodiment, it is nevertheless within the scope of this invention that one could remove the Outer Sleeve as there are no neural structures requiring protection, and insert the implants directly rather than through the Outer Sleeve.

As yet a further alternative of this method, where the height of the distracted interspace is such that the diameter of the implant required to span that height and to embed with sufficient depth into the opposed vertebral bodies is such that it is not possible to place two such implants side by side, then only a single implant which may be of significantly increased diameter, is used and placed centrally within the interspace rather than to either side. The placement of a singular central graft via the present invention method and instrumentation is in keeping with the methods already described and can be performed using either a drill or the "Trephine Method".

Referring to FIGS. **16-18**, a cylindrical embodiment of the spinal implant I of the present invention is shown. In FIG. **16** the implant I is shown attached to the insertion device **350**. In FIGS. **17** and **18** the implant I is shown installed in the disc space D, between the adjacent vertebrae.

The cylindrical implant I comprises a hollow tubular member which in the preferred embodiment is made of an ASTM surgically implantable material, preferably Titanium. The cylindrical implant I is closed at one end and open at the other end covered by a cap 394. The cylindrical implant I has a series of macro-sized openings 390 through the side walls of the cylindrical implant I. A series of external threads 392 are formed on the circumference of the cylindrical implant I. Any variety of threads may be used on the implant. The cap 374 has a hexagonal opening 394 for tightening the cap 374.

While the present invention has been described in association with the implant of a threaded spinal implant, it is recognized that other forms of implants may be used with the present method. For example, dowels, made from bone or artificial materials, knurled or irregularly shaped cylinders or spheres, or any other shaped implants that can be introduced through the outer sleeve may be used. Being able to perform

the procedure through the outer sleeve permits the procedure to be performed safely and quickly, and more accurately.

What is claimed is:

- 1. An apparatus for use in performing human interbody spinal surgery comprising:
  - a guard having a mid-longitudinal axis, and a first passage and a second passage for providing guided access to a disc space and a vertebral body adjacent the disc space, each of said first and second passages having a proximal end, an opposite distal end, a length from said proximal 10 end to said distal end, and being adapted to receive at least a portion of a spinal implant therethrough, said first and second passages being spaced apart from one another along a majority of their respective lengths, said first and second passages being located on first and sec- 15 ond sides, respectively, of a plane bisecting said guard along the mid-longitudinal axis thereof, said first and second passages having first and second cross sections, respectively, in a plane perpendicular to the mid-longitudinal axis, said first and second cross sections of said 20 first and second passages having identical dimensions;
  - a distracting element configured to distract two adjacent vertebral bodies, said distracting element having a blunt leading end, an opposite trailing end, upper and lower load-bearing surfaces between said leading and trailing ends, and a length from said leading end to said trailing end, said distracting element being fixed to a shaft having a portion configured to contact the posterior surfaces of the adjacent vertebral bodies to limit the depth of insertion of said distracting element into the disc space; 30
  - at least one orientation pin for holding said guard to the spine, said orientation pin having a distal end portion configured to extend parallel to and along a majority of the length of said distracting element when said distracting element and said orientation pin are inserted into the spine; and
  - a bone removal device for forming through each of said passages in turn an implantation space.
- 2. The apparatus of claim 1, wherein each of said passages forms a hollow tubular sleeve.
- 3. The apparatus of claim 1, wherein each of said passages has a wall, said walls of said passages being spaced apart the majority of the length of the respective passages.
- **4**. The apparatus of claim **1**, wherein each of said passages has a circular cross section.
- 5. The apparatus of claim 1, wherein said guard has an increased outer dimension portion at its proximal end.
- 6. The apparatus of claim 1, wherein said guard has a footplate.
- 7. The apparatus of claim 1, further comprising means for 50 cooperatively engaging to the proximal end of said guard an impaction end member for receiving an impaction force for driving said guard into the spine.
- **8**. The apparatus of claim **1**, further comprising a cap adapted to engage said proximal end of said guard.
- 9. The apparatus of claim 1, wherein, said passages are bridged in part proximate one of said ends of said passages.

- 10. The apparatus of claim 1, wherein at least one of said orientation pins has a tapered leading end to facilitate placement of at least one of said spine engaging portions into the spine.
- 11. The apparatus of claim 1, wherein at least one of said orientation pins has upper and lower surfaces that are parallel to each other.
- 12. The apparatus of claim 1, further comprising a removable inner guard.
- 13. The apparatus of claim 12, wherein said removable inner guard is a hollow tubular sleeve.
- **14**. The apparatus of claim **12**, wherein said removable inner guard is adapted to be inserted into said guard.
- 15. The apparatus of claim 14, wherein said inner guard has limiting means for limiting the travel of said inner guard within said guard when said inner guard is inserted into said guard.
- 16. The apparatus of claim 12, wherein said inner guard has a collar at one end larger than said passage of said guard.
- 17. The apparatus of claim 1, wherein said bone removal device is a milling instrument.
- **18**. The apparatus of claim **1**, wherein each of said passages has a mid-longitudinal axis, the mid-longitudinal axes of said passages being in a convergent relationship to each other.
- 19. The apparatus of claim 1, wherein each of said passages has a mid-longitudinal axis, the mid-longitudinal axes of said passages being in a divergent relationship to each other.
- 20. The apparatus of claim 1, further comprising an implant driver sized in part for passage through one of said passages of said guard for passing an implant through said guard and into an implantation space.
- 21. The apparatus of claim 20, wherein said implant driver comprises an elongated shaft having means for engaging an implant at one end and means for manipulating said implant driver at the other end of said elongated shaft.
- 22. The apparatus of claim 21, further comprising means for limiting the depth of insertion of said implant driver into said guard
- 23. The apparatus of claim 22, wherein said limiting means includes a portion of said implant driver cooperating with said guard to limit the depth of insertion of said implant driver into said guard.
  - **24**. The apparatus of claim **1**, wherein said upper and lower load bearing surfaces of said distracting element are in a parallel relationship to each other.
  - 25. The apparatus of claim 1, further comprising a tap for insertion through said guard for tapping the two adjacent vertebral bodies.
  - **26**. The apparatus of claim **1**, wherein said apparatus is in combination with a spinal implant.
  - 27. The apparatus of claim 26, wherein said implant is in combination with a fusion promoting material.
  - **28**. The apparatus of claim **27**, wherein said fusion promoting material is harvested bone.

\* \* \* \* \*

#### UNITED STATES PATENT AND TRADEMARK OFFICE

## **CERTIFICATE OF CORRECTION**

PATENT NO. : 7,993,347 B1 Page 1 of 1

APPLICATION NO. : 09/626636 DATED : August 9, 2011

INVENTOR(S) : Gary Karlin Michelson

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

#### Title Page 2, Item (56), U.S. Patents:

Column 1, line 30: change "Entherly" to -- Eatherly --; Column 2, line 58: change "Vichl" to -- Vich --; and Column 2, line 66: change "Babgy" to -- Bagby --.

#### Title Page 3, Item (56), Foreign Patent Documents, Column 2:

Line 6: change "10/1991" to -- 4/1991 --; and

Line 11: delete in its entirety.

#### Title Page 4, Item (56), Other Publications:

Column 1, line 6: change "pencott" to -- pincott --;

Column 1, line 9: change "Advance" to -- Advances --; and

Column 2, after line 50: insert the following missing reference listings:

-- Adams, et al.; Outline of Orthopaedics, Eleventh Edition; Trunk and Spine, p, 194.

Herkowitz, et al.; Principles of Bone Fusion; The Spine, Third Edition; Chapter 44, p. 1739.

Muschler, et al.; The Biology of Spinal Fusion; Spinal Fusion Science and Technique, Cotler and Cotler, pp. 9-13.

Zindrick, et al.; Lumbar Spine Fusion: Different Types and Indications; The Lumbar Spine, Vol. 1, Second Edition, pp. 588-593 (1996).

Gillingham, F.J., et al.; Automatic Patient Monitoring in the Ward; Brit. J. Surg., Vol. 53, No. 10, pp. 864-866 (October 1966).

Maloney, A.F.J., et al.; Clinical and Pathological Observations in Fatal Head injuries, Brit. J. Surg., Vol. 56, No. 1, pp. 23-31 (January 1969). --.

Signed and Sealed this Eighteenth Day of September, 2012

David J. Kappos

Director of the United States Patent and Trademark Office