

US007166129B2

(12) United States Patent

Michelson

(10) Patent No.: US 7,166,129 B2

(45) Date of Patent: *Ja

*Jan. 23, 2007

(54) METHOD FOR FORMING A SPINAL IMPLANT SURFACE CONFIGURATION

(75) Inventor: Gary K. Michelson, Venice, CA (US)

(73) Assignee: Warsaw Orthopedic, Inc., Warsaw, IN

(US)

(*) Notice: Subject to any disclaimer, the term of this

patent is extended or adjusted under 35

U.S.C. 154(b) by 0 days.

This patent is subject to a terminal dis-

claimer.

(21) Appl. No.: 09/921,851

(22) Filed: Aug. 3, 2001

(65) Prior Publication Data

US 2002/0013624 A1 Jan. 31, 2002

Related U.S. Application Data

- (62) Division of application No. 09/457,228, filed on Dec. 8, 1999, now Pat. No. 6,827,740.
- (51) **Int. Cl.** *A61F 2/02* (2006.01)
- (52) U.S. Cl. 623/17.11

See application file for complete search history.

(56) References Cited

U.S. PATENT DOCUMENTS

| 4,349,921 | A | 9/1982 | Kuntz |
|-----------|---|--------|---------------------|
| 4,531,244 | A | 7/1985 | Hamas |
| 4,673,409 | A | 6/1987 | Van Kampen |
| 4,795,472 | A | 1/1989 | Crowinshield et al. |
| 4,834,757 | A | 5/1989 | Brantigan |
| 4,865,603 | A | 9/1989 | Noiles |
| 4,944,763 | A | 7/1990 | Willert et al. |
| | | | |

 4,955,907 A
 9/1990 Ledergerber

 5,019,107 A
 5/1991 Schelhas

 5,306,308 A
 4/1994 Gross et al.

(Continued)

FOREIGN PATENT DOCUMENTS

EP 0 637 440 2/1995


(Continued)


Primary Examiner—Bruce Snow (74) Attorney, Agent, or Firm—Martin & Ferraro, LLP

(57) ABSTRACT

The present invention is a specialized implant having opposed surfaces for engaging each of the vertebral bodies adjacent a disc space into which the implant is implanted. The surface comprises arrayed projections having at least one forward facing facet directed at least in part toward the leading end of the implant and at least one rearward portion directed at least in part toward the opposite trailing end of the implant. Each of the forward facet and rearward portion has a length and a slope. The length of the forward facet is longer than the length of the rearward facet. The slope of the rearward facet is steeper than the slope of the forward facet. The surface projections also have opposed side facets directed generally toward the sides of the implant. The side facets are located between the forward facet and rearward facet and converge toward each other in a direction away from the base of the surface projections. The surface may also include projections having left and right forward side facets and a rearward facet. The surface further may include projections having left and right rearward side facets and a forward facet.


57 Claims, 7 Drawing Sheets

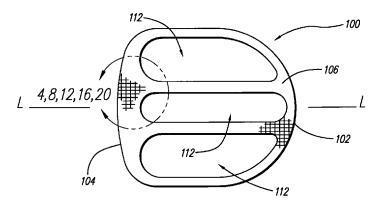


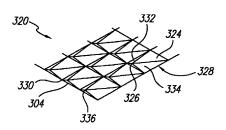



US 7,166,129 B2Page 2


| U.S. PA | ATENT | DOCUMENTS | | ,143,033 | | 11/2000 1/2001 | Paul et al. Suddaby | | |
|---------------------------------------|--------|---------------------|--------------------------|----------|-------|-------------------|------------------------|-----------|--|
| 5,425,772 A | 6/1995 | Brantigan | | .190.414 | | | Young et al. | | |
| 5,522,899 A | 6/1996 | Michelson | | ,258,125 | | | Paul et al | 623/17.11 | |
| 5,553,476 A | 9/1996 | Oehy et al. | | ,325,827 | | | Lin | 020/1//11 | |
| 5,593,409 A | 1/1997 | Michelson | | .425.920 | | | Hamada | | |
| 5,609,635 A | 3/1997 | Michelson | 6 | ,432,106 | В1 | 8/2002 | Fraser | | |
| | | Wagner et al. | 6 | ,511,509 | В1 | 1/2003 | Ford et al. | | |
| 5,713,899 A | 2/1998 | Marnay et al. | 6 | ,592,624 | В1 | 7/2003 | Fraser et al. | | |
| 5,716,412 A | 2/1998 | DeCarlo, Jr. et al. | 6 | ,660,038 | В1 | 12/2003 | Boyer, II et al. | | |
| | | Oehy et al. | 2002/ | 0068978 | A1 | 6/2002 | Camino et al. | | |
| , , | | Michelson | | | | | | | |
| , , , , , , , , , , , , , , , , , , , | | Thalgott | FOREIGN PATENT DOCUMENTS | | | | | | |
| , , | 3/1999 | | | | | | | | |
| 5,899,941 A | 5/1999 | Nishijima et al. | SU | | 1 107 | | 8/1984 | | |
| 5,989,289 A 1 | 1/1999 | Coates et al. | WO | WC | 98/58 | 3604 ' | * 12/1998 | | |
| 6,019,793 A | 2/2000 | Perren et al. | | | | | | | |
| 6,080,158 A | 6/2000 | Lin | * cited | l by exa | miner | | | | |

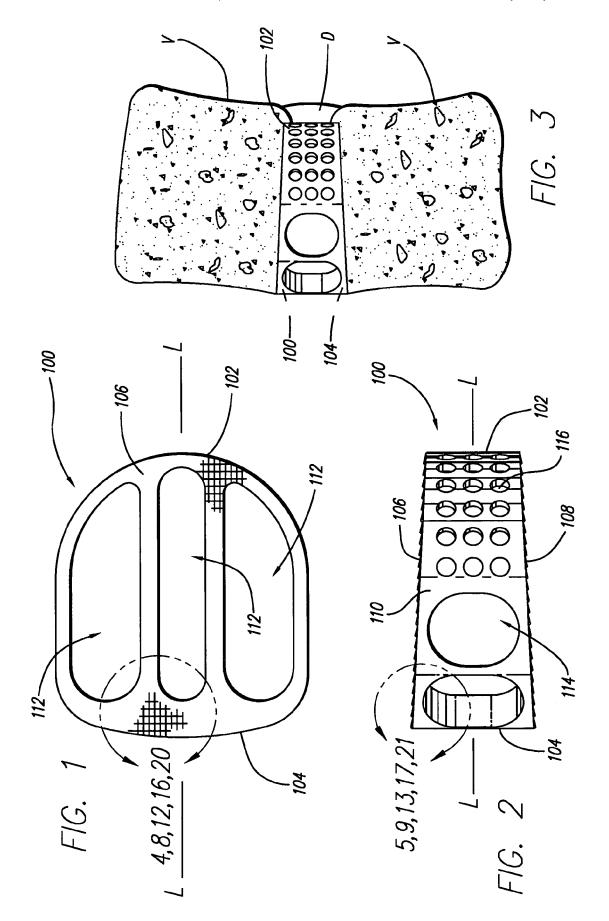


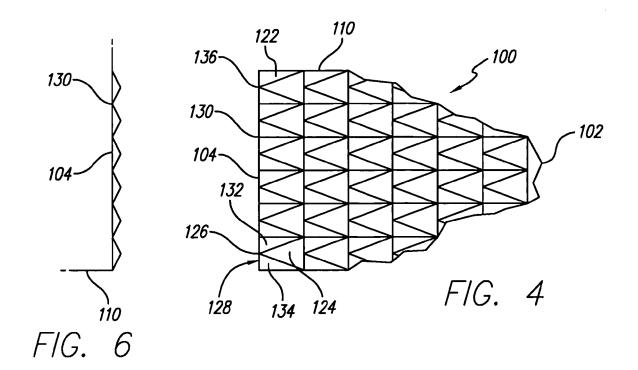


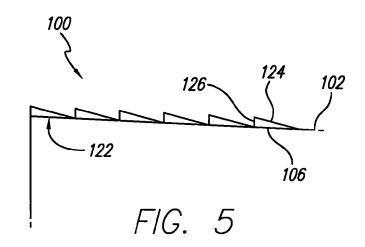

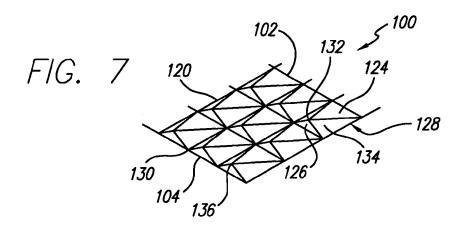

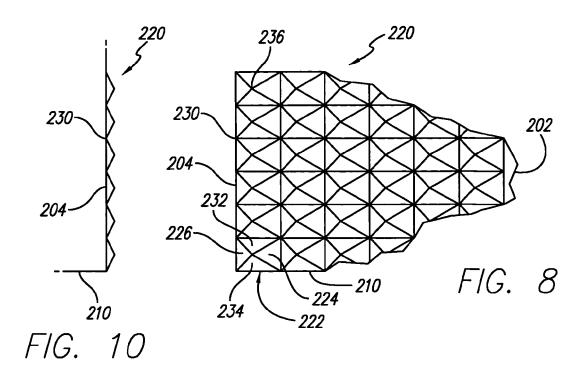

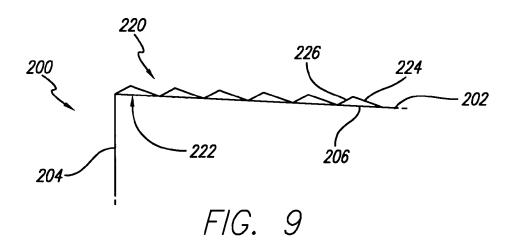

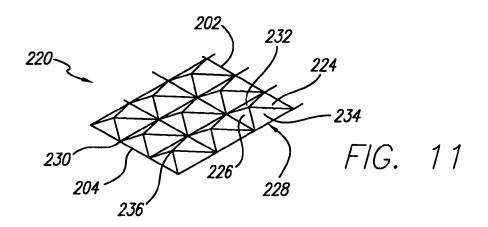





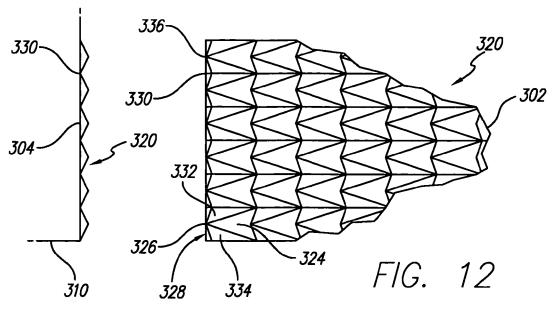



FIG. 14

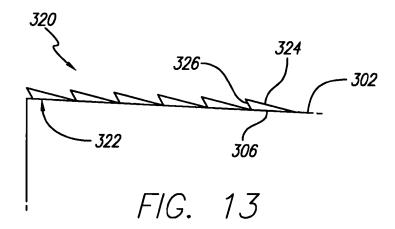


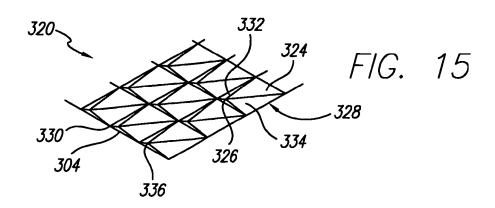

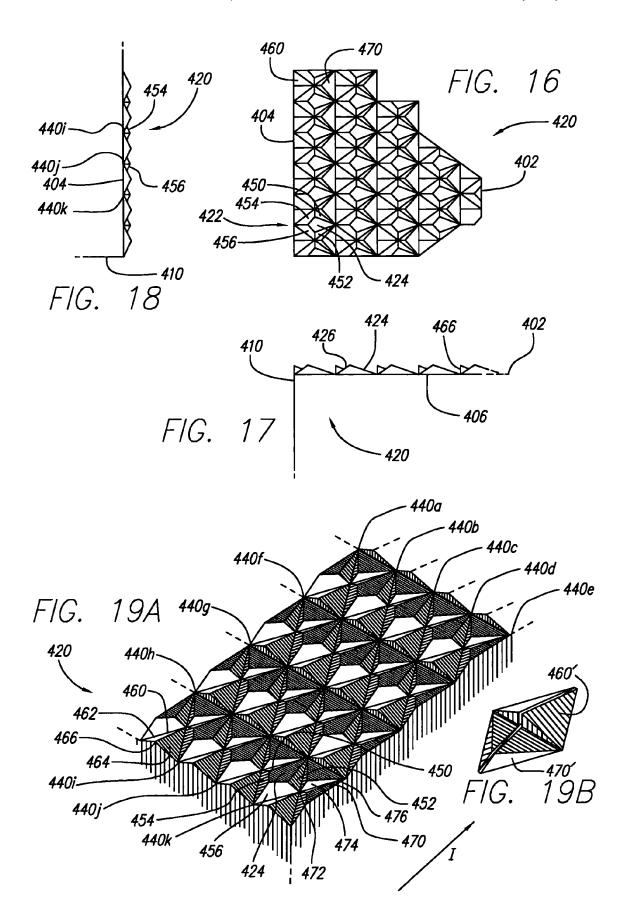

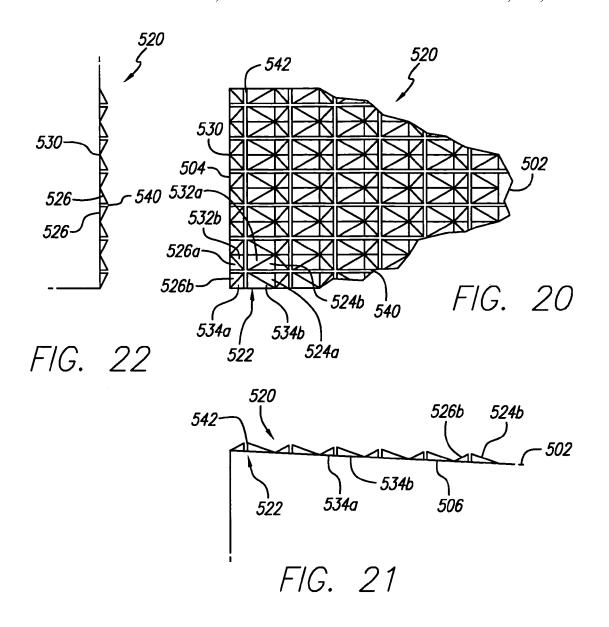

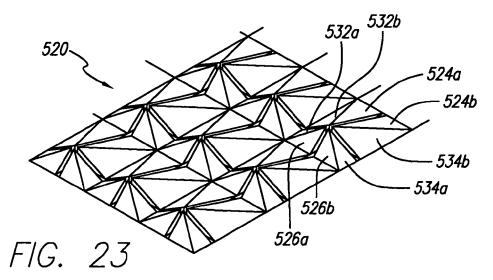

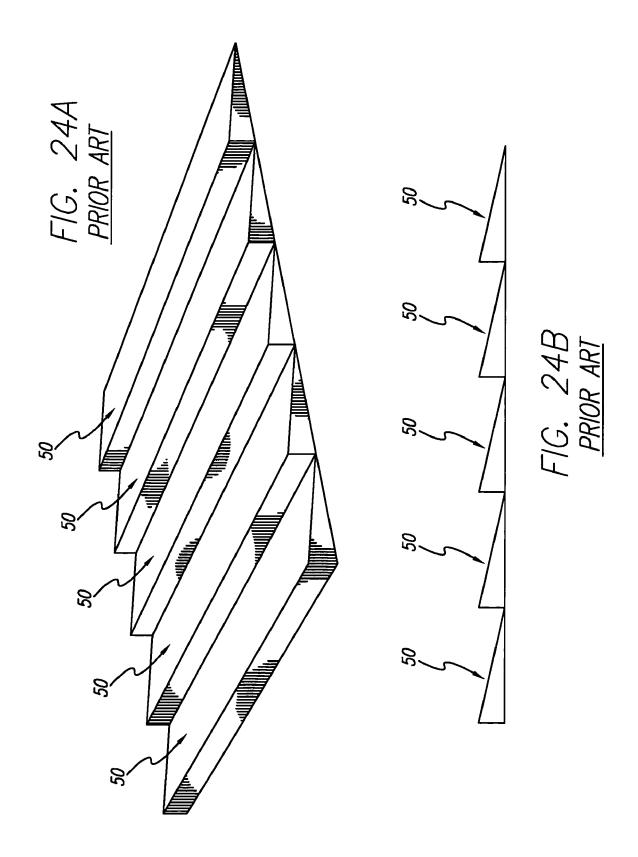

METHOD FOR FORMING A SPINAL IMPLANT SURFACE CONFIGURATION

BACKGROUND OF THE INVENTION

This application is a divisional of application 09/457,228, filed Dec. 8, 1999, now U.S. Pat. No. 6,827,740.

DESCRIPTION OF THE RELATED ART

The present invention relates to interbody spinal implants for placement between adjacent vertebral bodies of a human spine, and more specifically to a specialized surface for such interbody implants, for engaging the adjacent vertebral bodies. Vital to the functioning of all interbody spinal 15 implants is their ability to remain properly located within the spine after installation. In U.S. Pat. Nos. 5,593,409 and 5,609,635 Michelson described the use of surface roughenings, such as knurling or ratcheting on the opposed upper and lower vertebral body engaging surfaces of interbody 20 spinal fusion implants. Knurling has been particularly beneficial for increasing the grip of the implant surface to the adjacent vertebral bodies in a rather uniform manner without a directional bias. Spinal implants have a propensity to move in a particular direction, which is opposite to their path of 25 insertion, because this is the path of least resistance. Such propensity to move is further increased when the opposed upper and lower vertebral body engaging surfaces are in angular relationship to each other, such that they are spaced further apart at the implant's trailing end than at the 30 implant's leading end. In such circumstances where it is desirable then to gain stability in resistance to a particular direction of movement of the interbody spinal implant, the use of a plurality of forward facing ratchetings on the implant's vertebral body engaging surfaces has been pref- 35 erable to the previously described knurling for that purpose.

The term "ratcheting" as used herein is defined as a plurality of angular teeth or ridges or protrusions projecting from the surface of an implant to resist motion of the implant at least in one direction. The phrase "forward facing ratch- 40 etings" as used herein is defined as a ratcheting having at least one forward facing facet having a length greater than a rearward facing facet and an angle from the implant surface from which the forward facing facet arises that is less steep than the angle of the rearward facet. On an implant 45 surface, forward facing ratchetings facilitate the insertion of the implant in one direction and after insertion, resisting movement of the implant in a direction opposite to the direction of insertion. An example of forward facing ratchetings of the prior art is shown in partial fragmentary view 50 in FIGS. 24A and 24B, generally referred to by the reference numeral 50.

Knurled surfaces of the related art provide some stability in all directions, but lack the ability to resist a particular direction of motion preferentially. The above-described 55 ratcheted surface best resists motion in a particular direction. There exists a need for an improved interbody spinal implant surface configuration, wherein the opposed upper and lower vertebral body engaging surfaces of the implant are configured to be resistant to implant movement in all directions, 60 and preferentially or in particularly in one direction.

SUMMARY OF THE INVENTION

The present invention relates to interbody spinal implants 65 having a specialized surface configuration on the opposed exterior surfaces adapted for engaging the vertebral bodies

2

adjacent a disc space into which the interbody implant is to be implanted. Such an implant surface configuration has utility with a wide variety of shapes of interbody spinal implants where enhanced skeletal fixation is desired. Such an implant surface configuration can provide for enhanced stability, increased surface area, and a surface for the delivery of fusion promoting substances other than bone. In a preferred embodiment, the implant surface can provide for resisting motion in all directions, and particularly in at least one direction, such as counter to the direction of insertion of the implant.

While various embodiments of the present invention are presented by way of example only and not limitation, common to each of them is that the surface configuration incorporates a plurality of spatially integrated surface projections having at least one forward facing facet directed at least in part toward the leading end of the implant and at least one rearward portion directed at least in part toward the opposite trailing end of the implant. By way of example and not limitation, the rearward portion may be a facet, a line, or an edge of the rearward aspect of the surface projection formed where two facets come together. Each of the forward and rearward facets have a length and a slope. The length of the forward facet is longer than the length of the rearward facet. The slope of the rearward facet is steeper than the slope of the forward facet. In various embodiments, the surface projections also have opposed side facets directed generally toward the sides of the implant. The side facets are located between the forward facet and rearward facet and may converge toward each other in a direction away from the base of the surface projections. The surface comprises multifaceted ratcheted projections that are organized in geometrically disposed fields or arrays which are at a minimum located on at least a portion of the opposed vertebral body engaging surfaces of the implant. From the teachings disclosed herein, it is appreciated that the surface projections can be geometrically arranged in a pattern wherein at least a portion of the projection is aligned along a longitudinal, horizontal, diagonal, or curved line. The upper and lower surfaces of the implant can be at least in part arcuate or planar and can converge along a portion or all of the length of the implant.

In various preferred embodiments of the present invention, the rearward facets of the surface projections can be perpendicular or at angles greater or less than 90 degrees to at least one of the upper or lower surfaces of the implant from which the projections arise. The opposed side facets of the surface projections can have at least a first portion in a plane at an angle to the longitudinal axis of the implant. The opposed side facets can intersect each other, and can converge to form a peak at the top of each of the surface projections. The peaks can be aligned along lines that are perpendicular, parallel, or diagonal to the longitudinal axis of the implant. The surface projections can be cleaved such as by longitudinal and/or horizontal cuts to increase the number of exposed sides of the projections and thus increase the available surface area to contact and engage the bone of the adjacent vertebral bodies and increase the number of recessed areas to contain fusion promoting substances. Alternatively, the peaks of each surface projection can be cleaved, truncated, or flattened at least in part.

The surface projections can include a left forward side facet and a right forward side facet directed toward the leading end and sides, respectively, of the implant. Similarly, the surface projections can include a left rearward side facet and a right rearward side facet directed toward the trailing end and sides, respectively, of the implant. The side facets of

adjacent surface projections can be spaced apart to define a groove therebetween. A plurality of adjacent surface projections can be spaced apart to form a plurality of grooves that can be parallel or at an angle to the longitudinal axis of the implant, wherein the angle can be less than 90 degrees. The 5 grooves can have a horizontal cross section that is a V-shape, U-shape, or a box-like shape, for example.

Sequential projections can be positioned on an implant wherein each surface projection has forward facing facets facing the same direction, such that consecutive projections 10 are oriented forward facing facet to rearward facing facet. The lower most portion of the slope of the forward facing facet of a first surface projection in a sequence can be coincident with the rearward facet of the next surface projection in the sequence. Alternatively, the forward facet 15 of a first surface projection and the rearward facet of the next surface projection in a sequence can be spaced apart and the space can be at least in part flat, curved, or any other surface contour suitable for the intended use. The surface projections can be oriented relative to one another to form fields or 20 arrays that further can be geometrically disposed relative to one another, preferably in a pattern wherein at least a portion of the projection is aligned along a longitudinal, horizontal, diagonal, or curved line.

The surface configuration of the present invention can be 25 formed by casting, machining, or any other techniques known to one of ordinary skill in the art. The present surface configuration may readily be machined by milling the implant surface from side to side, across the upper and lower vertebral body engaging surfaces, to form ratchetings gen- 30 erally disposed perpendicular to the long axis of the implant and generally formed facing to the insertion end of the implant. The ratchetings may be cross machined with an angled cutting face to form grooves passing through the ratchetings. For example, a milling machine having a cutting 35 tool, with a V-shaped profile, can be run through the plurality of ratchetings parallel to the longitudinal axis of the implant to form the above-described surface. In a preferred embodiment, the V-shaped cutting tool of the milling machine has opposed cutting faces with an angle of approximately 90 40 degrees to each other, which faces are each at a 45-degree angle to the plane of the surfaces being machined. Without departing from the scope of the present invention, the angle of the cutting faces can be more or less than 90 degrees, and the angle of the cutting face to the surface to be cut can be 45 more or less than 45 degrees. It is appreciated that rather than the cutting element being run parallel to the longitudinal axis of the implant, the cutting element could be run at some other angle. By way of example only and not limitation, this angle could be at 45 degrees to the longitudinal axis 50 of the implant and to the projections. Each surface projection could then be formed by a cutter crossing in two passes to form two grooves at a 90 degree angle to each other.

The surface of the present invention for engaging each of the adjacent vertebral bodies may be incorporated into 55 various types of spinal implants. Such spinal implants may be for the purpose of achieving interbody spinal fusion, or for stabilizing a device to space apart and allow motion between the adjacent vertebral bodies. Such spinal implants appropriate for the intended purpose. Such materials would include, but are not limited to, implant quality metals, including, but not limited to, titanium and its alloys, surgical grade plastics and plastic composites which may or may not be bioresorbable, ceramics, and cortical bone. Some 65 examples of interbody spinal implants that may benefit from the present teaching, include but are not limited to the

following patents and applications by Michelson which are incorporated by reference herein: U.S. Pat. Nos. 5,015,247; 5,522,899; 5,593,409; 5,609,635; 5,860,973; and application Ser. No. 08/480,904.

BRIEF DESCRIPTION OF THE DRAWINGS

FIG. 1 is a top plan view of an interbody spinal implant having a surface configuration in accordance with the present invention.

FIG. 2 is a side elevation view of the spinal implant of FIG. 1.

FIG. 3 is a side elevation view of the interbody spinal implant of FIG. 1 installed in an implantation site formed across the disc space between two adjacent vertebral bodies of the spine shown in partial cross-section.

FIG. 4 is an enlarged fragmentary top plan view of an implant surface of one embodiment of the present invention from a view taken along area 4 of FIG. 1.

FIG. 5 is a fragmentary side elevation view of the implant surface of FIG. 4 from a view taken along area 5 of FIG. 2.

FIG. 6 is a fragmentary end elevation view of FIG. 4.

FIG. 7 is a fragmentary perspective view of the implant surface of FIG. 4.

FIG. 8 is an enlarged fragmentary top plan view of a second embodiment of the implant surface of the present invention from a view taken along area 8 of FIG. 1.

FIG. 9 is a fragmentary side elevation view of the implant surface of FIG. 8 from a view taken along area 9 of FIG. 2.

FIG. 10 is a fragmentary end view of the implant surface of FIG. 8.

FIG. 11 is a fragmentary perspective view of the implant surface of FIG. 8.

FIG. 12 is an enlarged fragmentary top plan view of a third embodiment of the implant surface of the present invention from a view taken along area 12 of FIG. 1.

FIG. 13 is a fragmentary side elevation view of the implant surface of FIG. 12 from a view taken along area 13 of FIG. 2.

FIG. 14 is a fragmentary end view of FIG. 12.

FIG. 15 is a fragmentary perspective view of the implant surface of FIG. 12.

FIG. 16 is an enlarged fragmentary top plan view of a fourth embodiment of the implant surface of the present invention from a view taken along area 16 of FIG. 1.

FIG. 17 is a fragmentary side elevation view of the implant surface of FIG. 16 from a view taken along area 17 of FIG. 2.

FIG. 18 is a fragmentary end view of FIG. 16.

FIG. 19A is an enlarged fragmentary perspective view of the implant surface of FIG. 16.

FIG. 19B is an enlarged fragmentary perspective view of a variation on the second and third surface projections of the fourth embodiment of the implant surface of the present invention with a cleave therethrough.

FIG. 20 is an enlarged fragmentary top plan view of a fifth embodiment of the implant surface of the present invention from a view taken along area 20 of FIG. 1.

FIG. 21 is a fragmentary side elevation view of the may comprise any artificial or naturally occurring material 60 implant surface of FIG. 20 from a view taken along line 21 of FIG. 2.

FIG. 22 is a fragmentary end view of FIG. 20.

FIG. 23 is an enlarged fragmentary perspective view of the implant surface of FIG. 20.

FIGS. 24A and 24B are perspective and side elevation views, respectively, of a prior art implant surface having forward facing ratchetings.

DETAILED DESCRIPTION OF THE DRAWINGS

As shown in FIGS. 1–7, an interbody spinal implant 100 has a leading end 102, a trailing end 104, an upper surface 106, a lower surface 108, and a side wall 110 between upper 5 and lower surfaces 106, 108. Upper and lower surfaces 106, 108 may converge from trailing end 104 to leading end 102 along a longitudinal axis L of implant 100 as shown, or may diverge, be parallel, or any combination thereof. Upper and lower surfaces 106, 108 are configured to be placed against 10 and in contact or engagement with the bone of vertebral bodies V of two adjacent vertebrae adjacent disc D of the human spine. Upper and lower surfaces 106,108 and side wall 110 may include large and/or small openings 112, 114, and 116, respectively, to permit bone growth into and 15 through implant 100. Upper and lower surfaces 106,108 of implant 100 can be generally planar as shown in the figures, or can be opposed arcuate surfaces as shown and described in U.S. Pat. No. 5,593,409, incorporated herein by reference, or any other configuration suitable for the desired use.

As shown in detail in FIGS. 4-7, at least a portion of upper and lower surfaces 106, 108 of implant 100 have a surface configuration generally referred to by the numeral 120. In accordance with a first embodiment of the present invention, surface configuration 120 includes surface pro- 25 jections 122 configured to facilitate insertion of implant 100 into an implantation site while resisting expulsion of implant 100 in a direction opposite to the direction of insertion. Each of surface projections 122 has an angled forward facet 124 directed at least in part toward leading end 102 of implant 30 100 and a rearward facet 126 directed at least in part toward trailing end 104 of implant 100. Forward facet 124 has a length greater than the length of rearward facet 126. Rearward facet 126 has a slope that is steeper than the slope of forward facet 124. In this embodiment, the base of rearward 35 facet 126 forms an angle of approximately 90 degrees with respect to upper and/or lower surfaces 106, 108 of implant 100. It is appreciated that the angle of the base of rearward facet 126 with respect to upper and/or lower surfaces 106, 108 of implant 100 may be perpendicular to, greater than 40 perpendicular to, or less than perpendicular to the base of the surface where the facet arises. Forward facet 124 forms an angle in the range of approximately 10 to 60 degrees, with 25-45 degrees being preferred, with respect to upper and/or lower surfaces 106, 108. Each one of surface projections 122 45 also has a left side facet 132 and a right side facet 134 directed toward the sides of implant 100.

In this embodiment of surface configuration 120, a plurality of surface projections 122 are spaced apart laterally (side to side) by longitudinal grooves 130 formed along the 50 longitudinal axis L of implant 100. In one embodiment, longitudinal grooves 130 have a V-shaped horizontal crosssection. The lower most portions of left and right side facets 132, 134 of consecutive side-by-side projections 122 can be coincident with each other or may be spaced apart, any space 55 therebetween can be at least in part flat, curved, sloped or otherwise configured. Each surface projection 122 has left and right side facets 132, 134 that converge to form a high point or peak 136 at the top of each surface projection 122. Each peak 136 can be aligned along lines that are perpen- 60 dicular, parallel, and/or diagonally oriented to longitudinal axis L of implant 100. The left and right side facets 132,134 resist side-to-side motion of implant 100 after it is inserted into the implantation space. Peaks 136 engage the bone of vertebral bodies V adjacent to implant 100 in the implanta- 65 tion site. It is appreciated that in a variation of the present invention, the peaks may be modified such as to be truncated

6

or cut off to have a broader rather than shaper upper most surface. Moreover, the peaks can be cleaved in one or more directions so as to increase the surface area useful for engaging the bone of the vertebral bodies. The relieved areas of the cleaved projections are useful for containing and carrying fusion promoting substances other than bone such as bone morphogenetic proteins and genetic materials coding for the production of bone, or bone itself which could by way of example be in the form of a paste. It is further appreciated that for all the various embodiments of the surface configuration of the present invention, longitudinal grooves 130 can have horizontal cross-sections in a variety of configurations such as, without limitation, square-shaped or U-shaped configurations.

Sequential projections can be positioned on an implant wherein each surface projection has forward facing facets facing the same direction such that consecutive projections are oriented forward facing facet to rearward facing facet. The lower most portion of the slope of the forward facing 20 facet of a first surface projection in a sequence can be coincident with the rearward facet of the next surface projection in the sequence. Alternatively, the forward facet of a first surface projection and the rearward facet of the next surface projection in a sequence can be spaced apart and the space can be at least in part flat, curved, or any other surface configuration suitable for the intended use. The surface projections can be oriented relative to one another to form an array and are preferably geometrically disposed relative to one another in a pattern wherein at least a portion of the projection is aligned along a longitudinal, horizontal, diagonal, or curved line. Further, it is appreciated that the surface of the present invention can be useful with spinal implants of various configurations, including configurations wherein at least one of leading end, trailing end, and sides of the spinal implant is curved. By way of example and not limitation, the leading end, trailing end, and sides of the spinal implant can form an oval, an oblong, or a circle. As shown in FIGS. 8-11, a second embodiment of the surface configuration of the present invention is generally referred to by the numeral 220. Surface configuration 220 includes surface projections 222 to facilitate insertion of implant 100 into an implantation site while resisting expulsion of implant 100 in a direction opposite to the direction of insertion. Each of surface projections 222 has an angled forward facet 224 directed at least in part toward leading end 202 of implant 100 and a rearward facet 226 directed at least in part toward trailing end 204 of implant 100. Forward facet 224 has a length greater than the length of rearward facet 226. Rearward facet 226 has a slope that is steeper than the slope of forward facet 224. In this embodiment, the base of rearward facet 226 forms an angle of approximately 45 degrees with respect to upper and/or lower surfaces 206, 208 of implant 100. Each one of surface projections 222 has a left side facet 232 and a right side facet 234 directed toward the sides of implant 100, and forward facet 224 and rearward facet 226.

In this embodiment of surface configuration 220, longitudinal grooves 230 have a V-shaped horizontal cross-section. The lower most portions of left and right side facets 232, 234 of consecutive side-by-side projections 222 can be coincident with each other or may be spaced apart, any space therebetween can be at least in part flat, curved, sloped or otherwise configured. Each surface projection has left and right side facets 232, 234 that converge to form a high point or peak 236 at the top of each surface projections 222. Each peak 236 can be aligned along lines that are perpendicular, parallel, and/or diagonally oriented to the longitudinal axis L of implant 100. The left and right side facets 232, 234

resist side-to-side motion of implant 100 after it is inserted into the implantation space. Peaks 236 engage the bone of the vertebral bodies adjacent to implant 100 in the implantation site.

As shown in FIGS. 12-15, a third embodiment of the 5 surface configuration of the present invention is generally referred to by the numeral 320 is shown. Surface configuration 320 includes surface projections 322 to facilitate insertion of implant 100 into an implantation site while resisting expulsion of implant 100 in a direction opposite to 10 the direction of insertion. Each of surface projections 322 has an angled forward facet 324 directed at least in part toward leading end 302 of implant 100 and a rearward facet 326 directed at least in part toward trailing end 304 of implant 100. Forward facet 324 has a length greater than the 15 length of rearward facet 326. Rearward facet 326 has a slope that is steeper than the slope of forward facet 324. In this embodiment, the base of rearward facet 326 is "back cut" to form an angle greater than 90 degrees with respect to upper and/or lower surfaces 306, 308 of implant 100. The con-20 figuration of rearward facet 326 further enhances resistance of motion of the implant in a direction opposite to the direction of insertion. It is appreciated that the angle of the base of rearward facet 326 with respect to upper and/or lower surfaces 306, 308 of implant 100 can be any other 25 angle suitable for the intended purpose of the present invention. Each one of surface projections 322 has a left side facet 332 and a right side facet 334 directed toward the sides of implant 100, and a forward facet 324 and a rearward facet

In this embodiment of surface configuration 320, longitudinal grooves 330 have a V-shaped horizontal cross section. The lower most portions of left and right side facets 332, 334 of consecutive side-by-side projections 322 can be coincident with each other or may be spaced apart, and any 35 space therebetween can be at least in part flat, curved, sloped or otherwise configured. Each surface projection 322 has left and right side facets 332, 334 that converge to form a high point or peak 336 at the top of each surface projection 322. Each peak 336 can be aligned along lines that are perpen- 40 dicular, parallel, and/or diagonally oriented to the longitudinal axis L of implant 100. The left and right side facets 332, 334 resist side-to-side motion of implant 100 after it is inserted into the implantation space. Peaks 336 engage the bone of vertebral bodies V adjacent to implant 100 in the 45 implantation site.

As shown in FIGS. 16–19B, a fourth embodiment of the surface configuration of the present invention is generally referred to by the numeral 420. Surface configuration 420 includes surface projections 422 configured to facilitate 50 insertion of implant 100 in the direction of insertion into an implantation site while resisting expulsion of implant 100 in a direction opposite to the direction of insertion. Each of surface projections 422 has an angled forward facet 424 directed toward leading end 402 of implant 100 and a 55 rearward portion 426 directed toward trailing end 404 of implant 100. Forward facet 424 has a length greater than the length of rearward portion 426. Rearward portion 426 has a slope that is steeper than the slope of forward facet 424. In this embodiment, the base of rearward portion 426 forms an 60 angle of approximately 90 degrees with respect to upper and/or lower surfaces 406, 408 of implant 100. Rearward portion 426 can be a portion of surface projection 422, such as a facet, an edge, or a line for example. Each one of surface projections 422 has a left side forward facet 450, a right side 65 forward facet 452, a left side rearward facet 454, and a right side rearward facet 456 directed toward the front and sides,

8

and directed toward the rear and sides of implant 100, respectively, and forward facet 424 and rearward portion 426

Surface configuration 420 can further include a second plurality of surface projections 460 having at least a left forward side facet 462 and a right forward side facet 464 directed at least in part toward leading end 402 and sides of implant 100, respectively, and at least one rearward facet 466 directed at least in part toward trailing end 400. Left and right forward side facets 462, 464 have at least a first portion in a plane at an angle to the longitudinal axis of implant 100. Second surface projections 460 can be interspersed with surface projections 422.

Surface configuration 420 can further comprise a third plurality of surface projections 470 having at least a left rearward side facet 472 and a right rearward side facet 474 directed at least in part toward trailing end 404 and sides of implant 100, respectively, and at least one forward facet 476 directed at least in part toward leading end 402. Left and right rearward side facets 472, 474 have at least a first portion in a plane at an angle to the longitudinal axis of implant 100. Third surface projections 470 can be interspersed with surface projections 422 and/or second surface projections 460. Surface projections 422 may have a length approximating the combined length of second surface projections 460 and third surface projections 470.

In this embodiment, surface configuration 420 has angled grooves 440a-k that form a plurality of surface projections **422**. In this example, angled grooves **440***a*–*k* are formed at an angle that is approximately 45 degrees to longitudinal axis L of spinal implant 100 and in this example, angled grooves 440a-k are approximately 90 degrees to one another. The angled grooves 440a-k can be formed, if machined, by first passing a cutting element at a 45 degree angle to the longitudinal axis L of implant 100 and then passing the cutting element at a 90 degree angle to the path of the first pass of the cutting element, or otherwise formed by casting, molding, and other methods for forming a surface configuration. It is appreciated that angled grooves 400a-k can be formed at various angles to the longitudinal axis L of implant 100 and to each other. For example, such angles can be less than 180 degrees.

In this embodiment of surface configuration 420, angled grooves 440a–k have a V-shaped horizontal cross-section. Each surface projection 422 has left and right side facets 432 and 434 that are convergent and form a high point or peak 436 at the top of each surface projections 422. Each peak 436 can be aligned along lines that are horizontally, longitudinally, and/or diagonally oriented along implant 100. The left and right side forward and rearward facets 450, 452, 454, 456 function to prevent side-to-side motion of implant 100after it is inserted into the implantation space. Peaks 436 may also function like teeth to engage the bone of vertebral bodies V adjacent to the implant in the implantation site.

FIG. 19B shows a variation of second and third surface projections 460', 470' that can be cleaved in one or more directions to increase the number of exposed sides of each projection and thus increase the surface area of the implant bone engaging surface available to contact the bone of the vertebral bodies. A preferred embodiment of this variation of the second and third surface projections 460', 470' are cleaved by a longitudinal groove.

As shown in FIGS. 20–23, a fifth embodiment of the surface configuration of the present invention is generally referred to by the numeral 520. Surface configuration 520 includes surface projections 522 to facilitate insertion of implant 100 into an implantation site while resisting expul-

sion of implant 100 in a direction opposite to the direction of insertion. Surface projections 522 can be cleaved in one or more directions to increase the number of exposed sides of each projection and thus increase the surface area of the implant bone engaging surface available to contact the bone 5 of the vertebral bodies. For example, the surface projections can be cleaved by a longitudinal cut 540 generally parallel to the longitudinal axis L of implant 100 to form a surface projection having nine exposed sides. The surface projections may further be cleaved by a horizontal cut 542 generally perpendicular to the longitudinal axis L of implant 100 to form a surface projection having eighteen exposed sides. The cuts can penetrate the surface projection at a depth substantially equal to that of the height of the surface projections as measured from the upper or lower surfaces of 15 the implant. The cuts can be oriented along at least one of the longitudinal axis of the implant, an axis perpendicular to the longitudinal axis of said implant, and an axis at an angle between the longitudinal axis and the axis perpendicular to the longitudinal axis of the implant. It is appreciated that 20 cuts 540 and 542 may be formed as part of the molding process for forming the surface projections.

When cleaved by longitudinal cut 540 and horizontal cut 542, each of surface projections 522 has angled forward facet **524***a*, **524***b* directed at least in part toward leading end 25 502 of implant 100 and rearward facets 526a, 526b directed at least in part toward trailing end 504 of implant 100. Forward facet 524 has a length greater than the length of rearward facet **526**. Rearward facets **526**a, **526**b have a slope that is steeper than the slope of forward facets **524***a*, **524***b*. 30 The cleaved portion of surface projection 522 can be spaced apart by a predetermined distance and the space can be at least in part flat, curved, or any other surface configuration suitable for the intended use. In this embodiment, the base of rearward facets 526a, 526b forms an angle of approxi- 35 mately 45 degrees with respect to upper and/or lower surfaces 506, 508 of implant 100. Each one of surface projections 522 has left side facets 532a, 532b and right side facets 534a, 534b directed toward the sides of implant 100, and forward facets 524a, 524b and rearward facet 526a, 40 **526***b*. In this embodiment of surface configuration **520**, longitudinal grooves 530 have a V-shaped horizontal crosssection and each surface projection 522 has left and right side facets 532a, 532b, 534a, 534b that converge toward one another. The left and right side facets 532a, 532b, 534a, 45 534b resist side-to-side motion of implant 100 after it is inserted into the implantation space. The surface configuration of the present invention can be formed by molding, machining or otherwise. A preferred surface configuration of the present invention may readily be machined by milling 50 from side to side, across the upper and lower vertebral body engaging surfaces, surface projections. A milling machine with a cutting tool having an angled cutting face such as a V-shaped profile can then be run through the plurality of surface projections parallel to the longitudinal axis of the 55 implant to form the above-described surface. In a preferred embodiment, the V-shaped cutting tool of the milling machine has faces with an angle of approximately 90 degrees, which faces are at a 45 degree angle to the plane of the surfaces being so machined. Without departing from the 60 present invention, the angle of the cutting faces can be more or less than 90 degrees, the angle of the cutting face to the surface to be cut can be more or less than 45 degrees, and rather than running the cutter element parallel to the longitudinal axis of the implant, the cutting element may be run 65 at an angle. By way of example only and not limitation, this angle may be at 45 degrees to the longitudinal axis of the

10

implant and each surface projection can be formed by two grooves crossing the projections at a 90 degree angle to each other.

The spinal implants of the present invention are made of artificial or naturally occurring materials suitable for implantation in the human spine. The implants can comprise bone including, but not limited to, cortical bone, materials other than bone, such as metals including, but not limited to, titanium and its alloys or ASTM material, surgical grade plastics, plastic composites, ceramics, or other materials suitable for use as a spinal implant. The implants of the present invention can further comprise or be combined with bone growth promoting materials, including but not limited to, bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone. The implants can be treated with a bone growth promoting substance, can be a source of osteogenesis, or can be bioabsorbable at least in part. The implants of the present invention can be formed of a porous material.

The spinal implants of the present invention can be for the purpose of achieving fusion. The upper and lower surfaces of the fusion implants can include at least one opening, each in communication with the other, to permit for the growth of bone from vertebral body to adjacent vertebral body through the implant. The implant can have an internal chamber and may also have an access opening for accessing the internal chamber, in which case the implant can further have a cover such as a cap to close the access opening at least in part. Openings in the upper and lower surfaces of the implant can communicate with the internal chamber to permit further growth of bone from vertebral body to adjacent vertebral body through the implant. The internal chamber can contain bone growth promoting materials, including but not limited to, bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone. The implants of the present invention can be formed of a material that intrinsically participates in the growth of bone from one of the adjacent vertebral bodies to the other of the adjacent verte-

While various embodiments of the present invention are presented by way of example only and not limitation, common to each of them, is that the configuration of the surface is based on a plurality of surface projections disposed in arrays, each surface projection comprising at least one leading facet and at least one opposing trailing facet, in which the leading facet has a length greater than the trailing facet and the trailing facet has a steeper slope than the slope of the leading facet. The surface configuration is located on at least a portion of one of the opposed vertebral body engaging surfaces of the spinal implant.

While the implant shown in FIGS. 1, 2, and 3 is an interbody spinal fusion implant, it is appreciated that the surface configuration of the present invention is applicable to any interbody spinal fusion implants, including but not limited to, an artificial disc or motion preserving device having opposed surfaces incorporating the present inventive teachings for engaging each of the adjacent vertebral bodies.

It is believed that the operation and construction of the present invention will be apparent from the foregoing description and, while the invention shown and described herein has been characterized as particular embodiments, changes and modifications may be made therein without departing from the spirit and scope of the invention as defined in the following claims.

What is claimed is:

1. A method for forming an interbody spinal implant having a plurality of bone engaging structures for insertion

between adjacent vertebral bodies of a human spine, the method comprising the steps of:

providing the implant comprising a leading end for introduction of the spinal implant into the spine, an opposite trailing end, spaced apart sides therebetween, a midlongitudinal axis through the leading and trailing ends, and opposite upper and lower surfaces between said leading and trailing ends and said spaced apart sides, said upper surface adapted for placement toward the bone of one of the vertebral bodies and said opposite lower surface adapted for placement toward the bone of the other of the vertebral bodies when said implant is placed between the adjacent vertebral bodies; and

forming a plurality of surface projections as part of the upper and lower surfaces of the implant, each of the surface projections having a base with a length and a width, each of the surface projections having a height perpendicular to the length and the width, at least two of the surface projections each having at least one forward facet directed at least in part toward the leading 20 end and at least one rearward facet directed at least in part toward the trailing end, said rearward facet terminating at a first location proximate the base, said rearward facet terminating at a second location proximate said forward facing facet, the first location being 25 closer than the second location to a vertical plane transverse to the mid-longitudinal axis representing the furthest distal extent of the leading end of the implant. each of said at least two surface projections being formed to have opposed side facets extending from the base and being directed generally toward said spaced apart sides of the implant, respectively, said side facets converging toward one another between said forward facet and said rearward facet of each of said at least two surface projections, the width of the base of each of 35 said at least two surface projections being greater than the respective height of each of said at least two surface projections.

- 2. The method of claim 1, wherein the step of forming includes one of the sub-steps of grinding, milling, burning, lasering, burnishing, electric discharge machining, and broaching to form said surface projections.
- 3. The method of claim 1, wherein the steps of providing and forming include the sub-step of casting to form said implant with said surface projections.
- **4**. The method of claim **1**, wherein said forming step includes the sub-step of orienting said projections relative to one another to form an array.
- 5. The method of claim 1, wherein said forming step includes the sub-step of orienting said projections to be geometrically disposed relative to one another.
- **6.** The method of claim **1**, wherein the step of forming includes forming the forward facets of at least two of said at least two surface projections to face the same direction.
- 7. The method of claim 1, wherein the step of farming includes forming said at least two surface projections to have substantially the same maximum height from the surface of said implant.
- 8. The method of claim 1, wherein the step of providing 60 the implant includes providing an implant having at least one opening in each of the upper and lower surfaces in communication with one another, the openings being configured to permit for the growth of bone from vertebral body to adjacent vertebral body through the implant.
- 9. The method of claim 8, further comprising the step of combining the implant with at least one of harvested bone,

12

bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone.

- 10. The method of claim 8, wherein the step of providing the implant includes providing an implant having an internal chamber between the upper and lower surfaces and in communication with the at least one opening in each of the upper and lower surfaces, the internal chamber being adapted to contain bone growth promoting materials.
- 11. The method of claim 10, further comprising the step of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone.
- 12. The method of claim 1, wherein the step of forming the plurality of surface projections includes using a milling instrument.
- 13. The method of claim 12, wherein the milling instrument includes a cutting tool with a V-shaped profile.
- 14. The method of claim 1, wherein said forward facet and said rearward facet have a length and a slope, the step of forming including forming the length of said forward facet to be longer than the length of said rearward facet, and forming the slope of said rearward facet to be steeper than the slope of said forward facet.
- 15. The method of claim 1, wherein the step of providing includes providing the implant with a width and a height, the width being greater than the height.
- 16. The method of claim 1, wherein the step of forming includes forming each of the forward facet, the rearward facet, and the opposed side facets to be in the shape of a triangle.
- 17. The method of claim 1, wherein the step of forming includes forming the forward facet to be in the shape of an isosceles triangle.
- 18. The method of claim 1, wherein the step of providing includes providing the implant so that the upper and lower surfaces each include a plurality of openings to permit bone growth from adjacent vertebral body to adjacent vertebral body through the implant.
- 19. A method for forming an interbody spinal implant having an exterior surface with a plurality of bone engaging structures for insertion between adjacent vertebral bodies of a human spine, the method comprising the steps of:
 - providing the implant comprising a leading end for introduction of the spinal implant into the spine, an opposite trailing end, spaced apart sides therebetween, and a mid-longitudinal axis passing through the leading and trailing ends, opposite upper and lower surfaces between said leading and trailing ends and said spaced apart sides, said upper surface adapted for placement toward the bone of one of the vertebral bodies and said opposite lower surface adapted for placement toward the bone of the other of the vertebral bodies when the implant is placed between the adjacent vertebral bodies; and
 - forming surface projections as part of the upper and lower surfaces of the implant:
 - at least a first and second of said surface projections each having a first facet configuration with at least one forward facet directed at least in part toward the leading end and at least one rearward facet directed at least in part toward the trailing end, said forward facet and said rearward facet having a length and a slope, the length of said forward facet being longer than the length of said rearward facet, the slope of said rearward facet being steeper than the slope of said forward facet, said first and second surface projections having opposed side facets between said forward facet and said rear-

ward facet, said side facets having at least a first portion in a plane passing through and being at an angle to the mid-longitudinal axis of the implant, said first and second surface projections each having a peak along a first line that is transverse to the mid-longitudinal axis 5 of said implant; and

- at least a third and fourth of said surface projections each having a second facet configuration with at least one forward facet directed at least in part toward the leading end and at least one rearward facet directed at least in 10 part toward the trailing end, said forward facet and said rearward facet of said second facet configuration having a length and a slope, the length of said forward facet of said second facet configuration being longer than the length of said rearward facet of said second facet 15 configuration, the slope of said rearward facet of said second facet configuration being steeper than the slope of said forward facet of said second facet configuration, said third and fourth surface projections each having a peak along a second line that is transverse to the 20 mid-longitudinal axis and off-set from the first line transverse to the mid-longitudinal axis, the second facet configuration of the third and fourth surface projections being different from the first facet configuration of the first and second surface projections.
- 20. The method of claim 19, wherein the step of forming includes one of the sub-steps of grinding, milling, burning, lasering, burnishing, electric discharge machining, broaching, and machining to form said surface projections.
- 21. The method of claim 19, wherein the steps of pro- 30 viding and forming include the sub-step of casting to form said implant with said surface projections.
- 22. The method of claim 19, wherein said forming step includes the sub-step of orienting said projections relative to one another to form an array.
- 23. The method of claim 19, wherein said forming step includes the cub-step of orienting said projections to be geometrically disposed relative to one another.
- 24. The method of claim 19, wherein the step of providing one opening in each of the upper and lower surfaces in communication with one another, the openings being configured to permit for the growth of bone from vertebral body to adjacent vertebral body through the implant.
- 25. The method of claim 24, further comprising the step 45 of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone.
- 26. The method of claim 24, wherein the step of providing the implant includes providing an implant having an internal 50 chamber between the upper and lower surfaces and in communication with the at least one opening in each of the upper and lower surfaces, the internal chamber being adapted to contain bone growth promoting materials.
- 27. The method of claim 26, further comprising the step 55 of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone.
- 28. The method of claim 19, wherein the step of forming the plurality of surface projections includes using a milling 60
- 29. The method of claim 19, wherein at least a fifth and sixth of said surface projections formed during the step of forming each have a third facet configuration with at least one forward facet directed at least in part toward the leading 65 end and at least one rearward facet directed at least in part toward the trailing end, said forward facet and said rearward

14

facet of said third facet configuration having a length and a slope, the length of said forward facet of said third facet configuration being longer than the length of said rearward facet of said third facet configuration, the slope of said rearward facet of said third facet configuration being steeper than the slope of said forward facet of said third facet configuration, said fifth and sixth surface projections each having a peak along a third fine that is transverse to the mid-longitudinal axis and off-set from the first and second lines, the third facet configuration of the fifth and sixth surface projections being different from the first facet configuration of the first and second surface projections and the second facet configuration of the third and fourth surface projections.

- 30. The method of claim 19, wherein the step of farming includes forming said forward facets of said first and second surface projections to face the same direction.
- 31. The method of claim 19, wherein the step of forming includes forming at feast one of the surface projections along the first line to have a maximum height from the surface of the implant that is substantially the same as the maximum height of one of the surface projections along the second line.
- 32. The method of claim 19, wherein the step of forming 25 includes forming at least three surface projections having the first facet configuration along the first line and forming at least three surface projections having the second facet configuration along the second line.
 - 33. The method of claim 19, wherein the step of forming includes forming at least four surface projections having the first facet configuration along the first line and forming at least four surface projections having the second facet configuration along the second line.
 - 34. The method of claim 19, wherein the step of forming includes forming at least five surface projections having the first facet configuration along the first line and forming at least five surface projections having the second facet configuration along the second line.
- 35. A method for forming an interbody spinal implant the implant includes providing an implant having at least 40 having an exterior surface with a plurality of bone engaging structures for insertion between adjacent vertebral bodies of a human spine, the method comprising the steps of:
 - providing the implant comprising a leading end for introduction of the spinal implant into the spine, an opposite trailing end, spaced apart sides therebetween, and a mid-longitudinal axis passing through the leading and trailing ends, opposite upper and lower surfaces between said leading and trailing ends and said spaced apart sides, said upper surface adapted for placement toward the bone of one of the vertebral bodies and said opposite lower surface adapted for placement toward the bone of the other of the vertebral bodies when the implant is placed between the adjacent vertebral bodies: and
 - forming surface projections as part of the upper and lower surfaces of the implant with a milling instrument, at least two of said surface projections each having at least one forward facet directed at least in part toward the leading end and at least one rearward facet directed at feast in part toward the trailing end, said forward facet and said rearward facet having a length and a slope, the length of said forward facet being longer than the length of said rearward facet the slope of said rearward facet being steeper than the slope of said forward facet said at least two of said surface projections having opposed side facets between said forward facet and said rearward facet, said side facets having at

least a first portion in a plane passing through and being at an angle to the mid-longitudinal axis of the implant, said forward facets of said at least two of said surface projections facing the same direction, wherein the milling instrument includes a culling tool with a 5 V-shaped profile.

- **36.** The method of claim **35,** wherein the step of providing the implant includes providing an implant having at least one opening in each of the upper and lower surfaces in communication with one another, the openings being configured to permit for the growth of bone from vertebral body to adjacent vertebral body through the implant.
- 37. The method of claim 36, further comprising the step of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and 15 genes coding for the production of bone.
- 38. The method of claim 36, wherein the step of providing the implant includes providing an implant having an internal chamber between the upper and lower surfaces and in communication with the at least one opening in each of the 20 upper and lower surfaces, the internal chamber being adapted to contain bone growth promoting materials.
- 39. The method of claim 38, further comprising the step of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and 25 genes coding for the production of bone.
- **40**. The method of claim **35**, wherein the forming step includes the sub-step of orienting said projections relative to one another to form an array.
- **41**. The method of claim **35**, wherein the forming step ³⁰ includes the sub-step of orienting said projections to be geometrically disposed relative to one another.
- **42**. A method for forming an Interbody spinal implant having an exterior surface with a plurality of bone engaging structures for insertion between adjacent vertebral bodies of 35 a human spine, the method comprising the steps of:
 - providing the implant comprising a leading end for introduction of the spinal implant into the spine, an opposite trailing end, spaced apart sides therebetween, and a mid-longitudinal axis passing through the leading and 40 trailing ends, opposite upper and lower surfaces between said leading and trailing ends and said spaced apart sides, said upper surface adapted for placement toward the bone of one of the vertebral bodies and said opposite lower surface adapted for placement toward 45 the bone of the other of the vertebral bodies when the implant is placed between the adjacent vertebral bodies: and

forming surface projections as part of the upper and lower surfaces of the implant

at least a first and second of said surface projections each having a first facet configuration with at least one forward facet directed at least in part toward the leading end and at least one rearward facet directed at least in part toward the trailing end, said forward facet and said 55 rearward facet having a length and a slope, the length of said forward facet being longer than the length of said rearward facet, the slope of said rearward facet being steeper than the slope of said forward facet, said first and second surface projections each having a peak 60 along a first line that is transverse to the mid-longitudinal axis of said implant; and

at least a third and fourth of said surface projections each having a second facet configuration with at least one forward facet directed at least in part toward the leading 65 end and at least one rearward facet directed at least in part toward the trailing end, said forward facet and said

16

rearward facet of said second facet configuration having a length and a slope, the length of said forward facet of said second facet configuration being longer than the length of said rearward facet of said second facet configuration, the slope of said rearward facet of said second facet configuration being steeper than the slope of said forward facet of said second facet configuration, said third and fourth surface projections each having a peak along a second line that is transverse to the mid-longitudinal axis and off-set from the first line transverse to the mid-longitudinal axis, the second facet configuration of the third and fourth surface projections being different from the first facet configuration of the first and second surface projections.

- **43**. The method of claim **42**, wherein the step of forming includes one of the sub-steps of grinding, milling, burning, lasering, burnishing, electric discharge machining, broaching, and machining to form said surface projections.
- **44**. The method of claim **42**, wherein the steps of providing and forming include the sub-step of casting to form said implant with said surface projections.
- **45**. The method of claim **42**, wherein the forming step includes the sub-step of orienting said projections relative to one another to form an array.
- **46**. The method of claim **42**, wherein the forming step includes the sub-step of orienting said projections to be geometrically disposed relative to one another.
- 47. The method of claim 42, wherein the step of providing the implant includes providing an implant having at least one opening in each of the upper and lower surfaces in communication with one another, the openings being configured to permit for the growth of bone from vertebral body to adjacent vertebral body through the implant.
- **48**. The method of claim **47**, further comprising the step of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone.
- 49. The method of claim 47, wherein the step of providing the implant includes providing an implant having an internal chamber between the upper and lower surfaces and in communication with the at least one opening in each of the upper and lower surfaces, the internal chamber being adapted to contain bone growth promoting materials.
- toward the bone of one of the vertebral bodies and said opposite lower surface adapted for placement toward the bone of the other of the vertebral bodies when the implant is placed between the adjacent vertebral bodies when the implant is placed between the adjacent vertebral bodies when the implant is placed between the adjacent vertebral bodies and said of combining the implant with at least one of harvested bone, bone morphogenetic proteins, hydroxyapatite, and genes coding for the production of bone.
 - **51**. The method of claim **42**, wherein the step of forming the plurality of surface projections includes using a milling 50 instrument.
 - 52. The method of claim 42, wherein at least a fifth and sixth of said surface projections formed during the step of forming each have a third facet configuration with at least one forward facet directed at least in part toward the leading end and at least one rearward facet directed at least in part toward the trailing end, said forward facet and said rearward facet of said third facet configuration having a length and a slope, the length of said forward facet of said third facet configuration being longer than the length of said rearward facet of said third facet configuration, the slope of said rearward facet of said third facet configuration being steeper than the slope of said forward facet of said third facet configuration, said fifth and sixth surface projections each having a peak along a third line that is transverse to the mid-longitudinal axis and off-set from the first and second lines, the third facet configuration of the fifth and sixth surface projections being different from the first facet con-

figuration of the first and second surface projections and the second facet configuration of the third and fourth surface projections.

- **53**. The method of claim **42**, wherein the step of forming includes forming said forward facets of said first and second 5 surface projections to face the same direction.
- **54**. The method of claim **42**, wherein the step of forming includes forming at least one of the surface projections along the first line to have a maximum height from the surface of the implant that is substantially the same as the maximum 10 height of one of the surface projections along the second line.
- 55. The method of claim 42, wherein the step of forming includes forming at least three surface projections having the first facet configuration along the first line and forming at

18

least three surface projections having the second facet configuration along the second line.

- **56**. The method of claim **42**, wherein the step of forming includes forming at least four surface projections having the first facet configuration along the first line and forming at least four surface projections having the second facet configuration along the second line.
- 57. The method of claim 42, wherein the step of forming includes forming at least five surface projections having the first facet configuration along the first line and forming at least five surface projections having the second facet configuration along the second line.

* * * * *

UNITED STATES PATENT AND TRADEMARK OFFICE CERTIFICATE OF CORRECTION

PATENT NO. : 7,166,129 B2 Page 1 of 1

APPLICATION NO.: 09/921851
DATED: January 23, 2007
INVENTOR(S): Gary K. Michelson

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

Title Page 2:

Item (56), U.S. Patent Documents, after "6,432,106 B1 8/2002 Fraser" insert the following:

-- 6,482,233B1 11/2002 Aebi --.

Column 11, line 55:

Change "farming" to -- forming --.

<u>Column 13, line 37</u>:

Change "cub-step" to -- sub-step --.

Column 14:

Line 8: change "fine" to -- line --;

Line 15: change "farming" to -- forming --; Lines 19 and 60: change "feast" to -- least --; and Lines 63 and 65: change "facet" to -- facet, --.

Column 15:

Line 5: change "culling" to -- cutting --;

Line 33: change "Interbody" to -- Interbody --; and

Line 50: change "implant" to -- implant: --.

Signed and Sealed this

Twenty-seventh Day of March, 2007

JON W. DUDAS
Director of the United States Patent and Trademark Office